



Minimalism with materials

Stress less

Practical

Efficiency

Dr.SM Chung



Pre-op (2023-09-06)



Post-op (2023-09-06)




Final prosthesis (2023-11-02)






#21

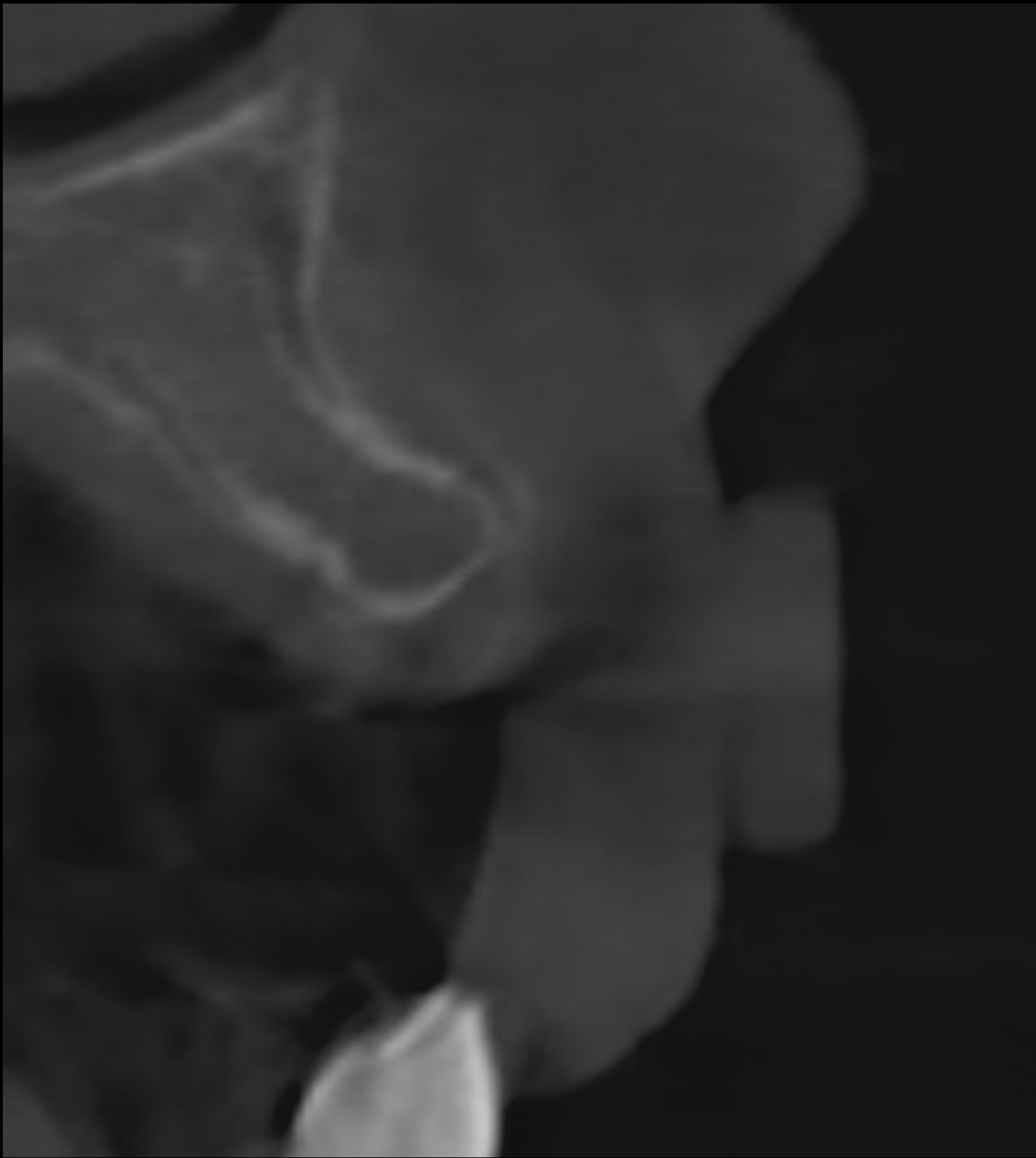

#22

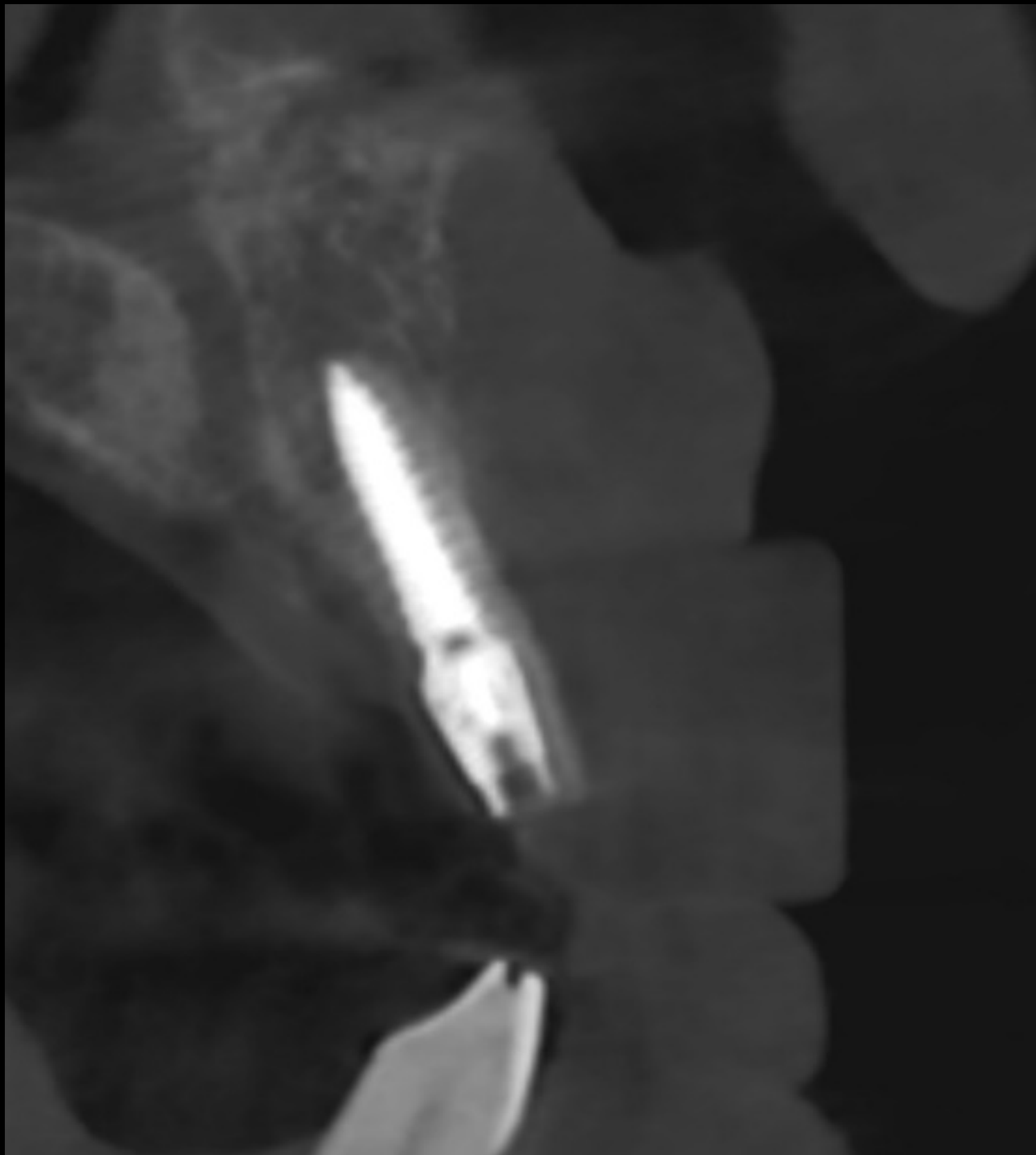

#21



#22





Pre-op (2023-09-06)






Surgery (2023-09-06)

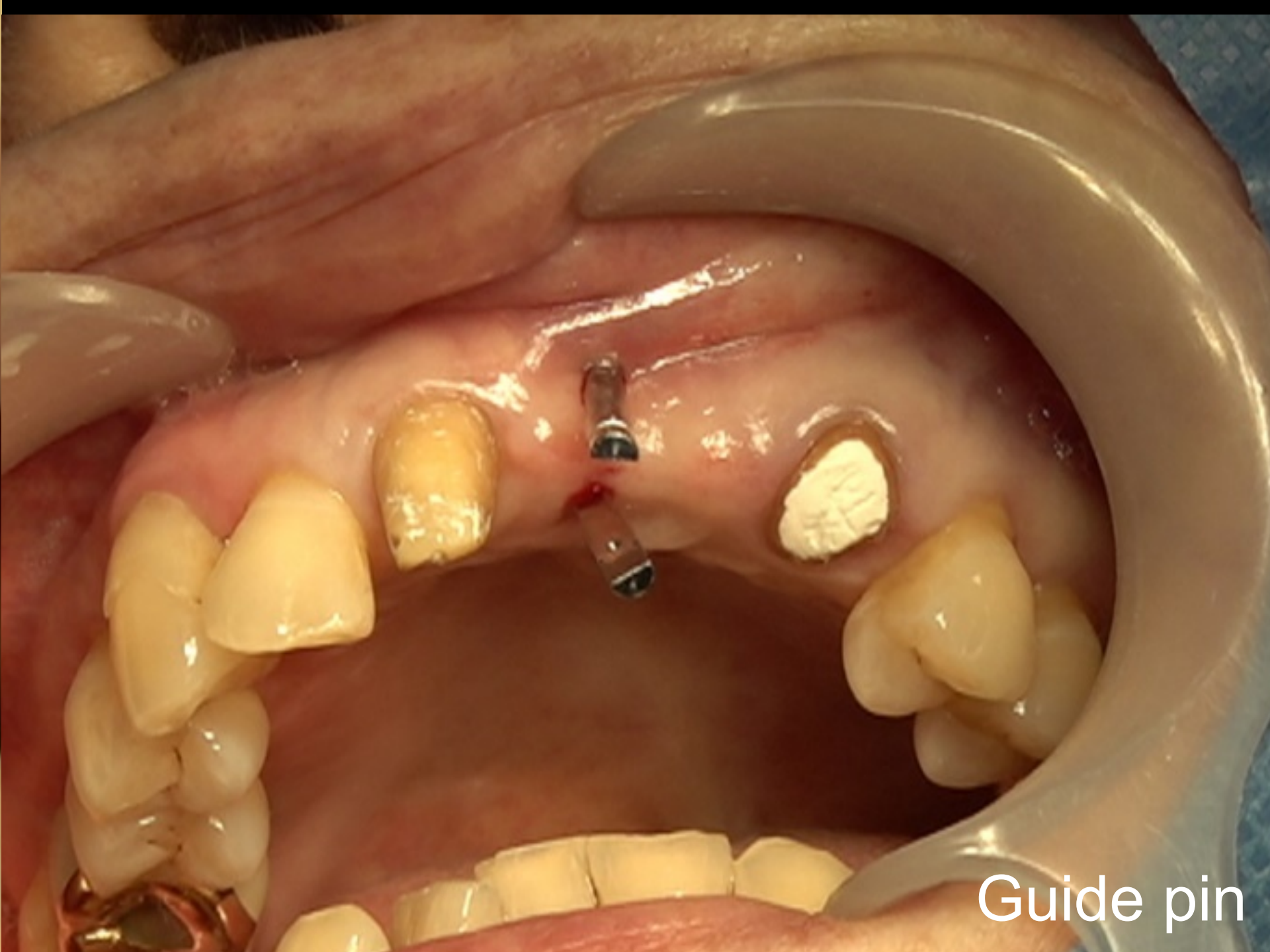
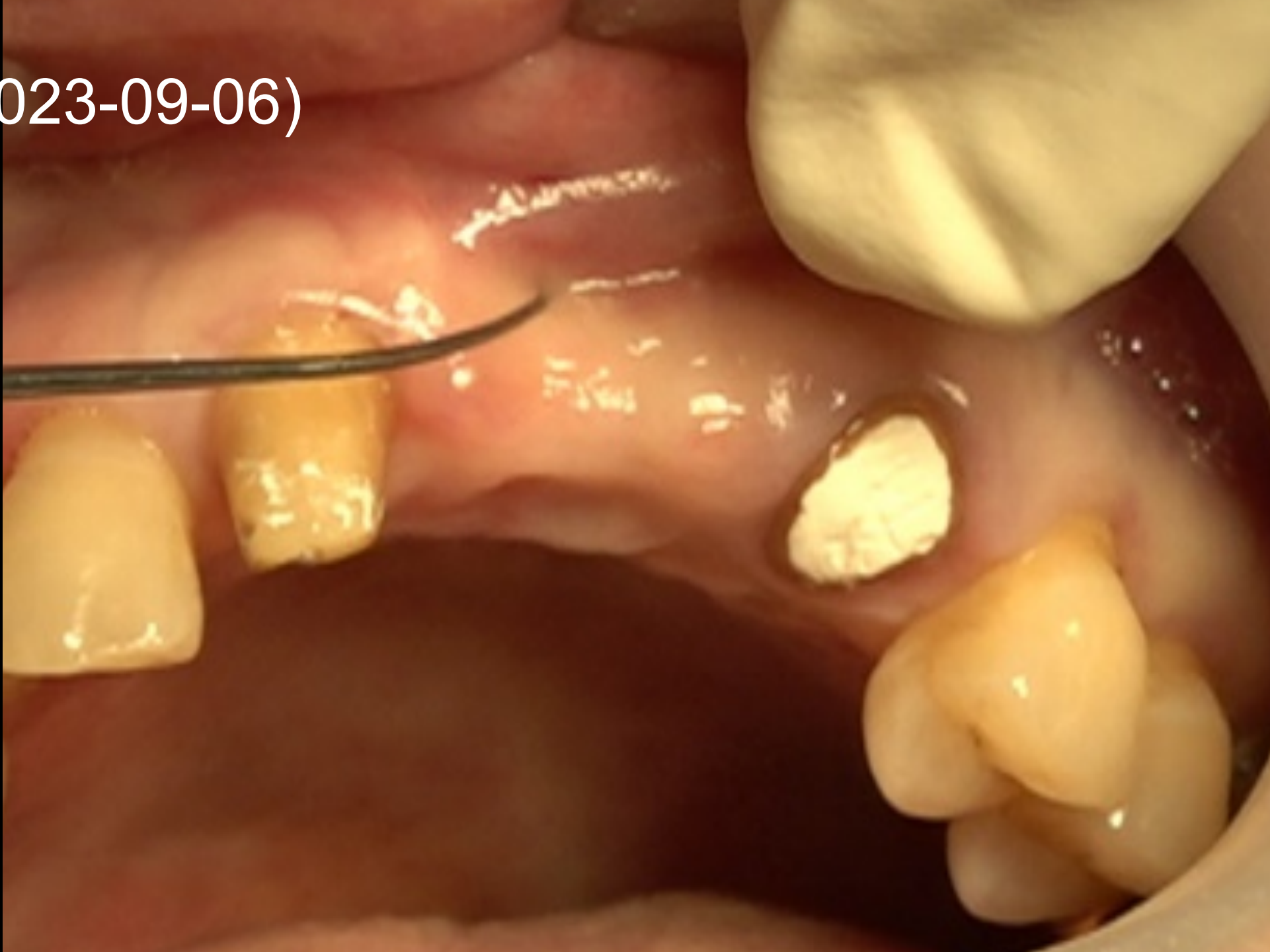


Guide pin

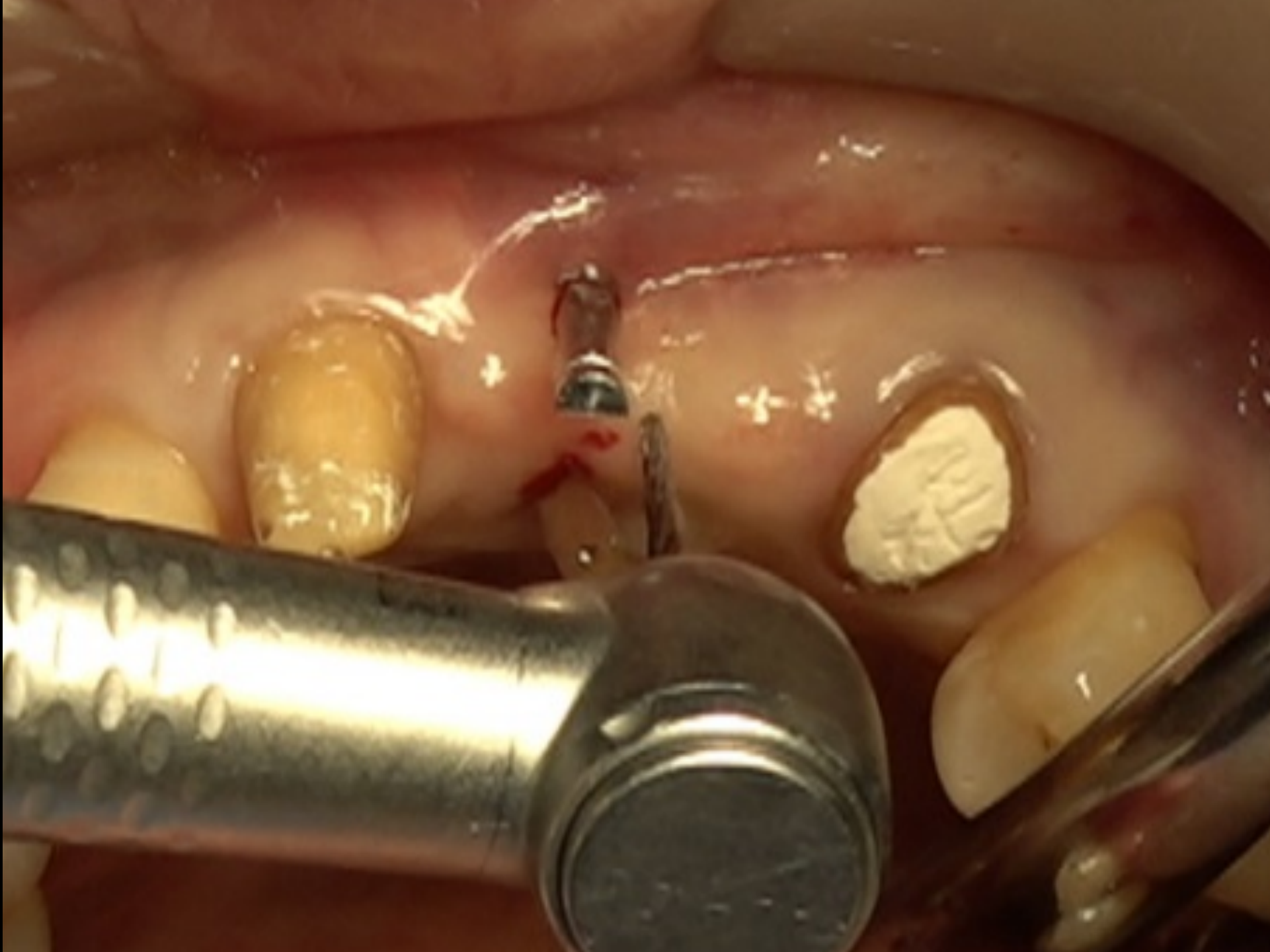




Surgery (2023-09-06)

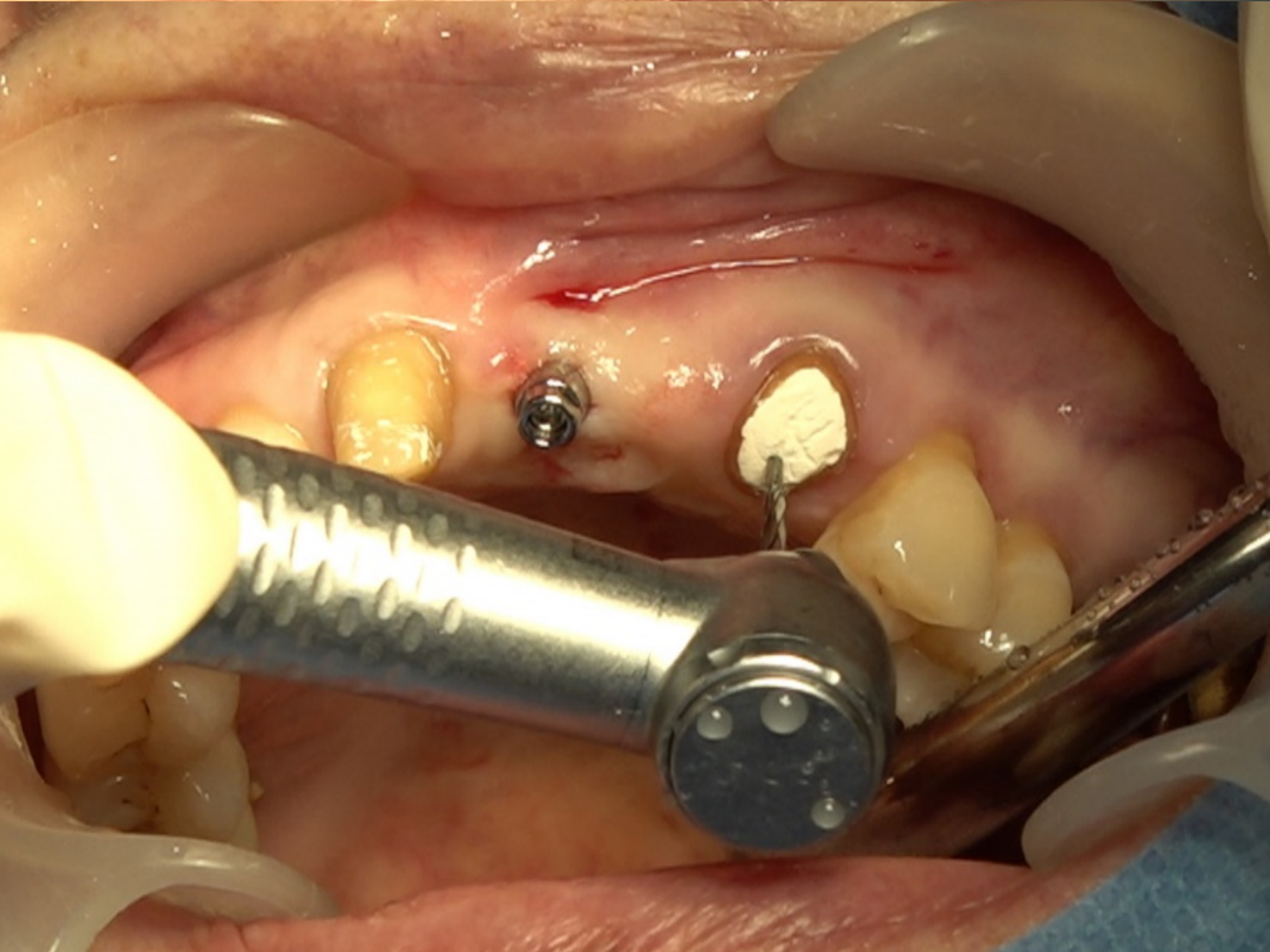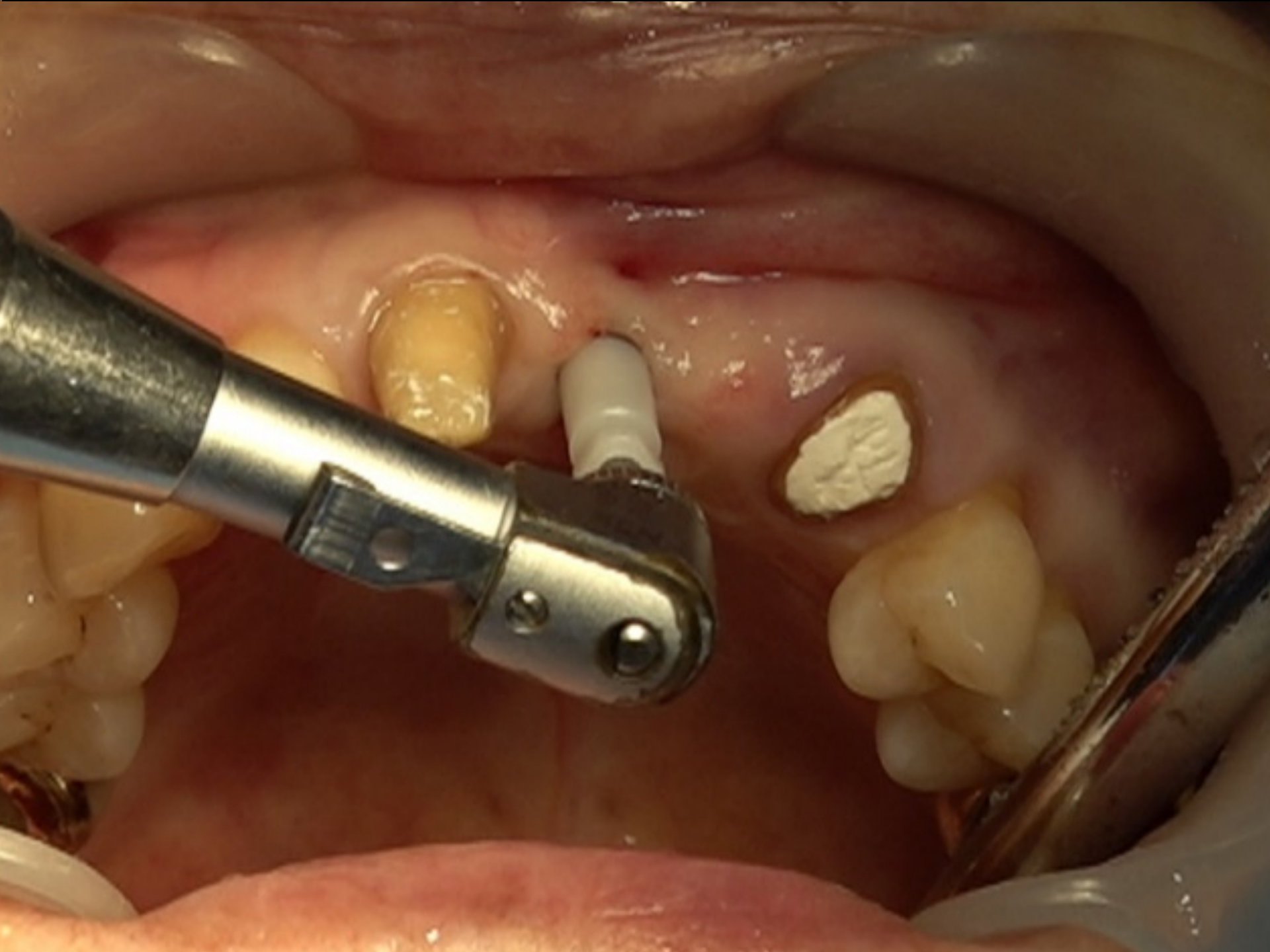




Surgery (2023-09-06)

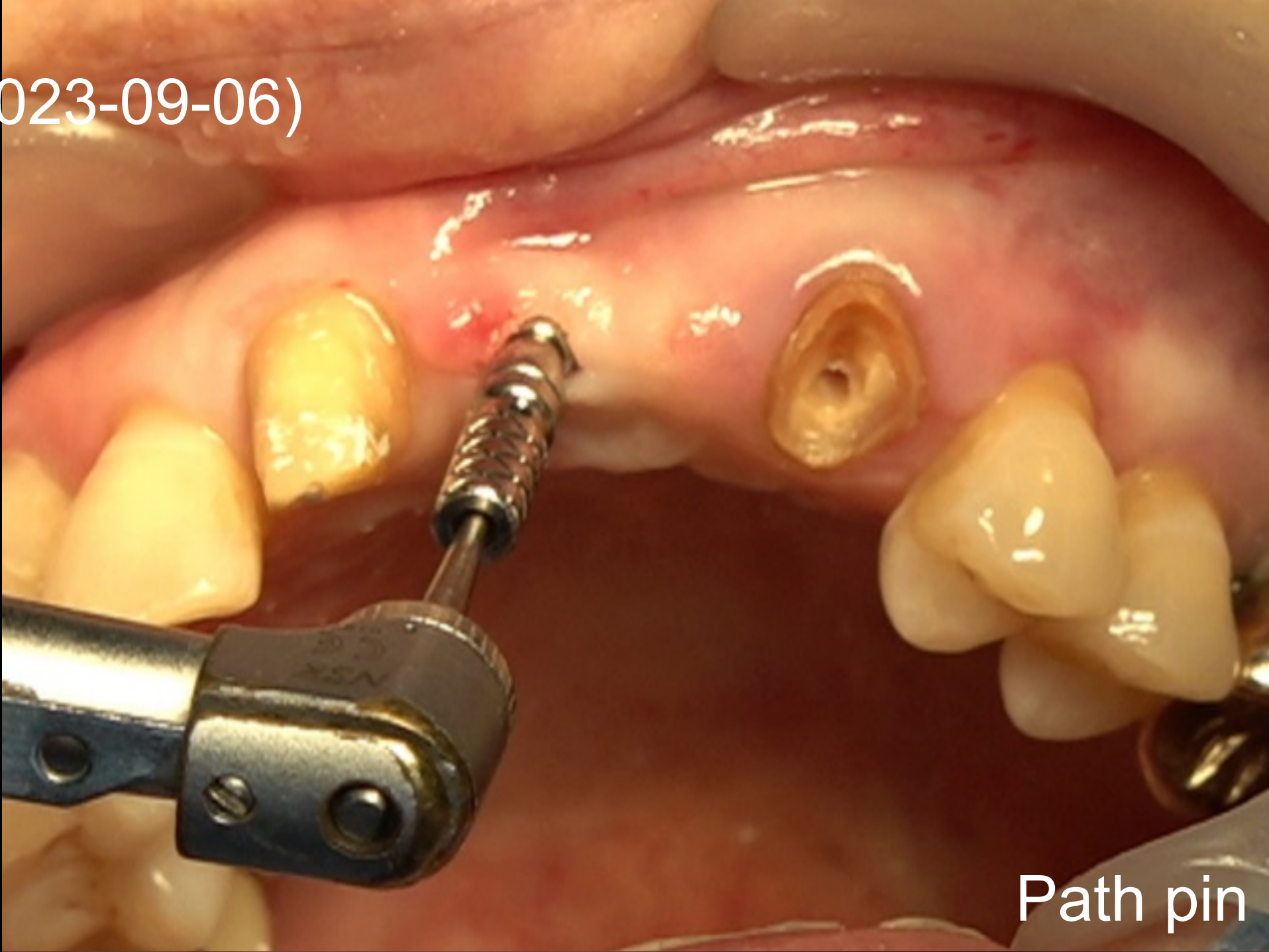

Path pin







Parallel pin

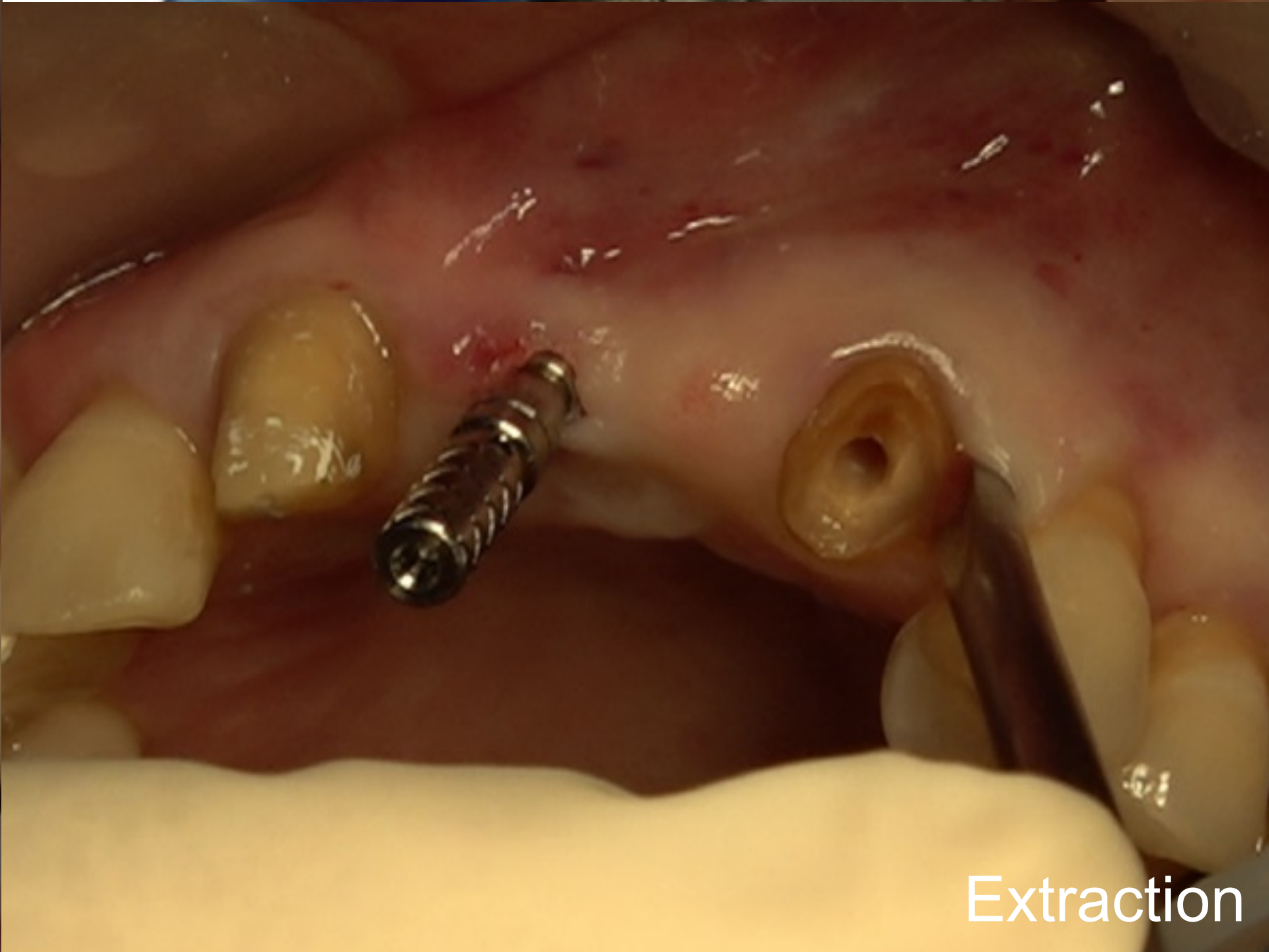

Extraction






Surgery (2023-09-06)

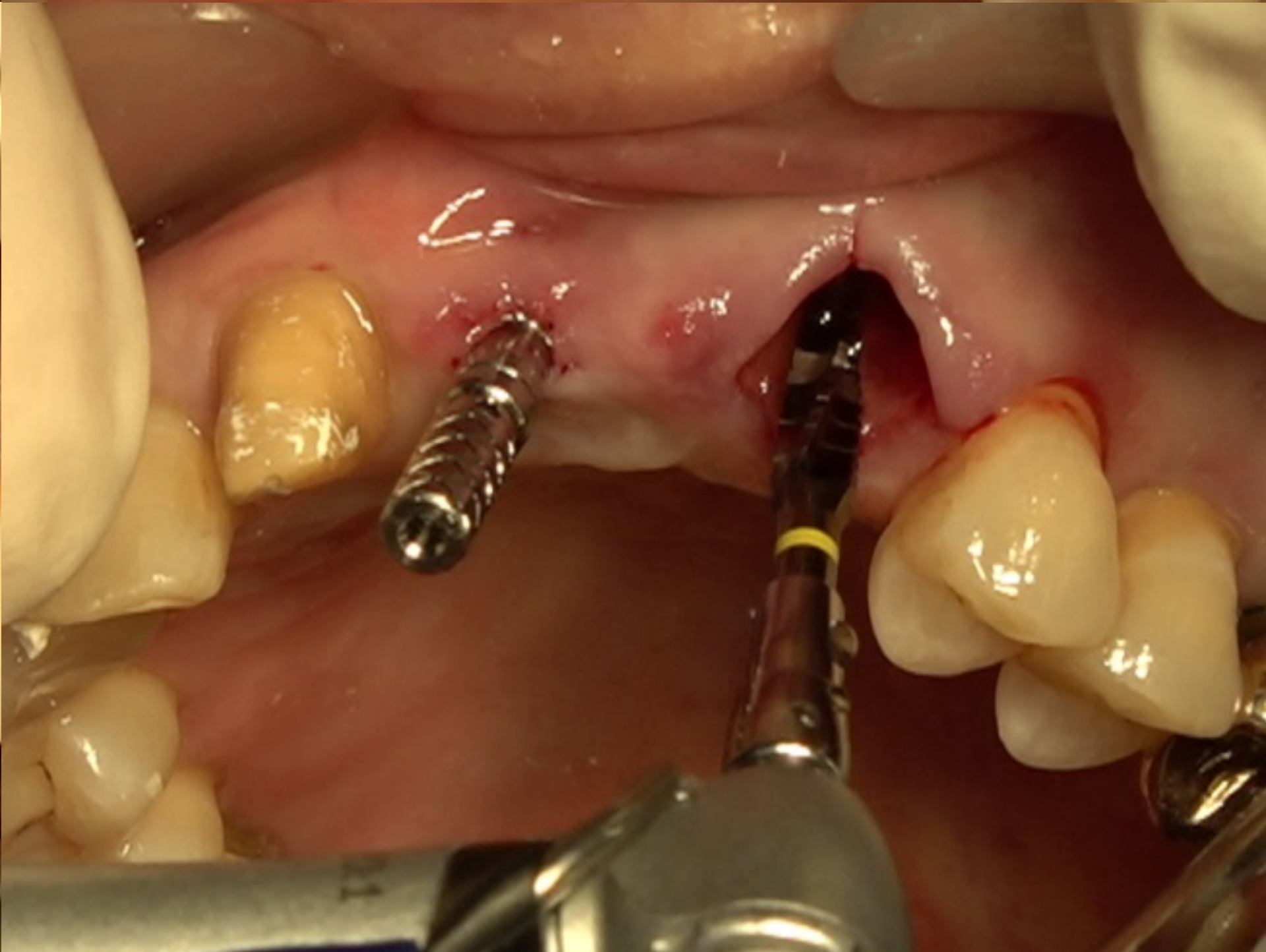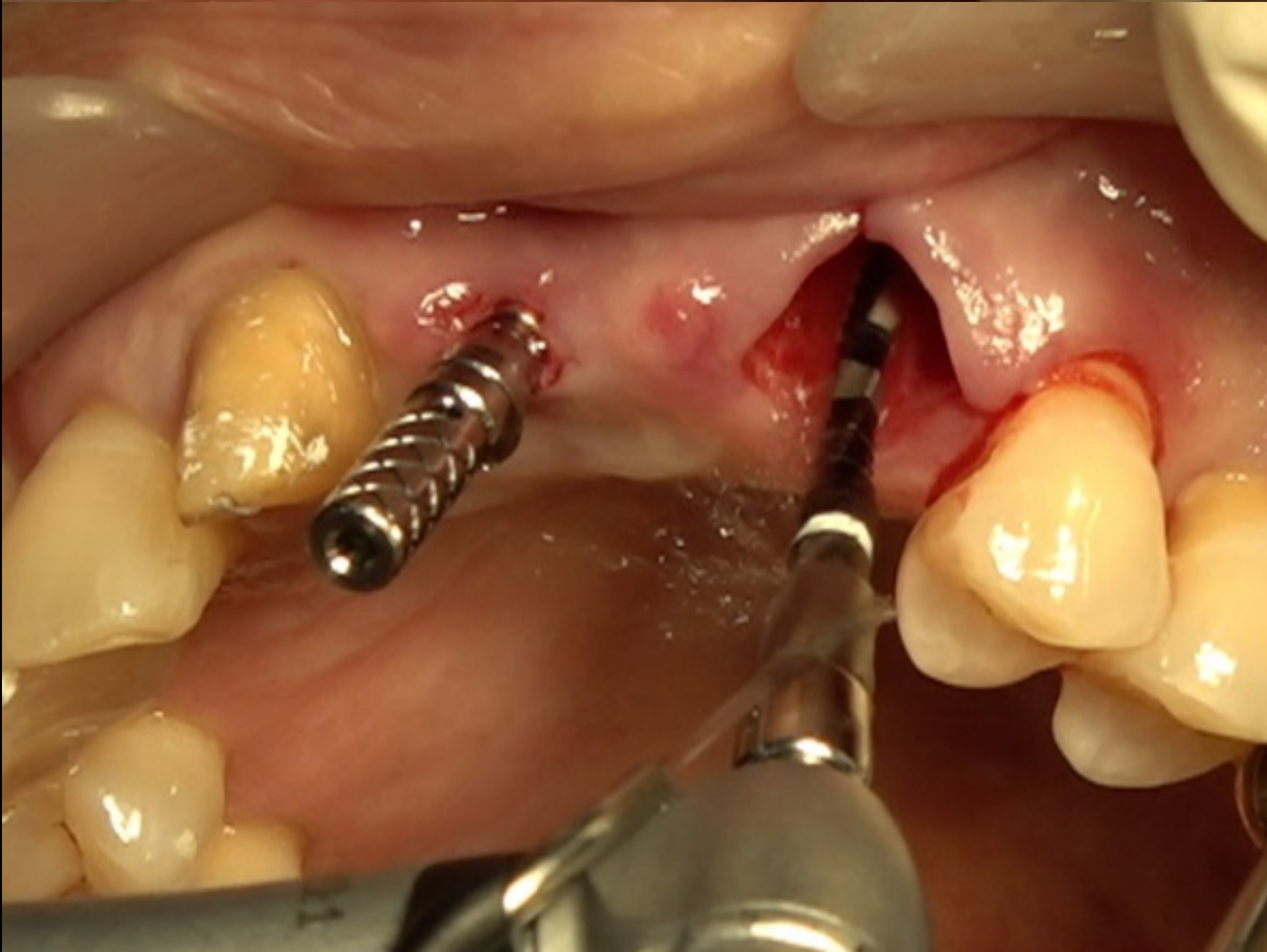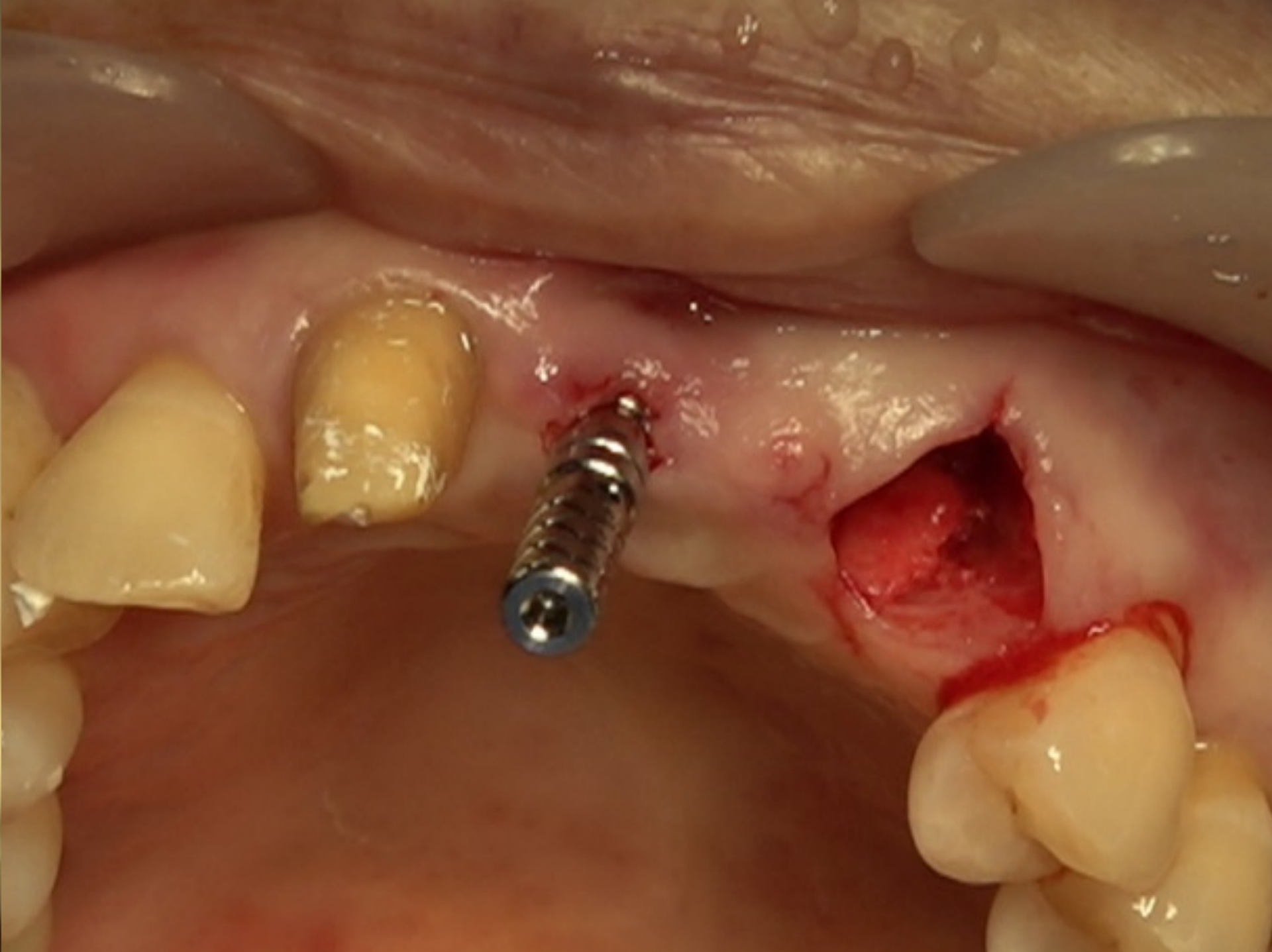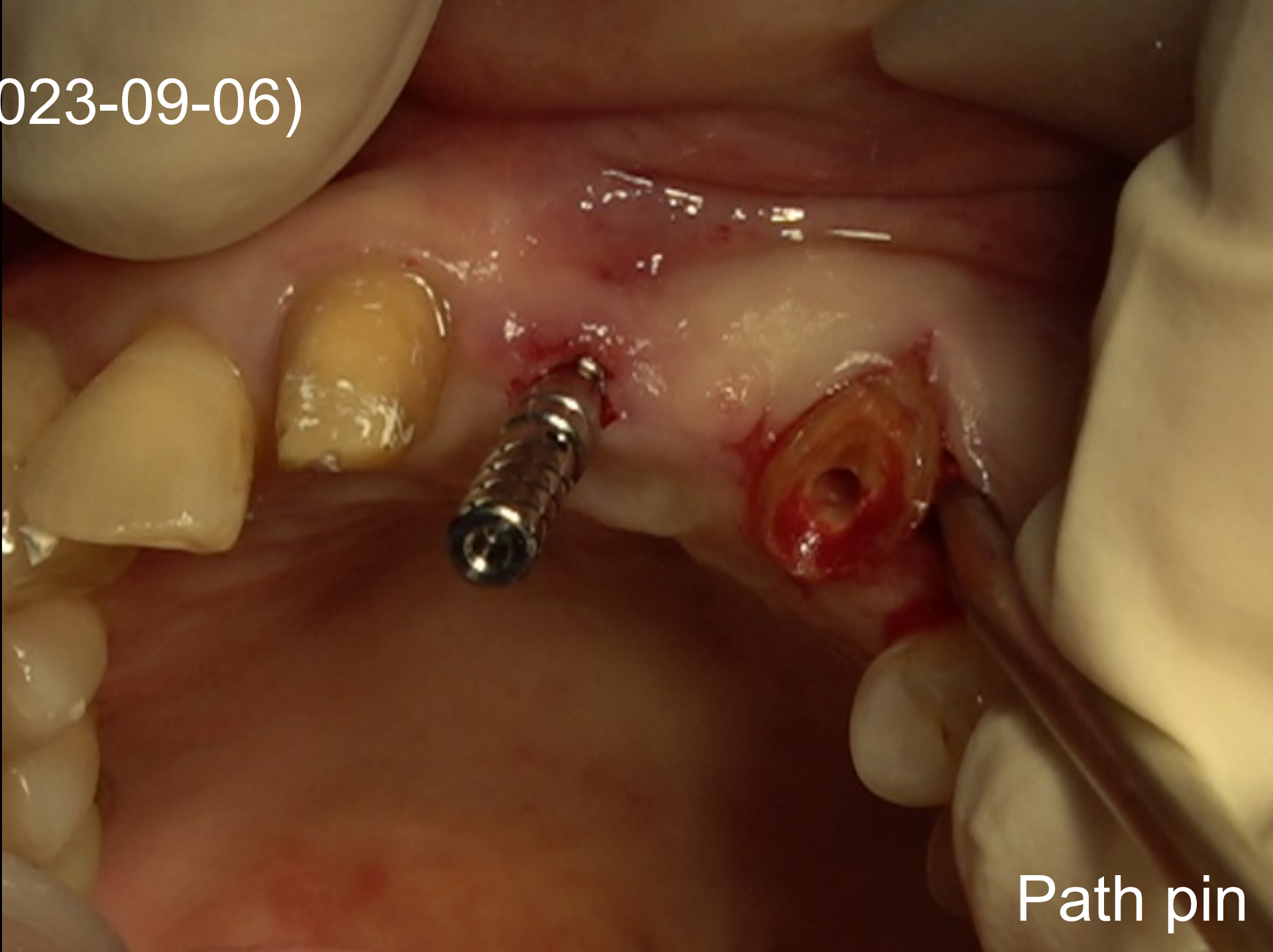





Surgery (2023-09-06)

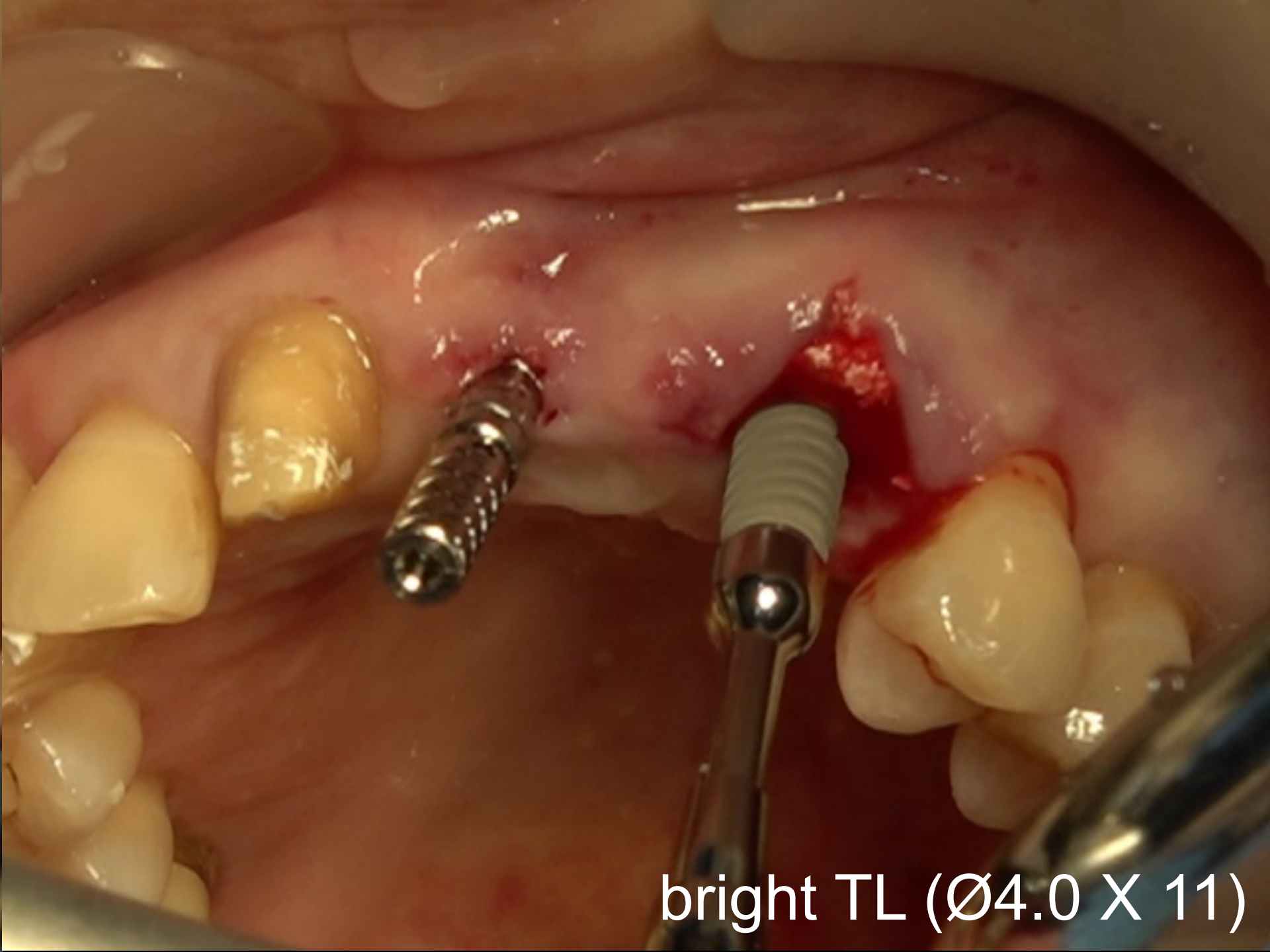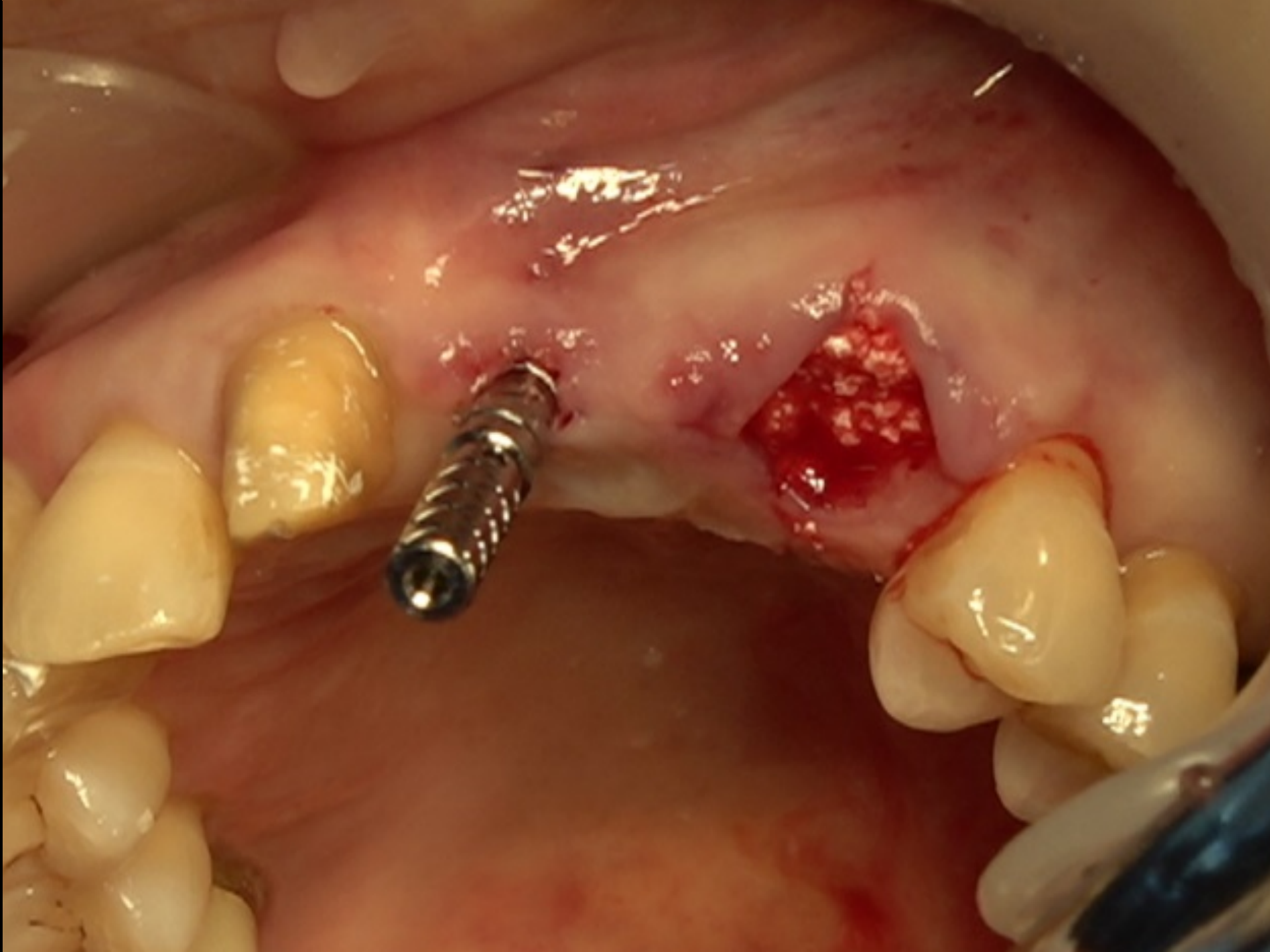
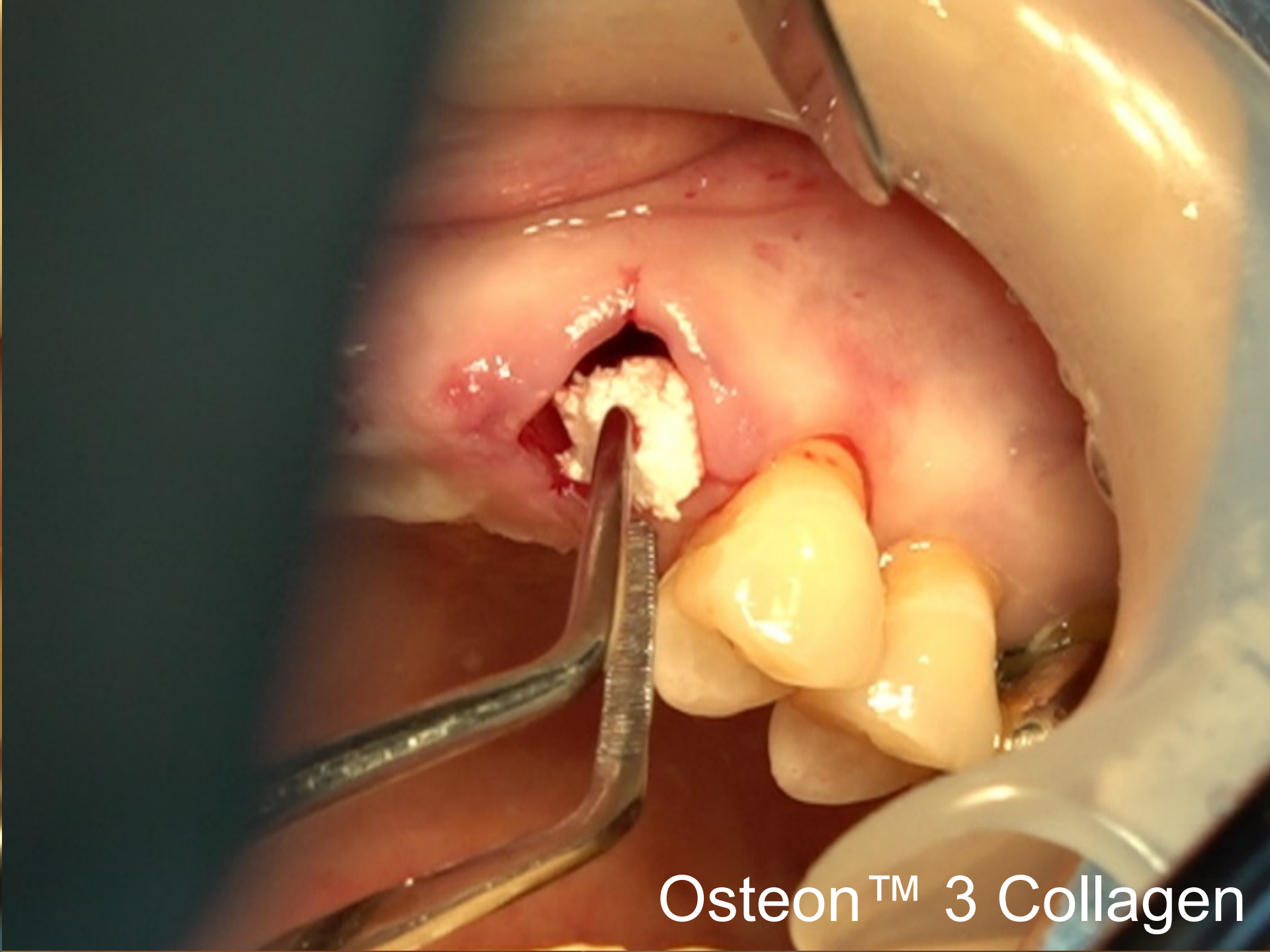




Surgery (2023-09-06)

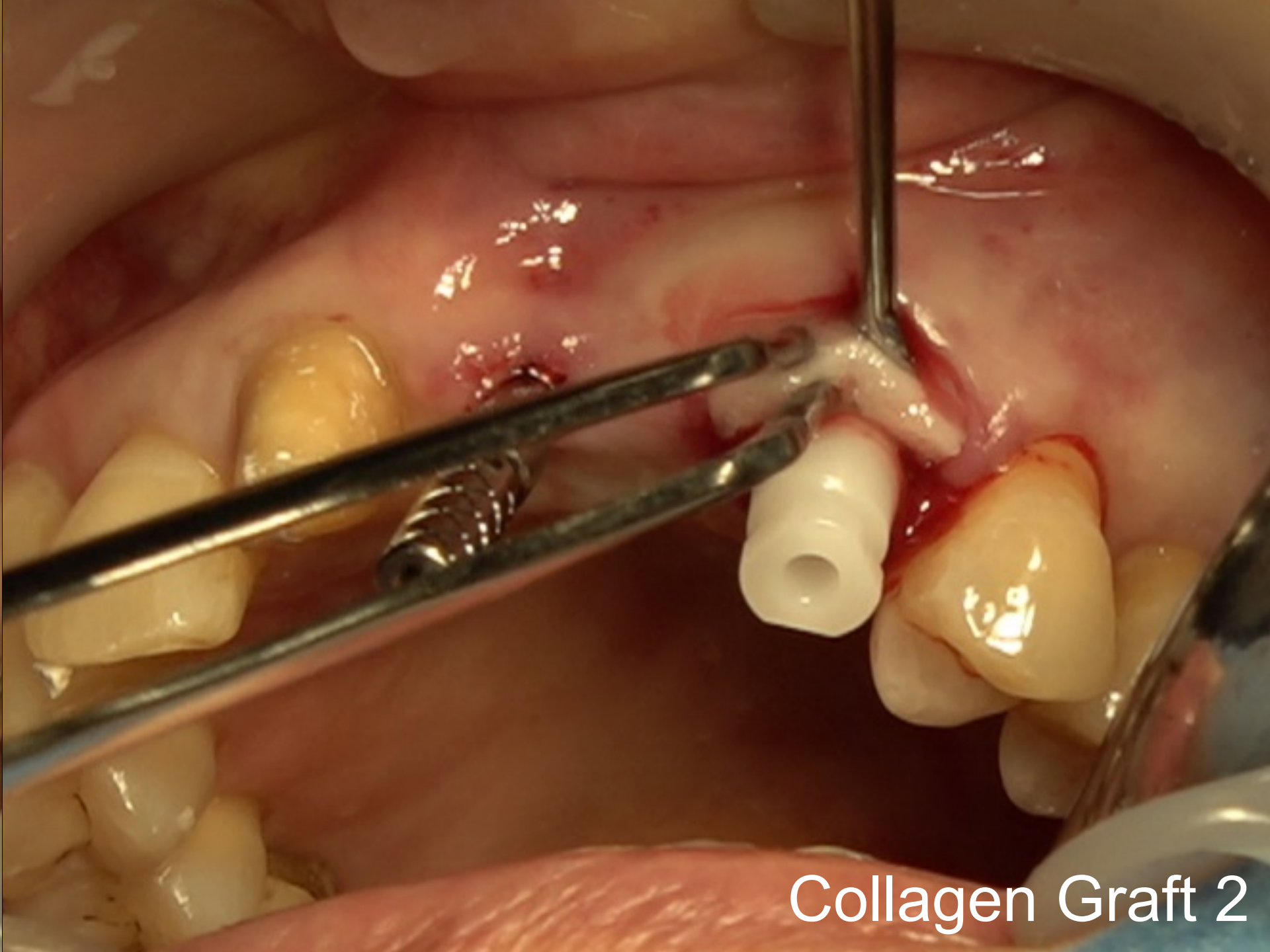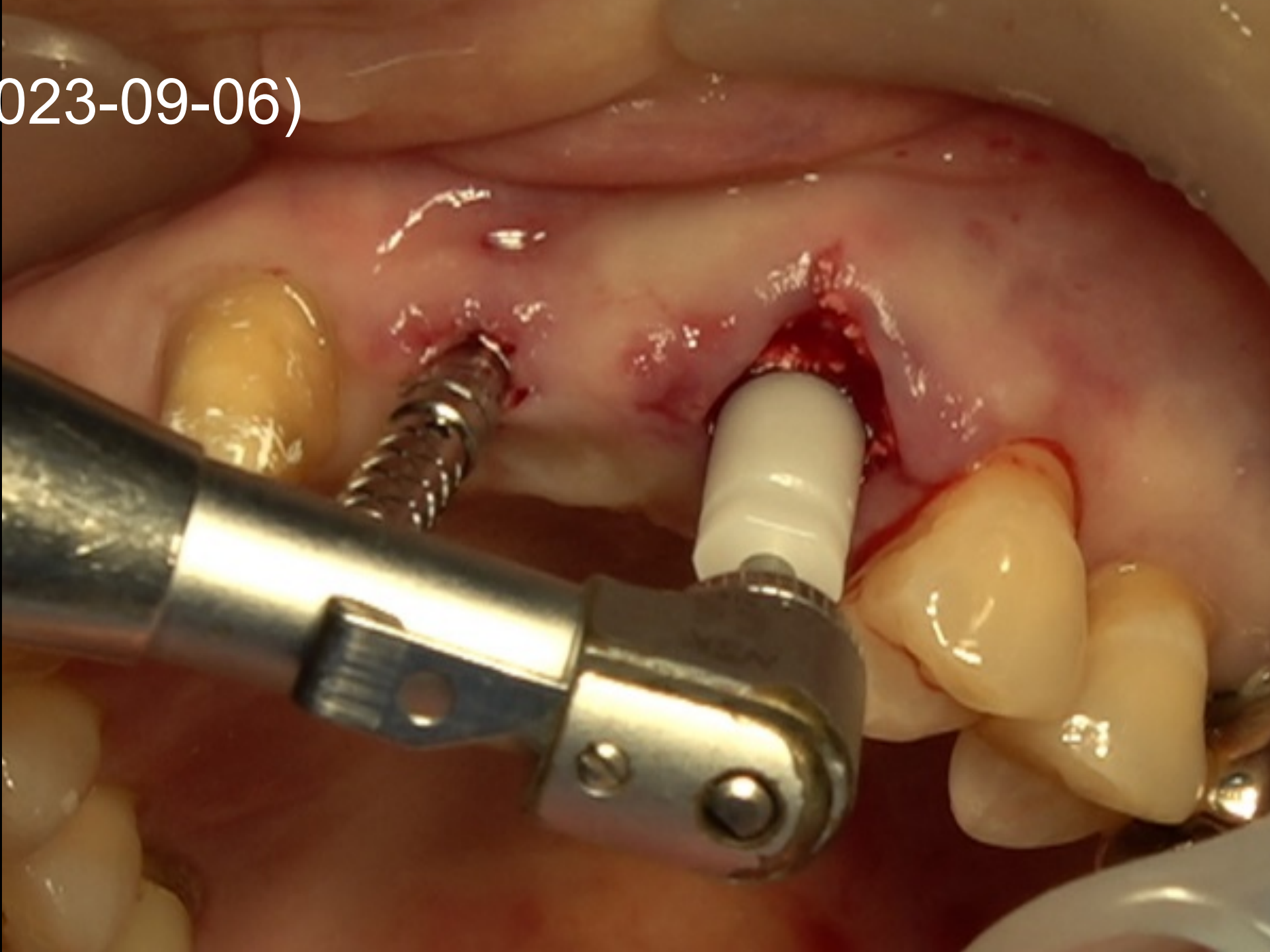


Collagen Graft 2

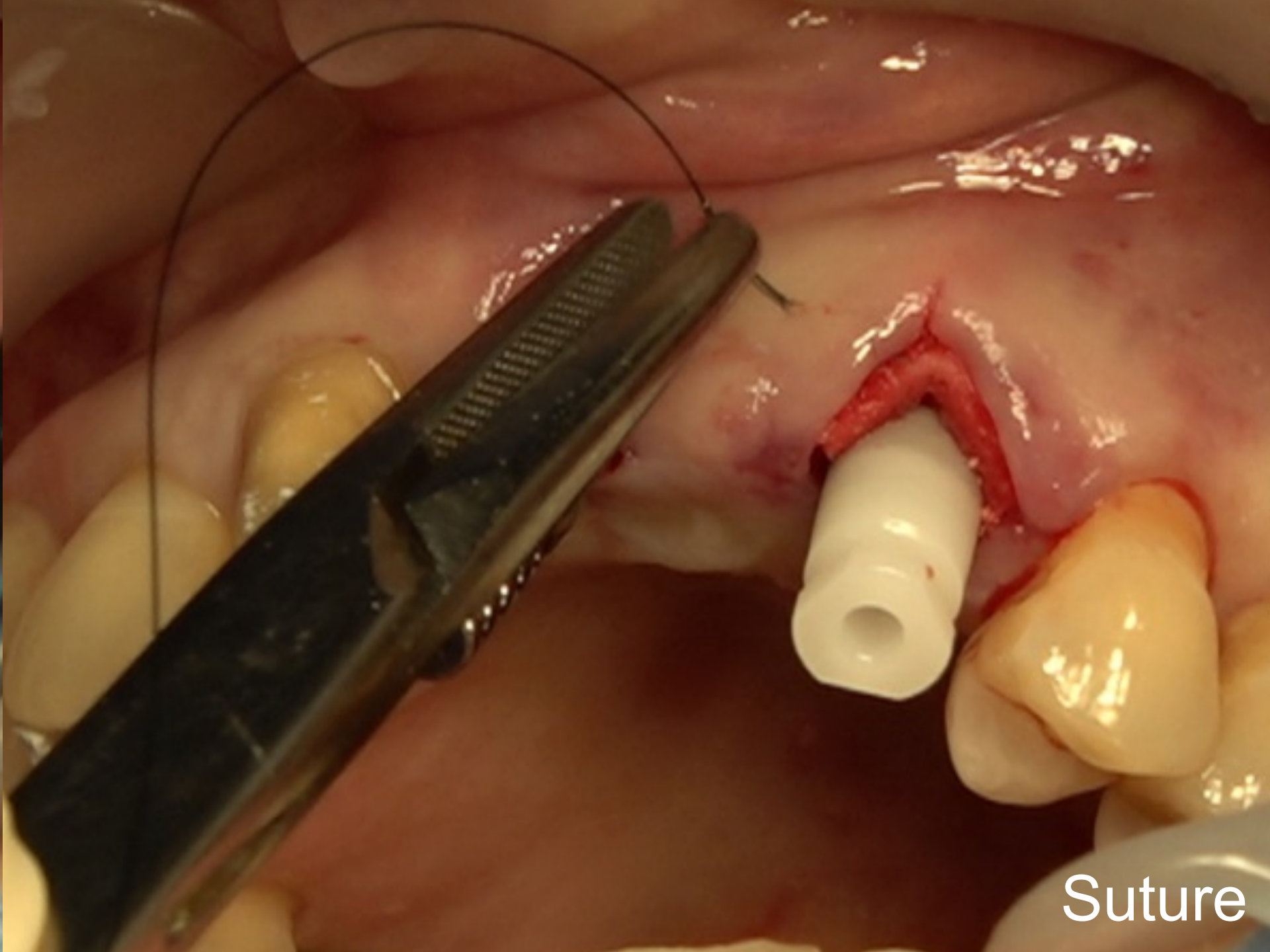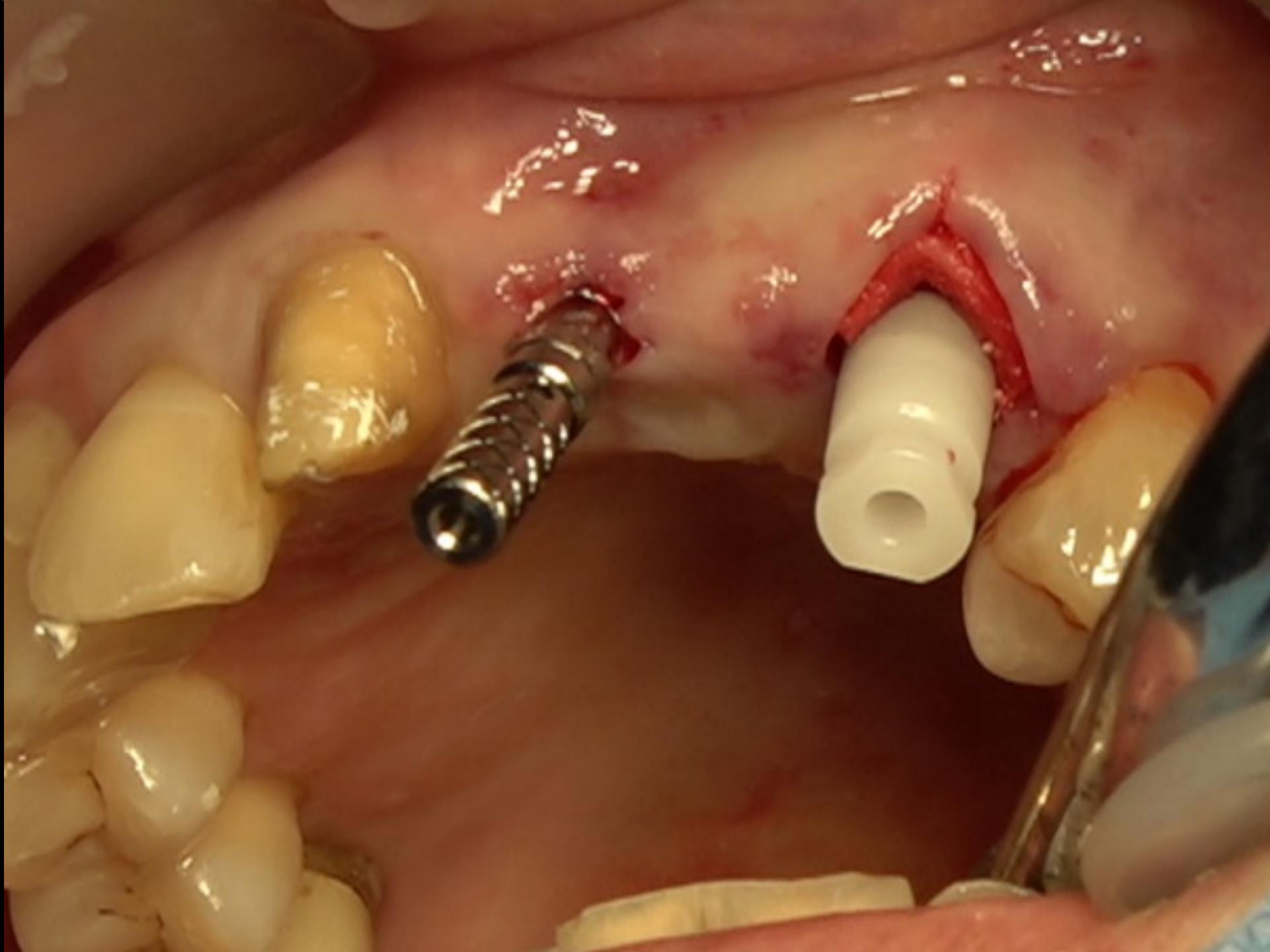

Suture



Post-op (2023-09-06)

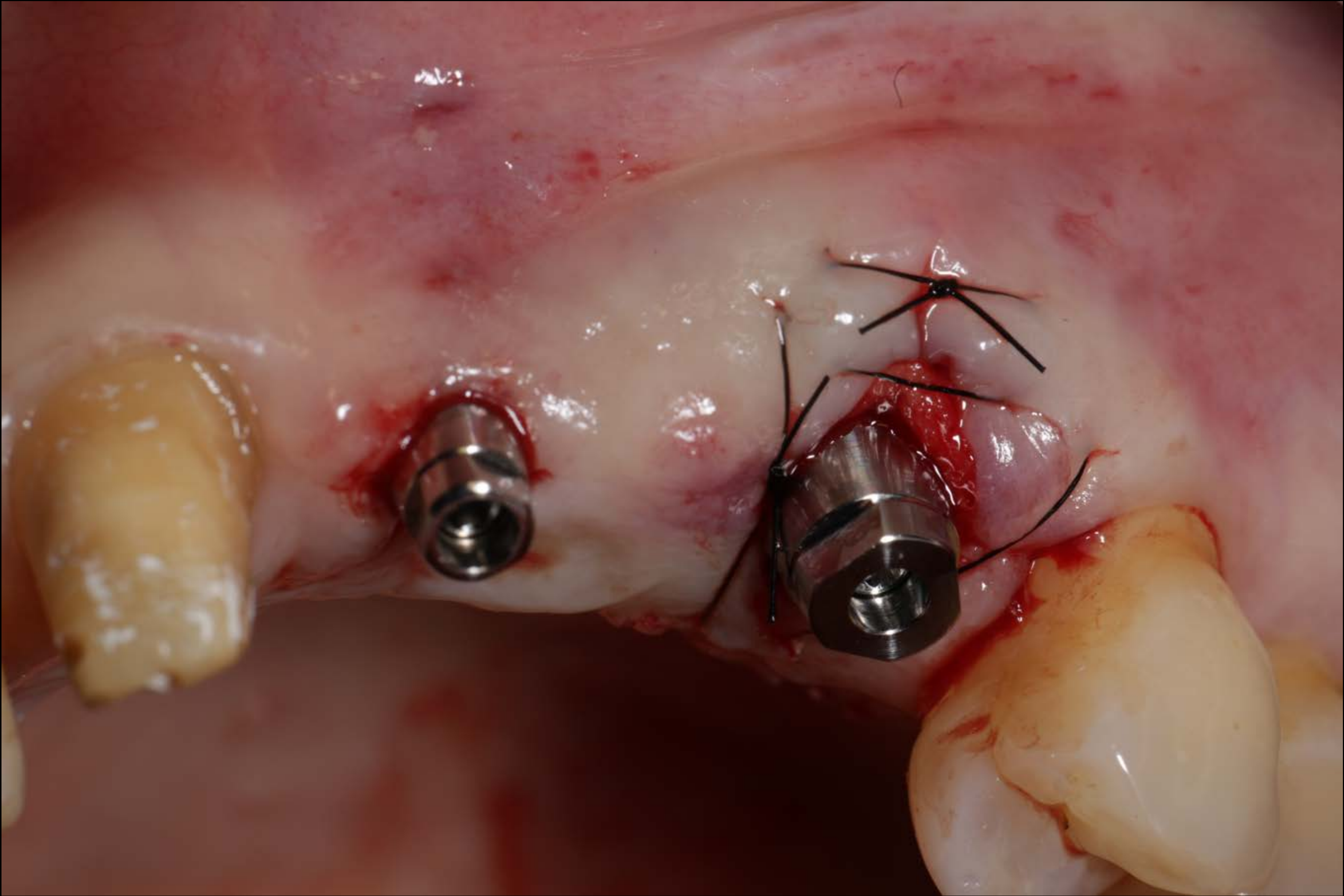







Intra-oral scanning





CAD (Temporary)

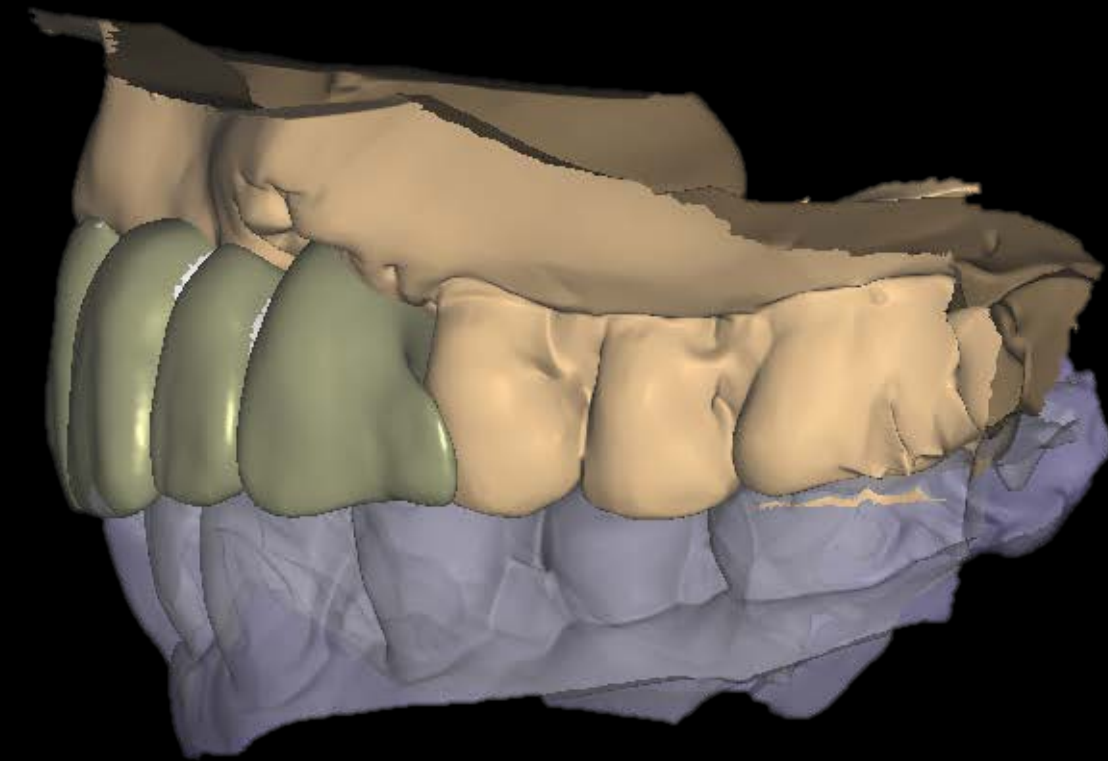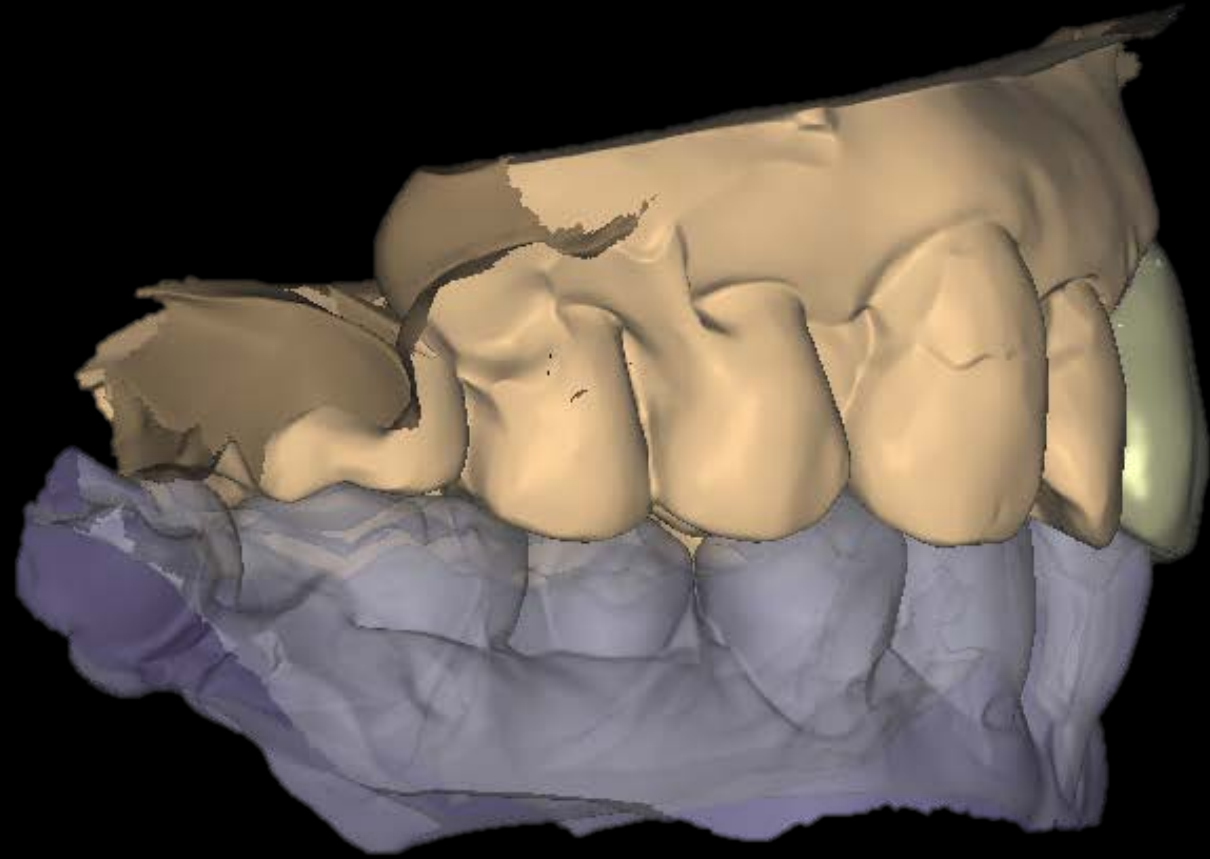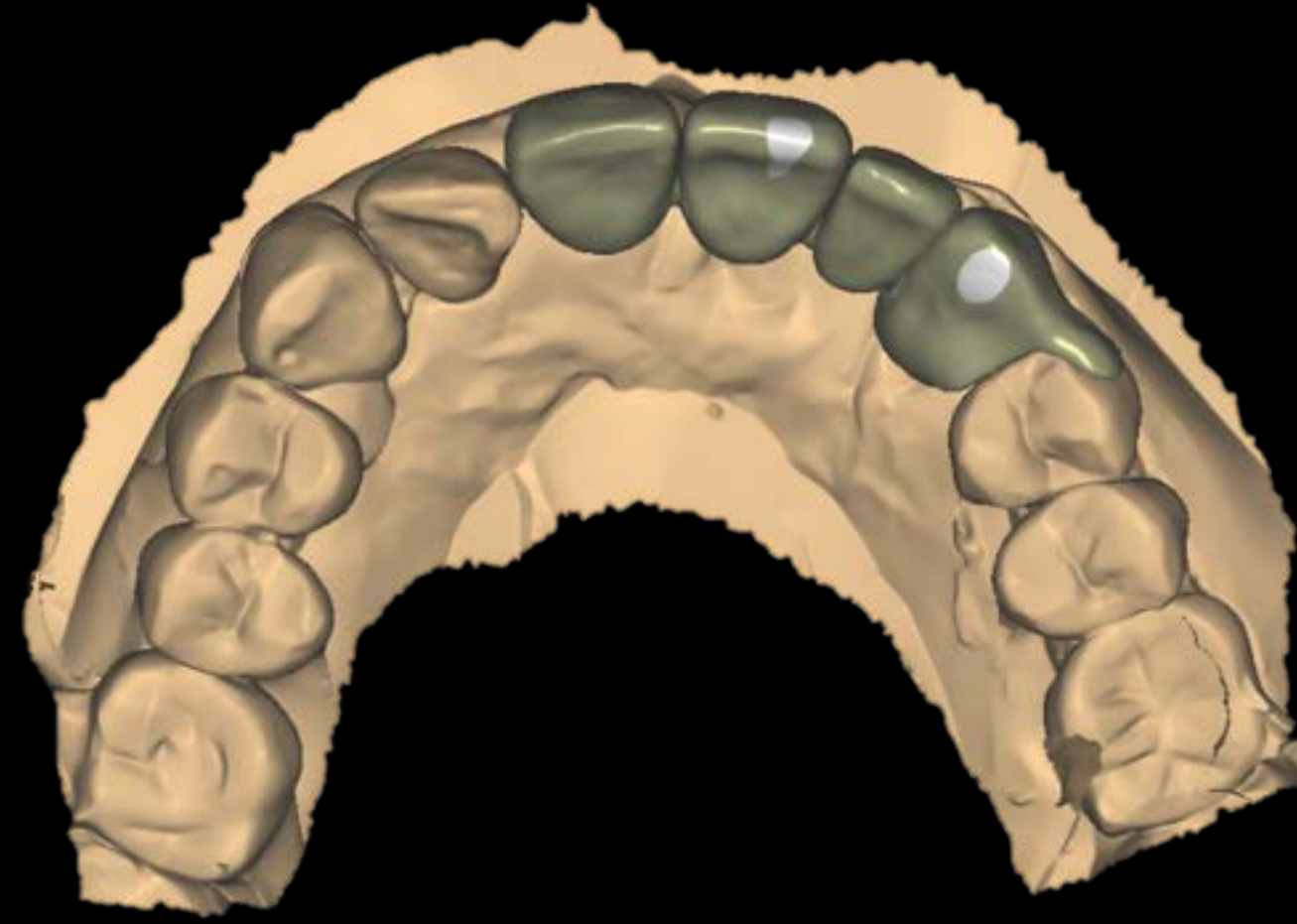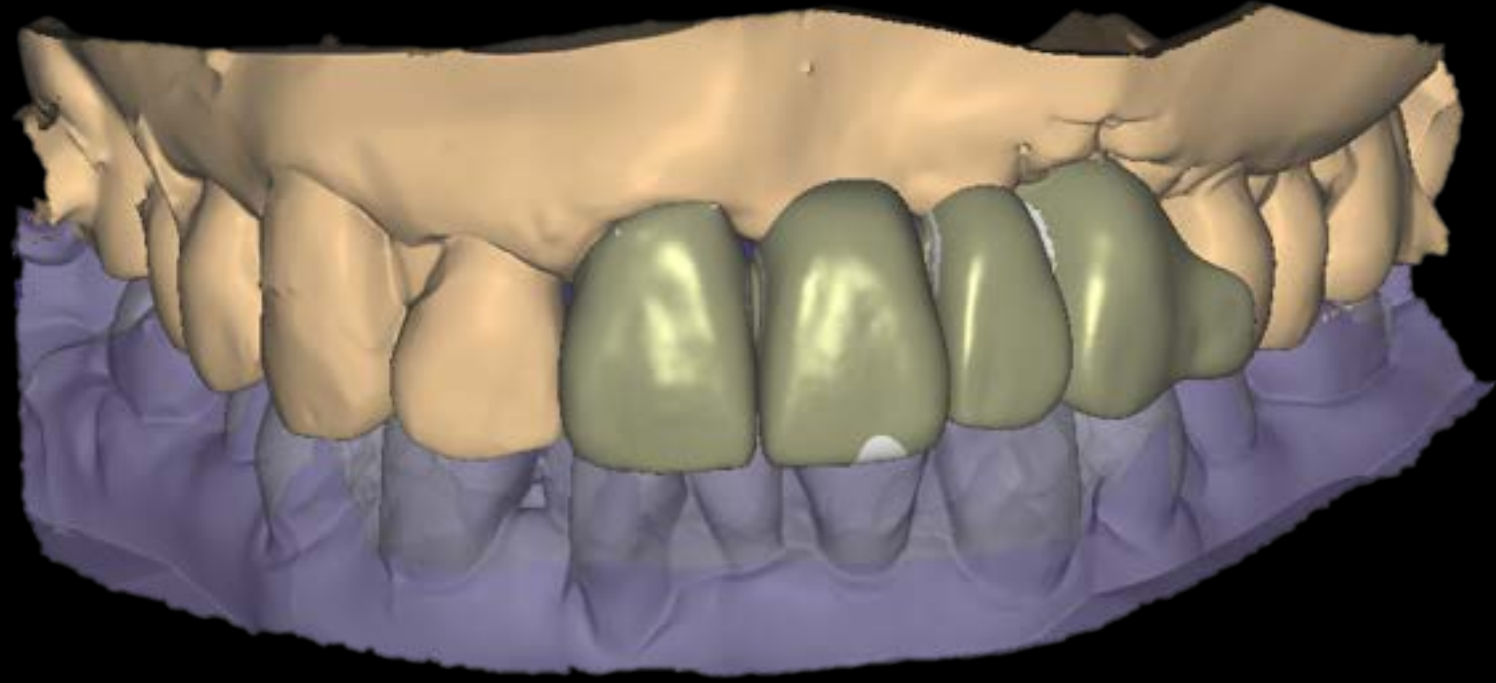



Provisional restoration (2023-11-02)







- **OSTEON™ Xeno**
- **Collagen Membrane**



- **OSTEON™ Xeno with PDRN injection**
- **Collagen Membrane**





- PDRN;  
polydeoxyribonucleotide sodium










- **Bone Patch**



- **OSTEON™ Xeno Collagen**
- **Collagen Graft 2**





Impression / Trimming / Scan spray

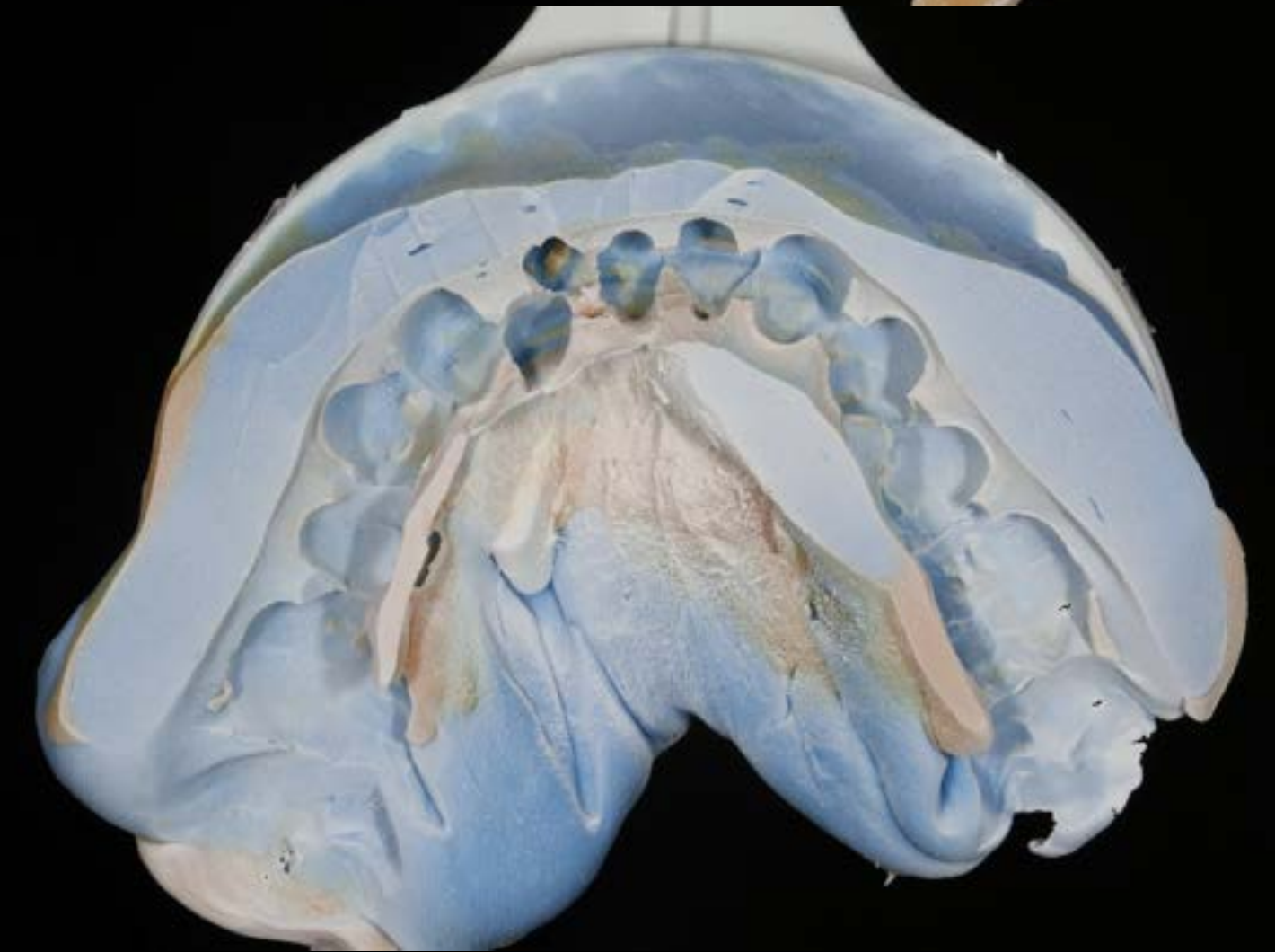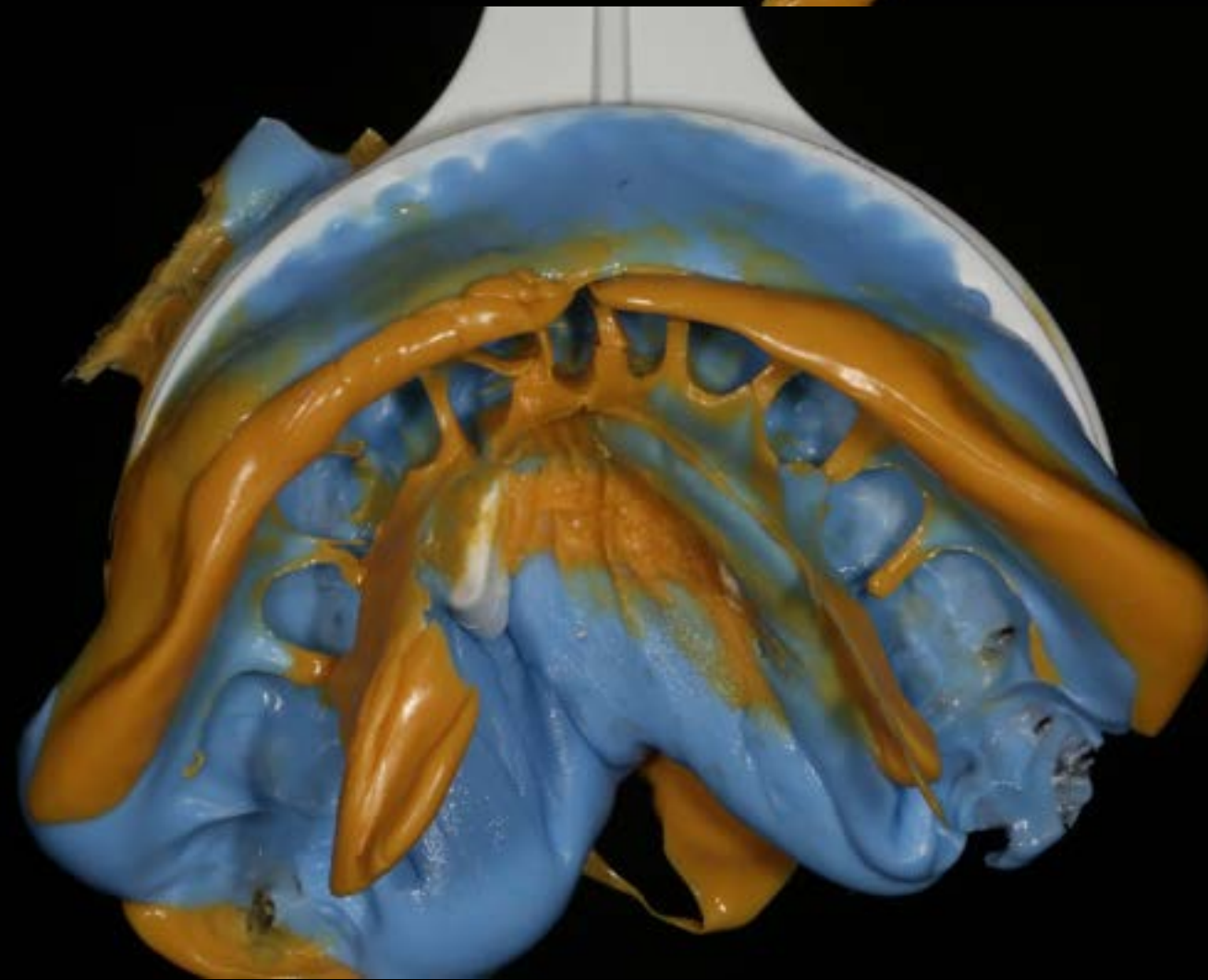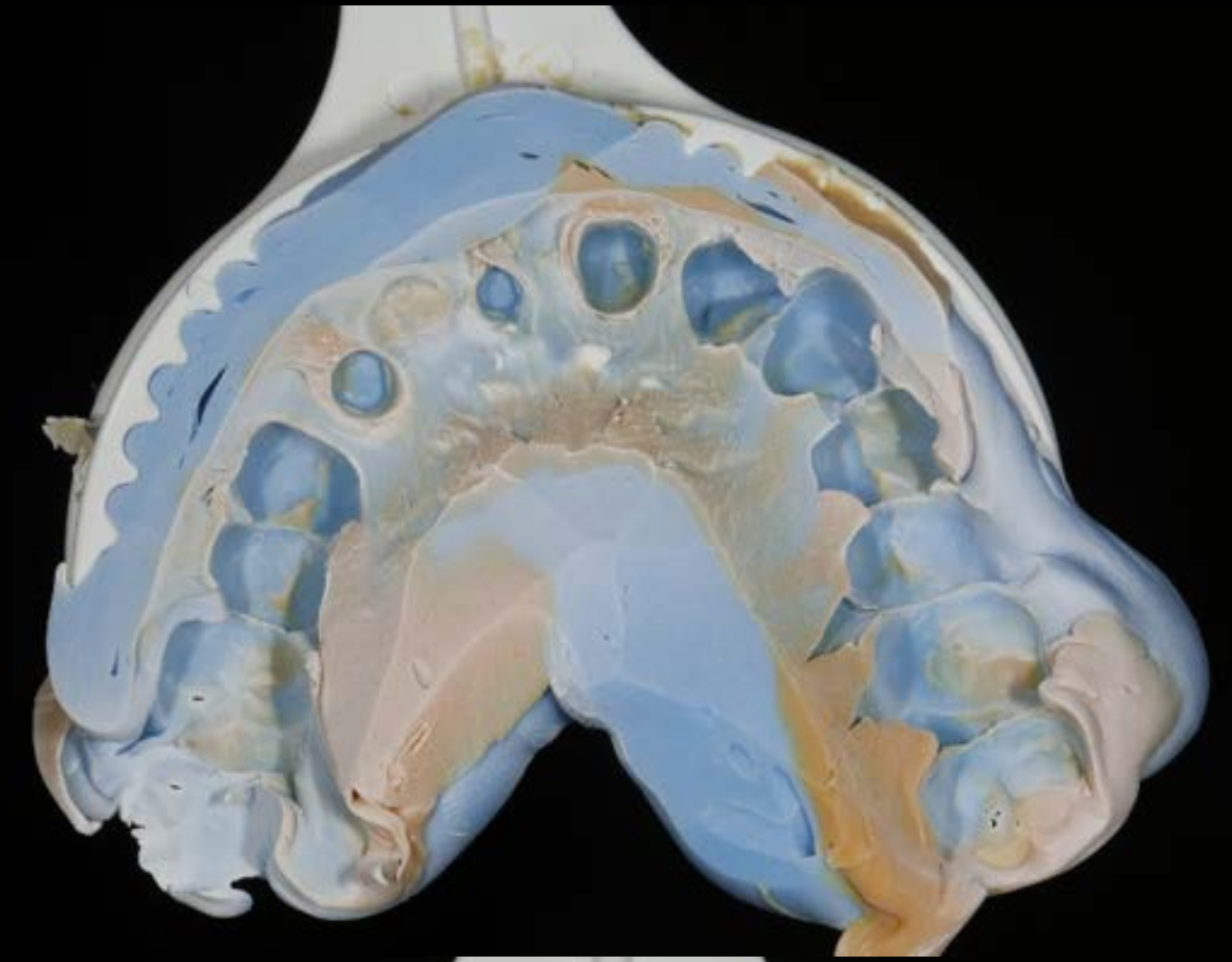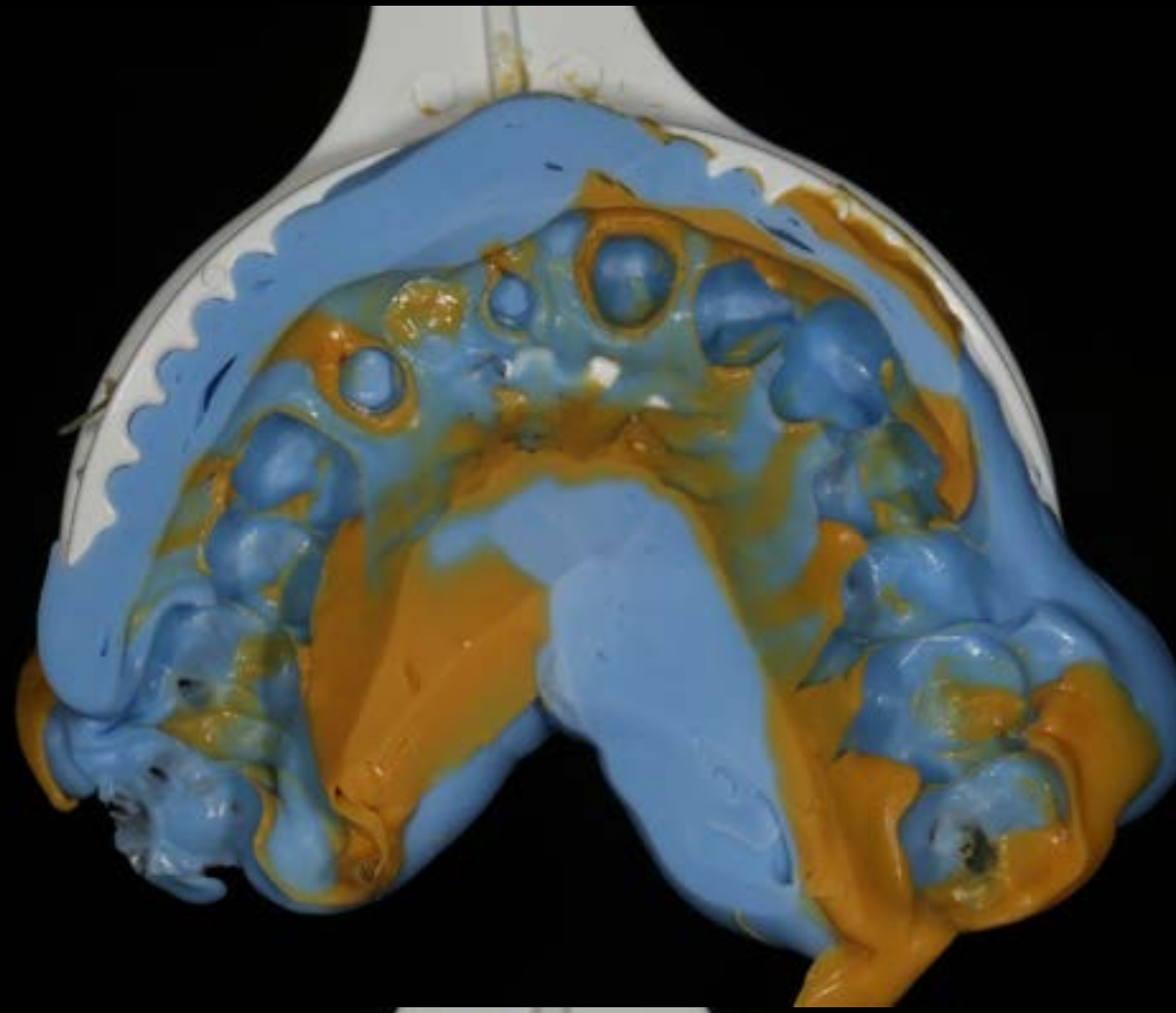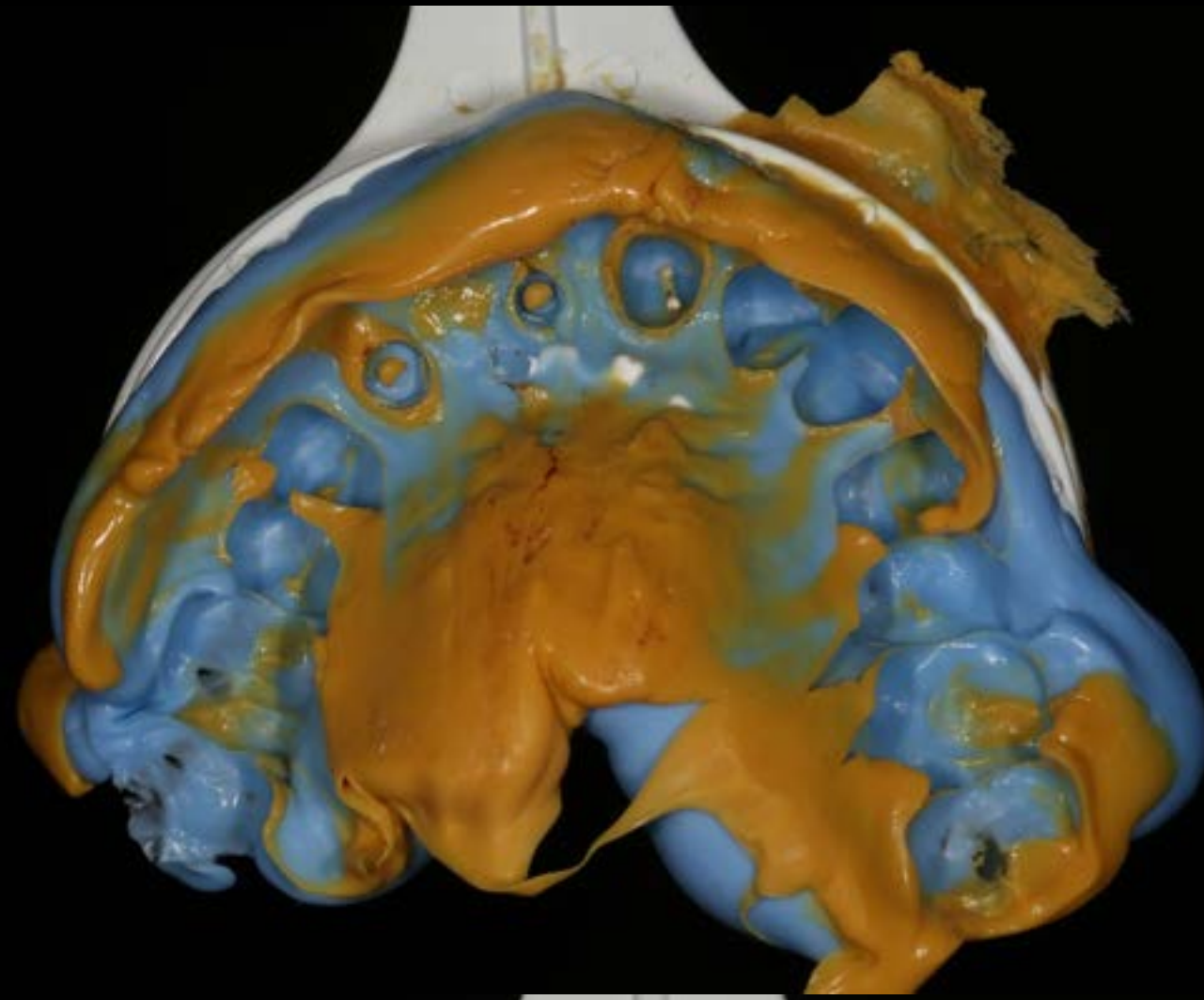




Working Model

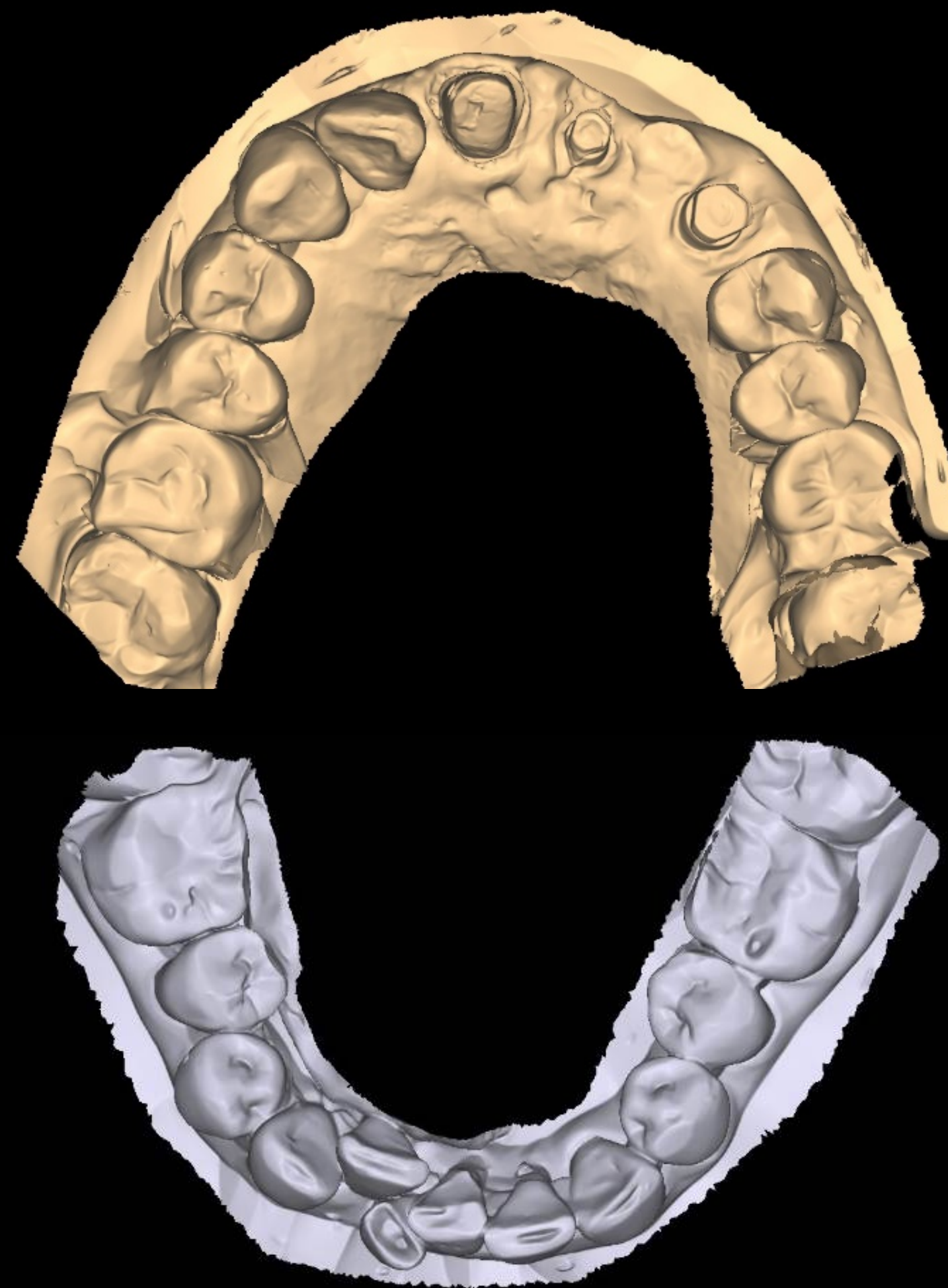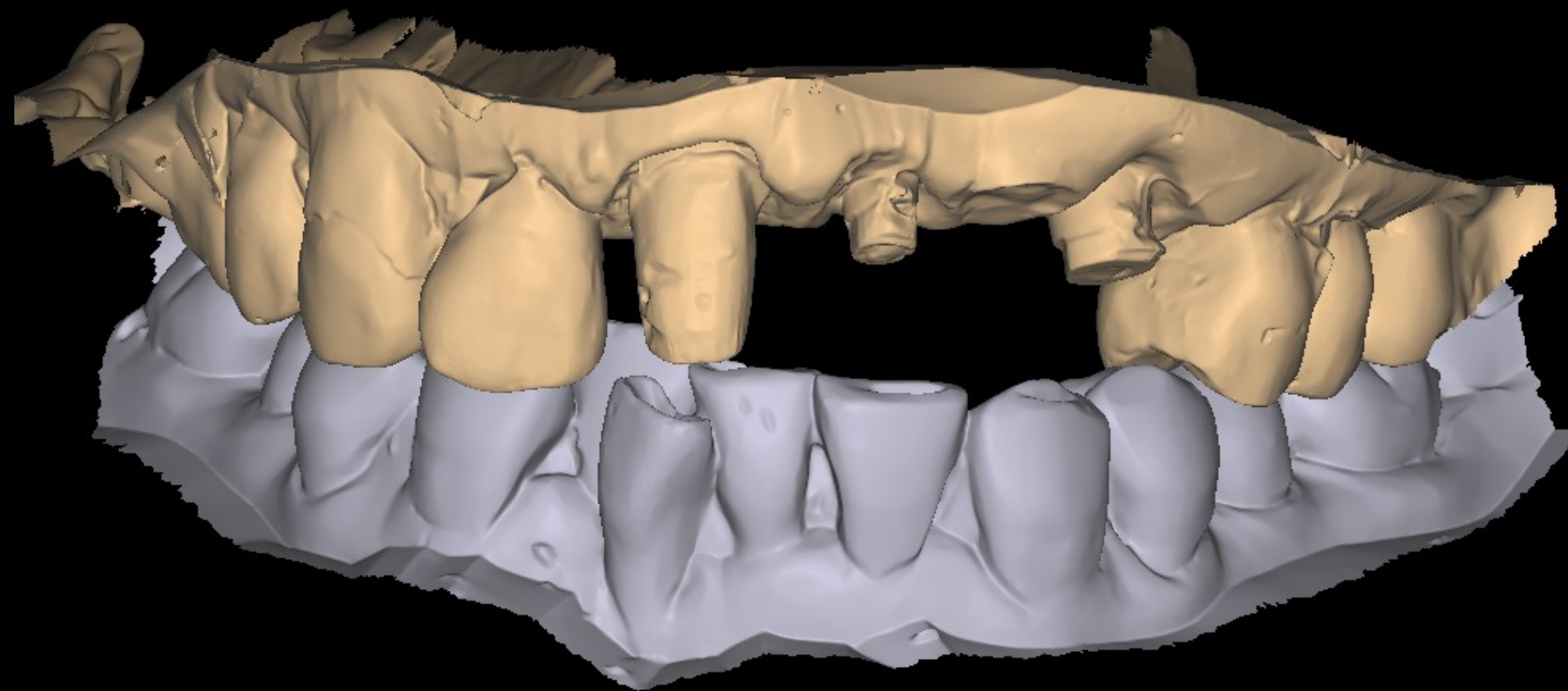



# CT STL export (rainbow 3D Viewer)

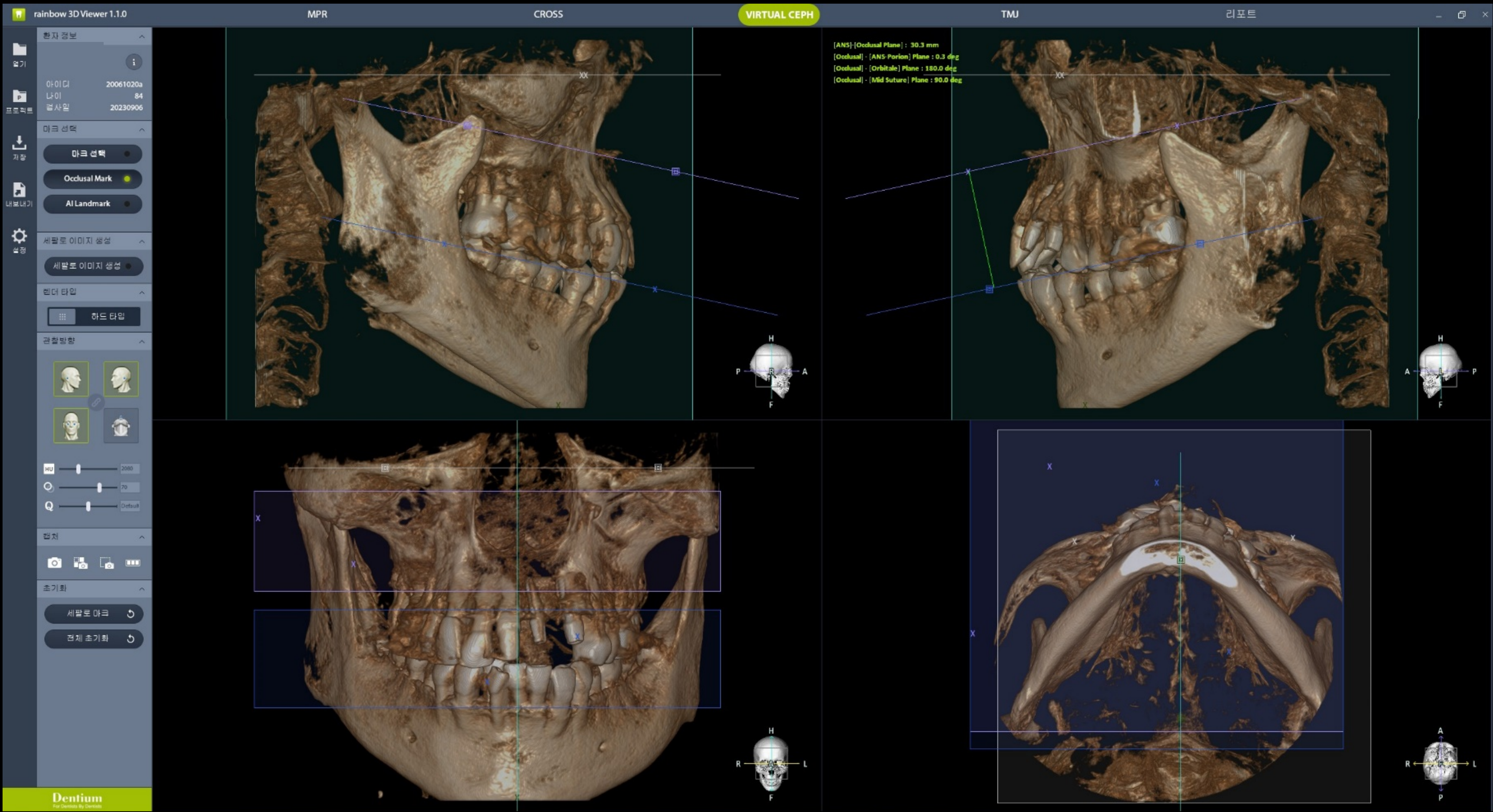



# bright CT 진단에서 치료로

Super Low dose & Large FOV (17.5X15)

100%



20% Low Dose




10% Low Dose



bright Metal Artifact Reduction : Clear Image



OFF

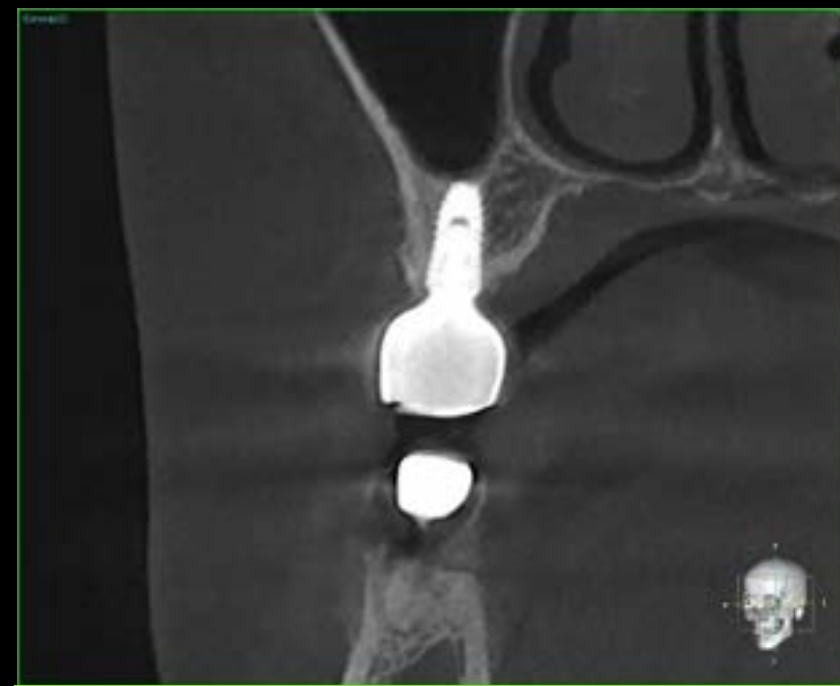

ON

Large FOV



17.5X15 (2Tile Stitching)







New Function





**AI** Occlusal Plane

Auto Set occlusal plane from anatomical landmark on CT

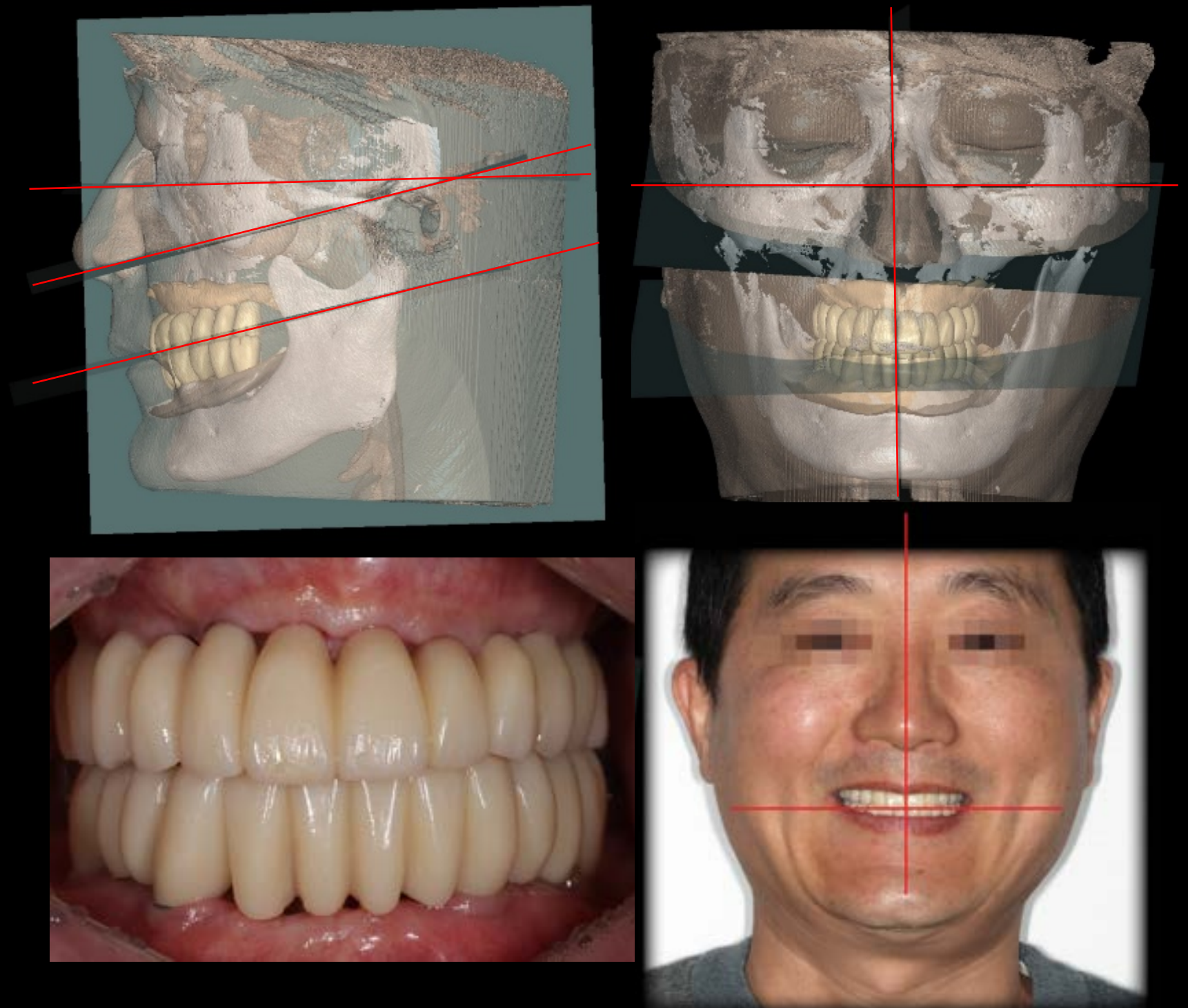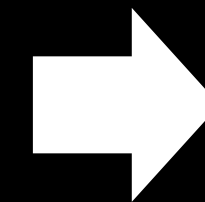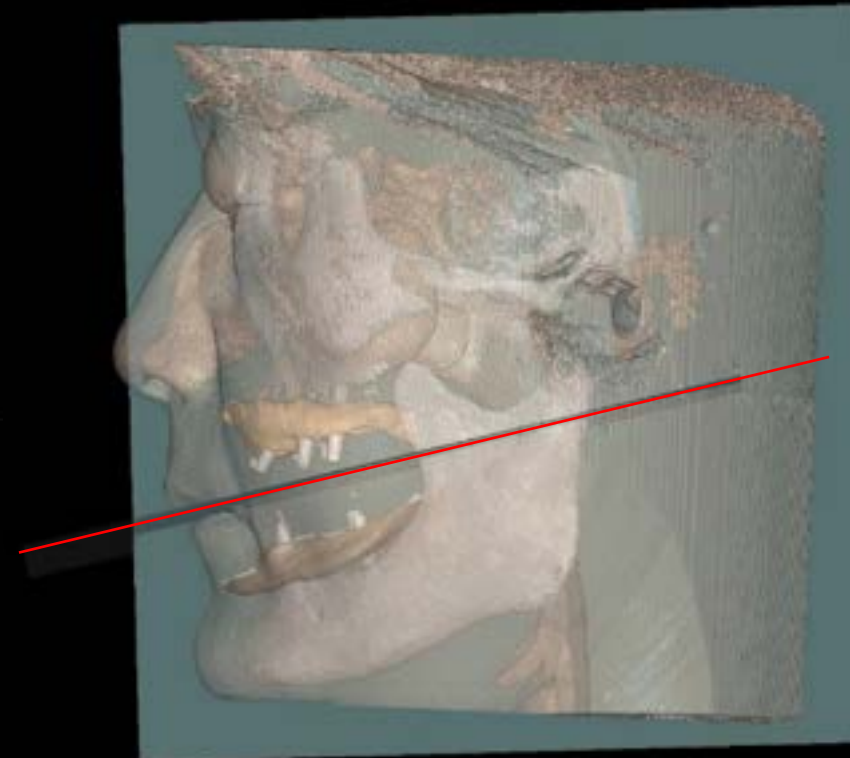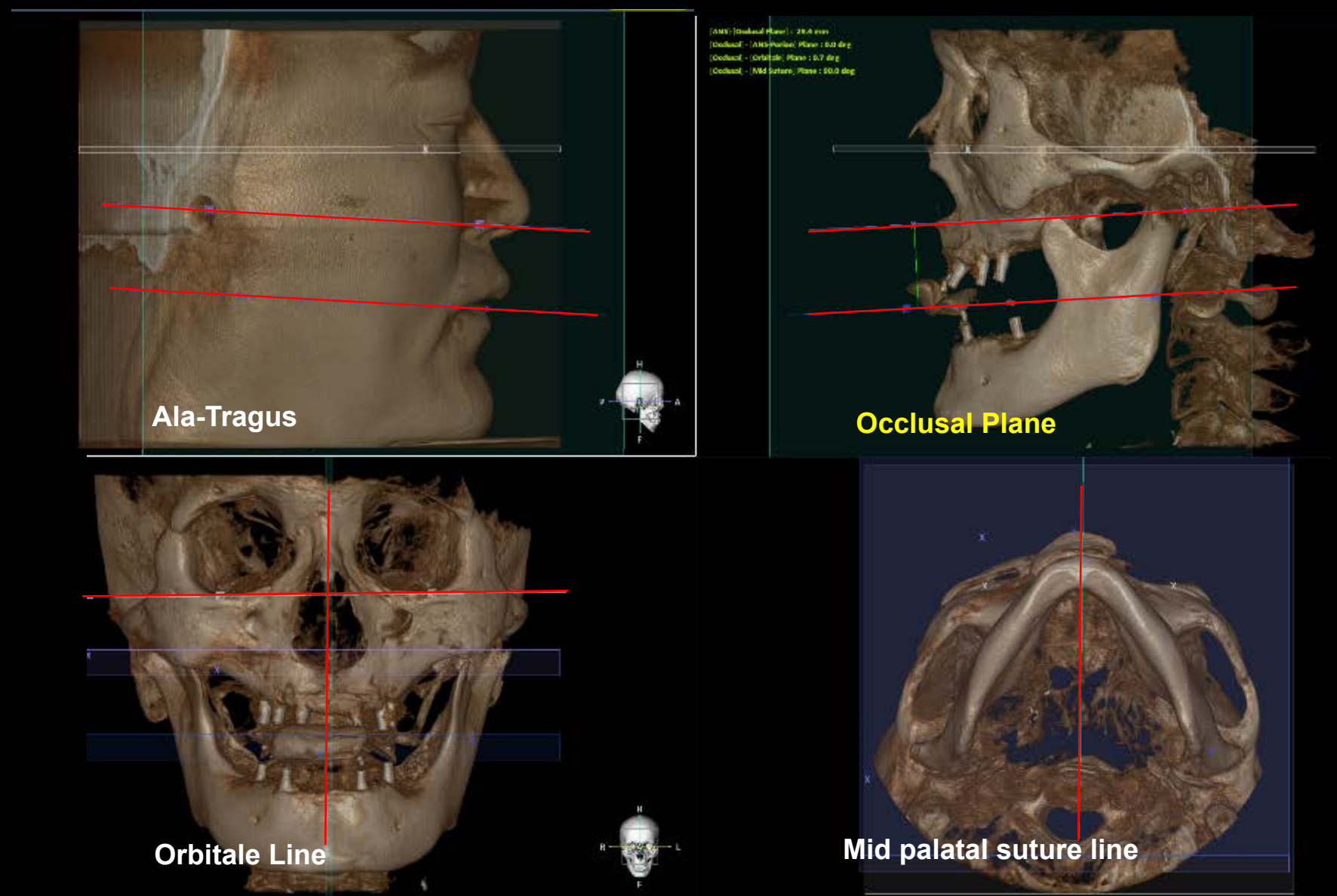






## CT STL & working model stitching

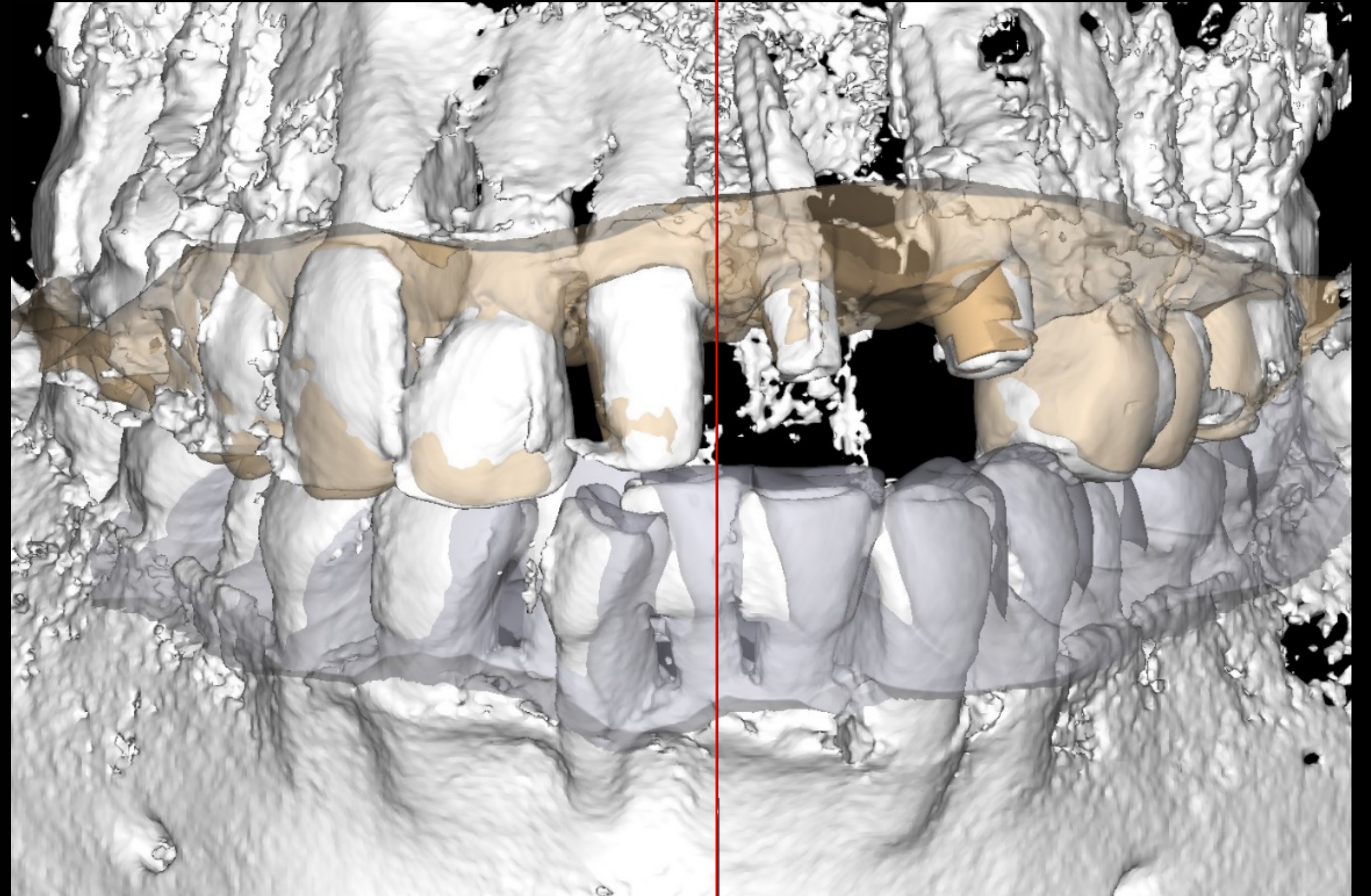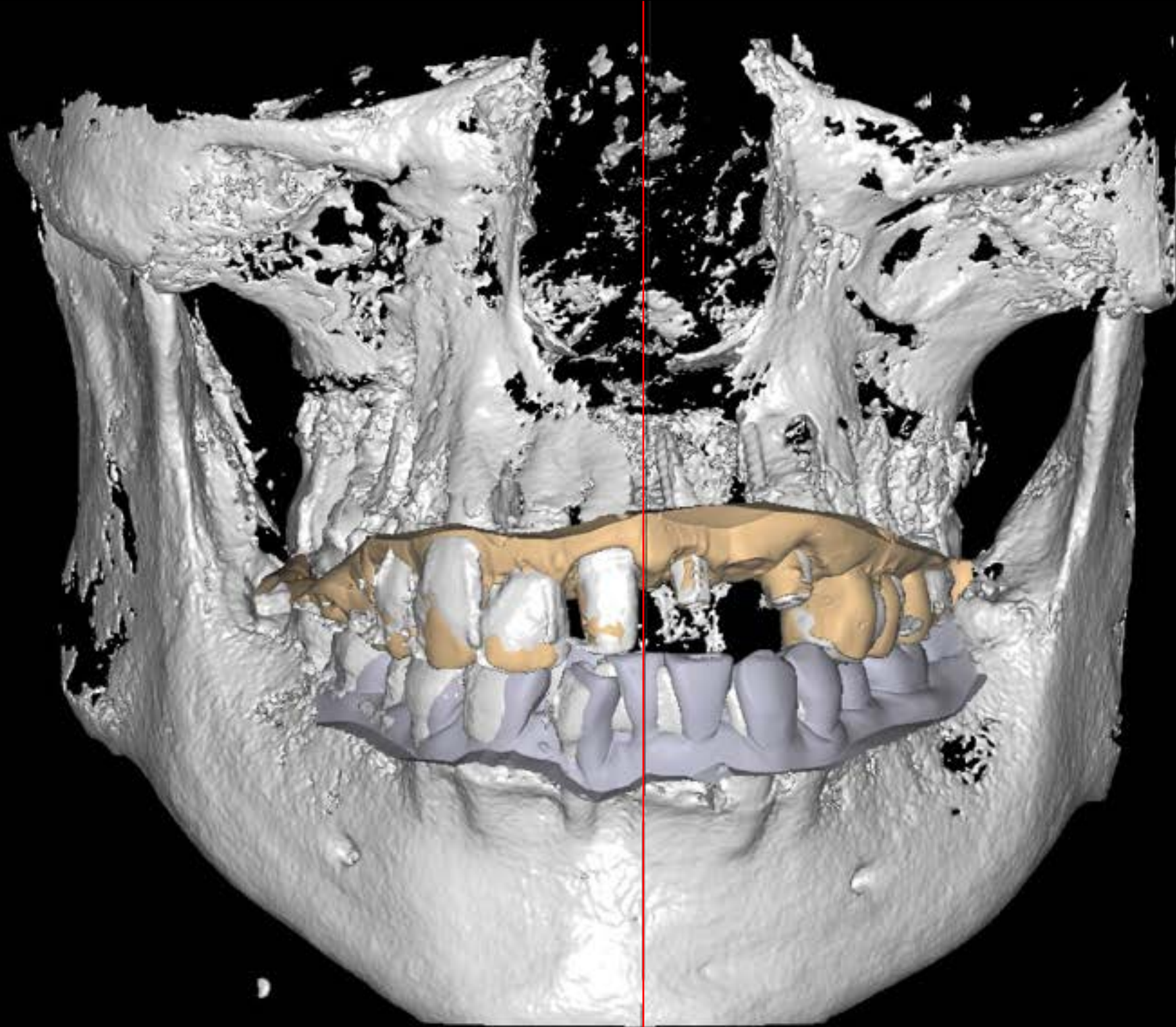




CAD (Final prosthesis)

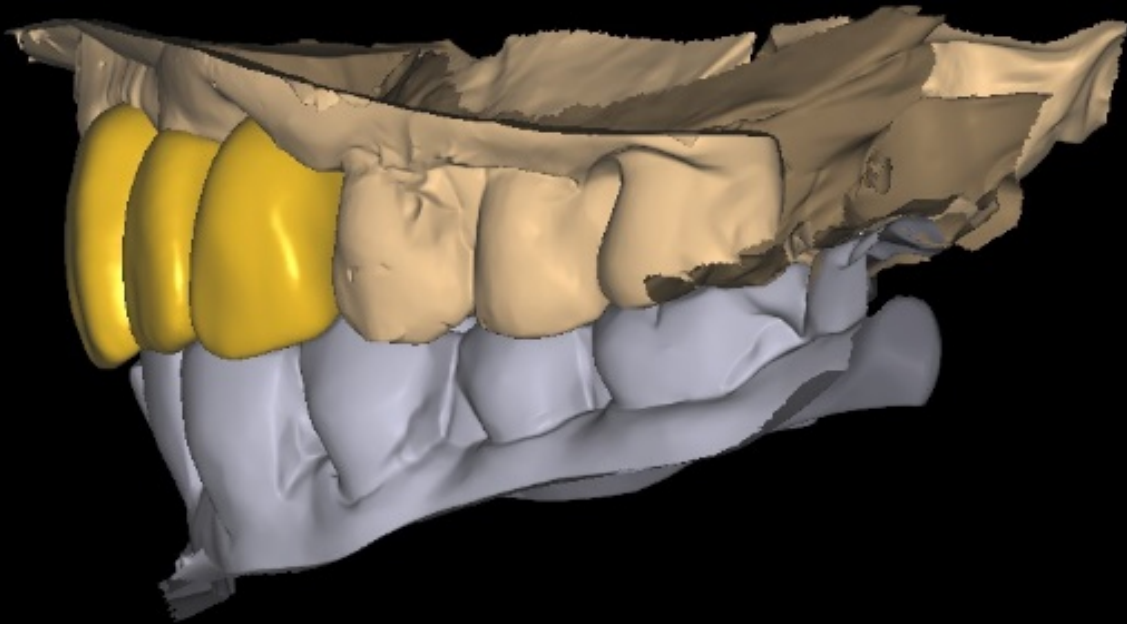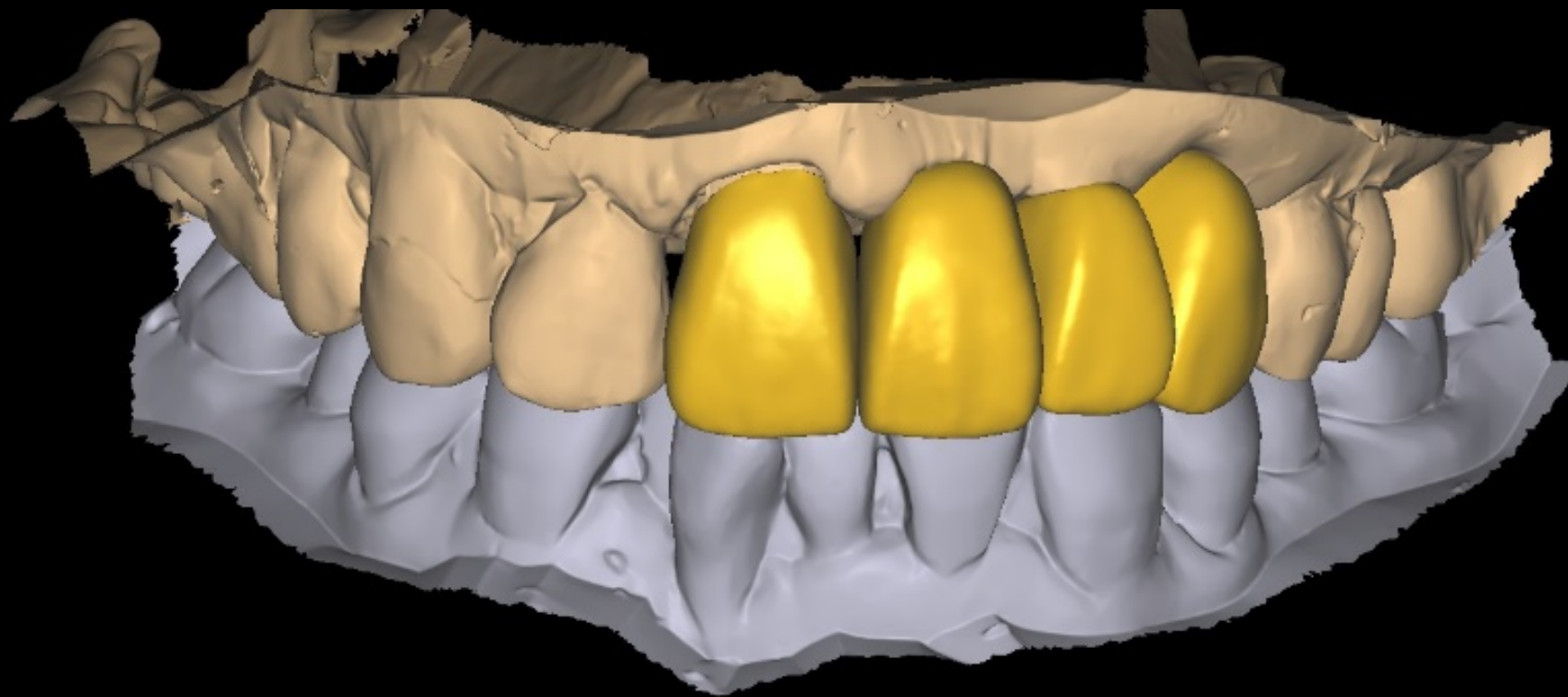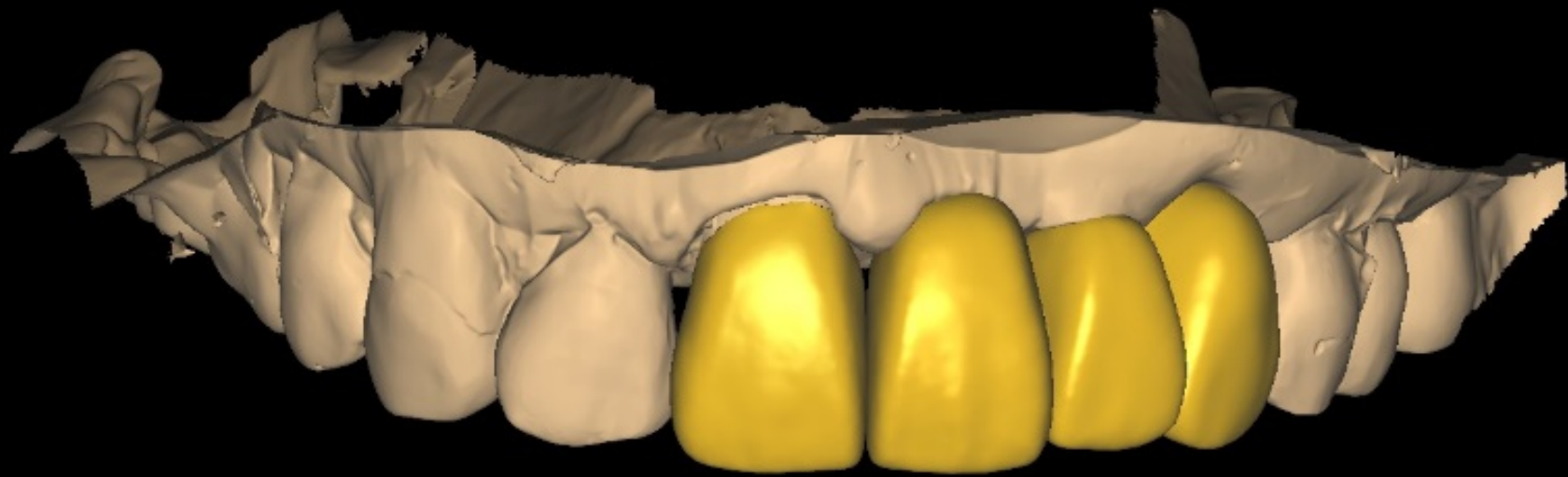





CT – Cad design stitching

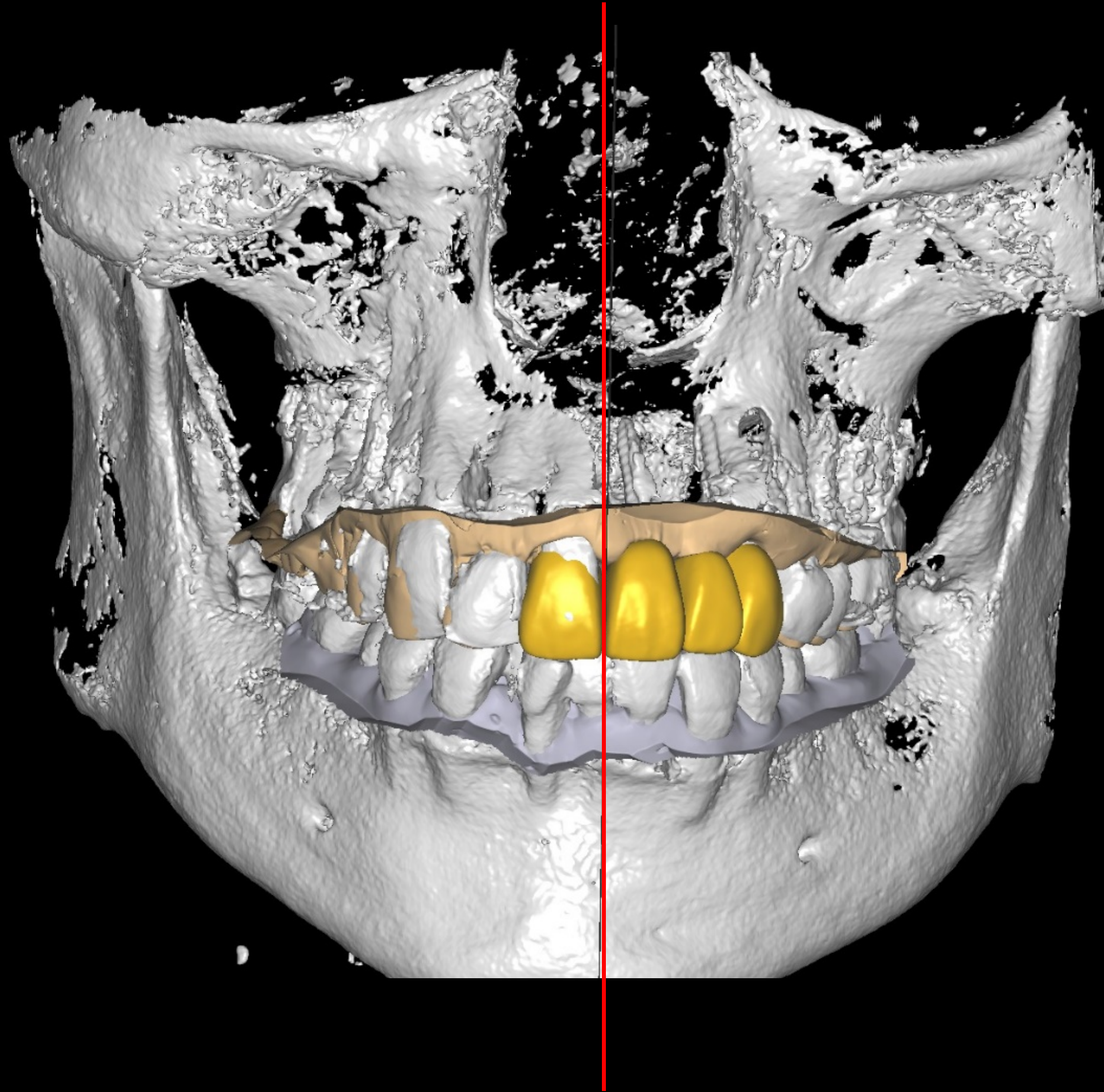



CAD : Mounting with virtual articulator

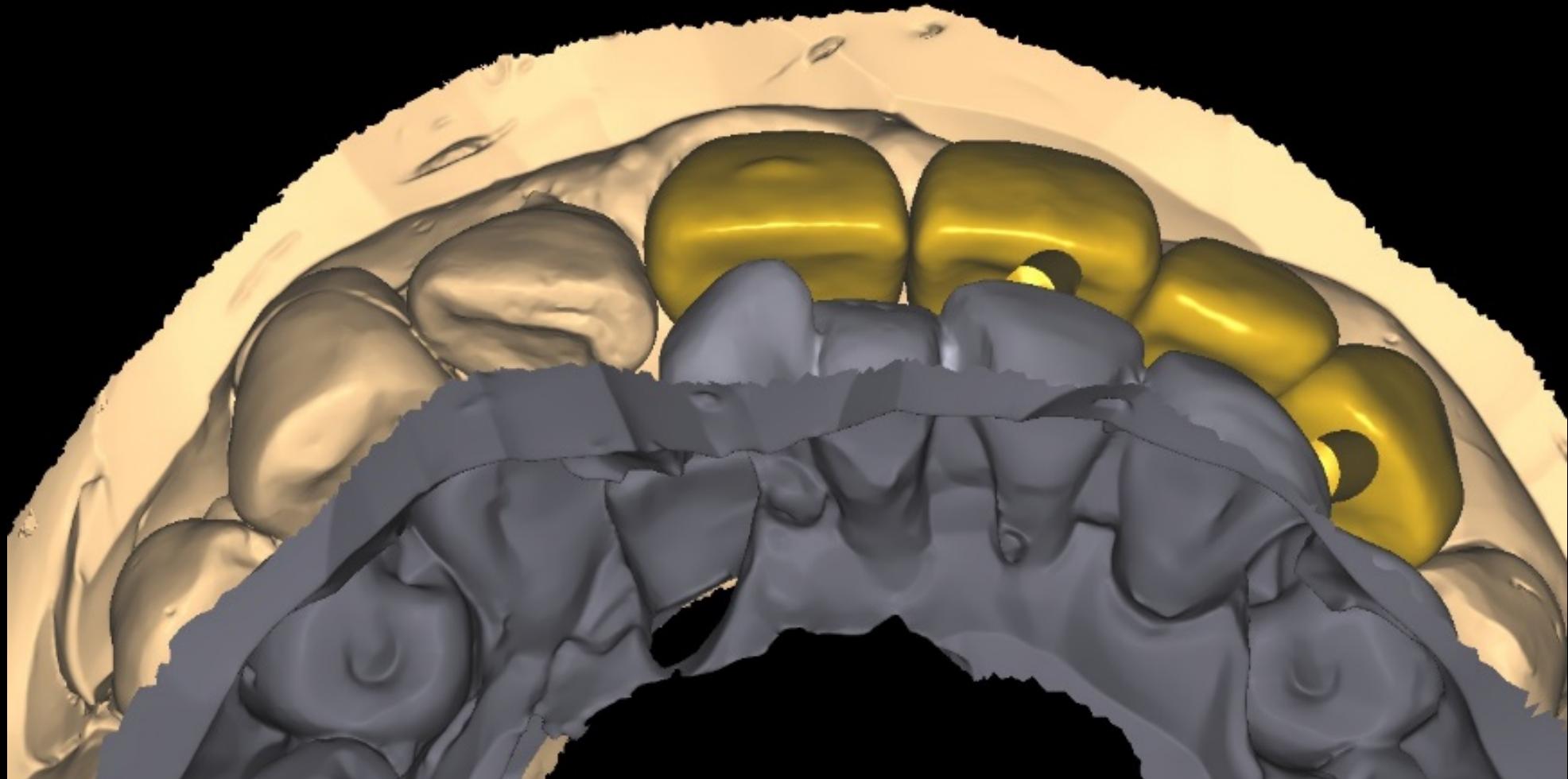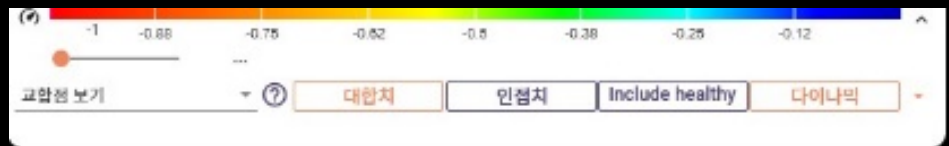




Shade taking(A3)





Hyper dent

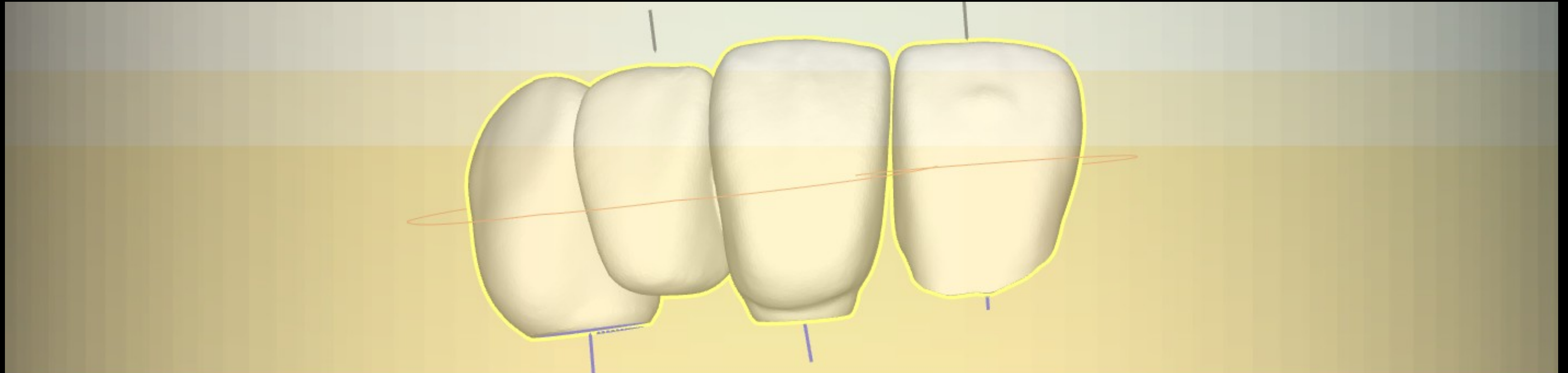

**A3 22T 3Layer Block**



Milling / Coloring / Sintering / Glazing



Inner – White Opauqe





Final prosthesis (2023-11-02)





# bright 3-Layer




**49% Transparency**  
**층간 자연스러운 Value변화**

**No Coloring & Staining**

직관적인 색 결정

절단층과 중간층 고정

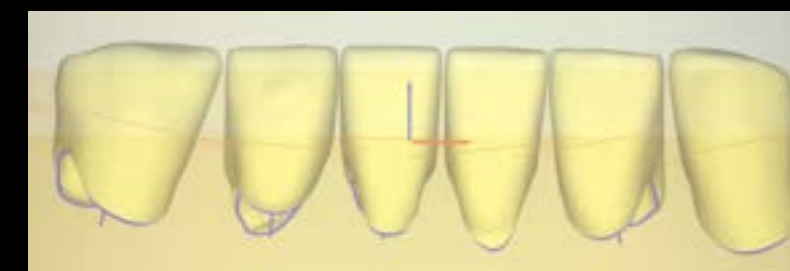




2 types of Zirconia

3 layer of different value

2 hour fast sintered mono block



# 3 layer Zr block

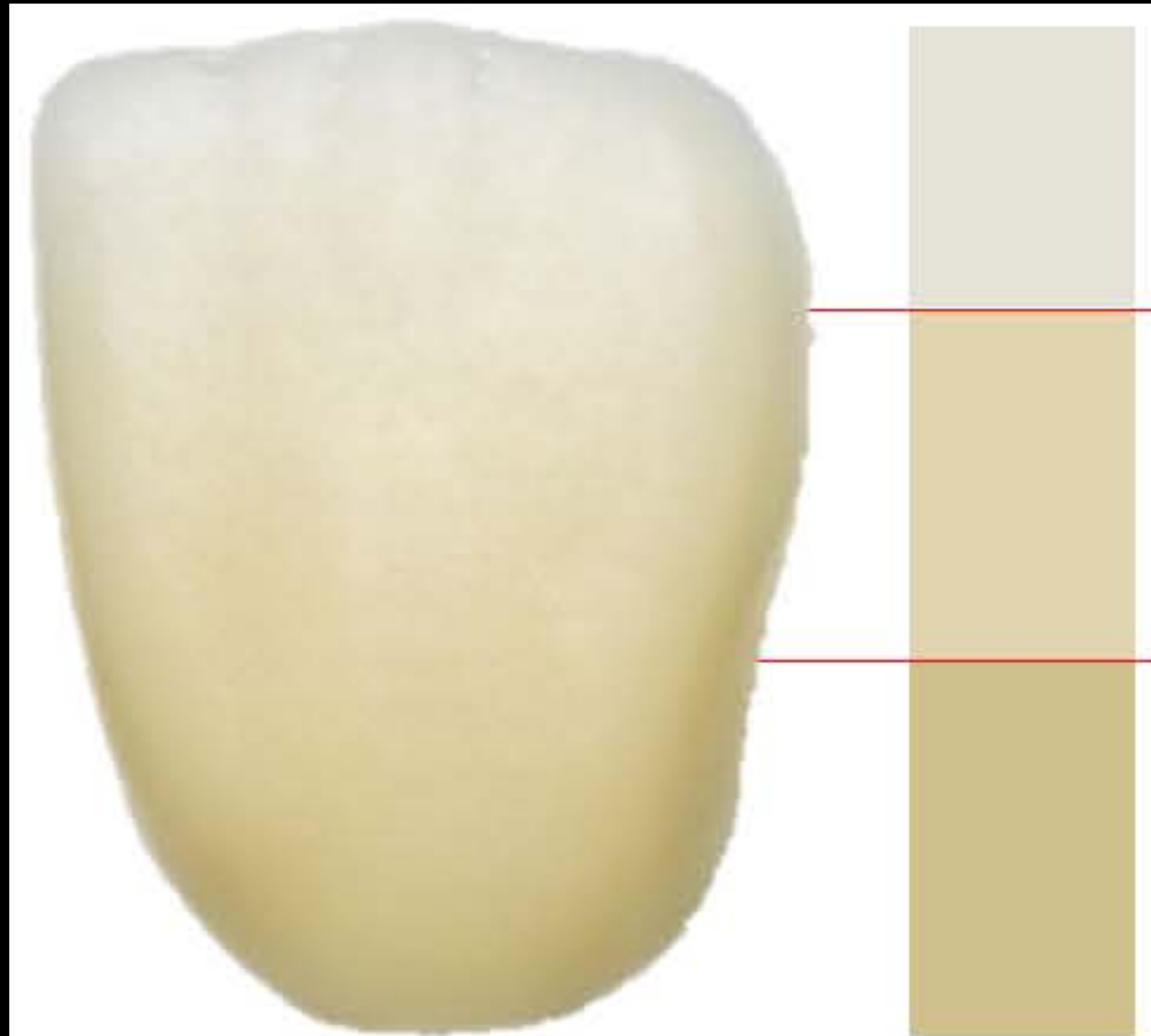

5 Y

Translucent 48  
Value -0.3

5Y + 4Y

Translucent 47  
Value 0

4Y

Trandlucent 46  
Value +0.3



6) Natural teeth prosthesis \_ upper Anterior

## Final prosthesis



내면 – White opaque



bright 3 Layer A3





Oral exam








- Bright Impress **Light**
- Bright Impress **Medium**
- Bright Impress **Heavy**
- Bright Impress **Bite**
- Bright Impress **Putty**
- Bright **Bite Tray**

Dentium



Impression / Trimming / Scan spray



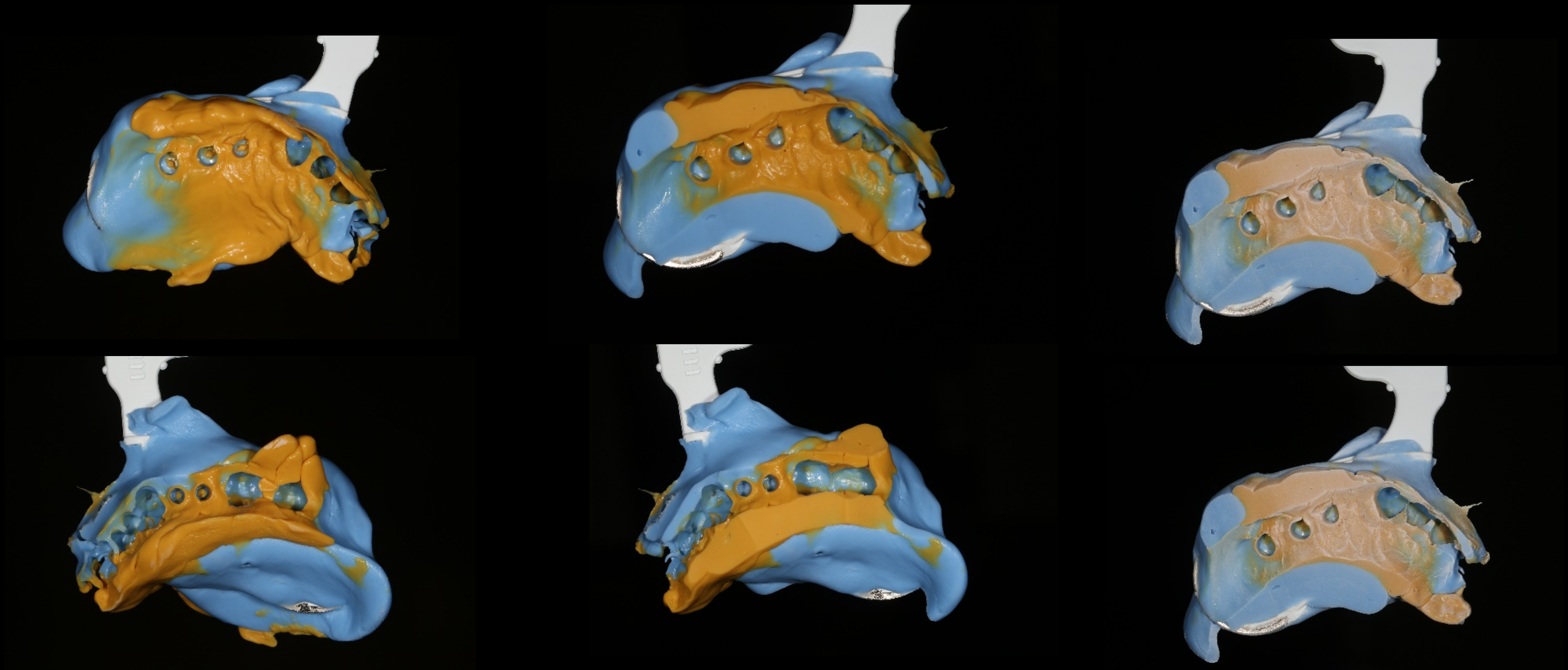




# Impression scan / Working Model






Shining 3D  
AutoScan-DS-EX Pro (H)

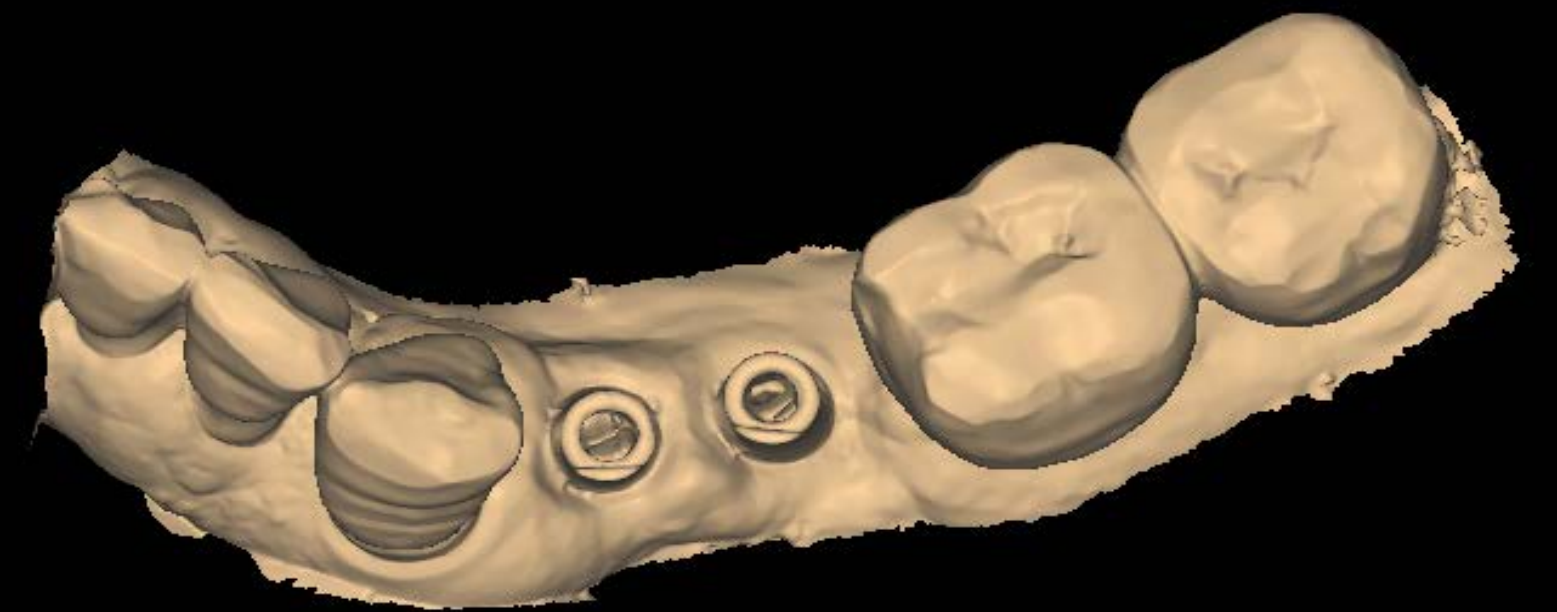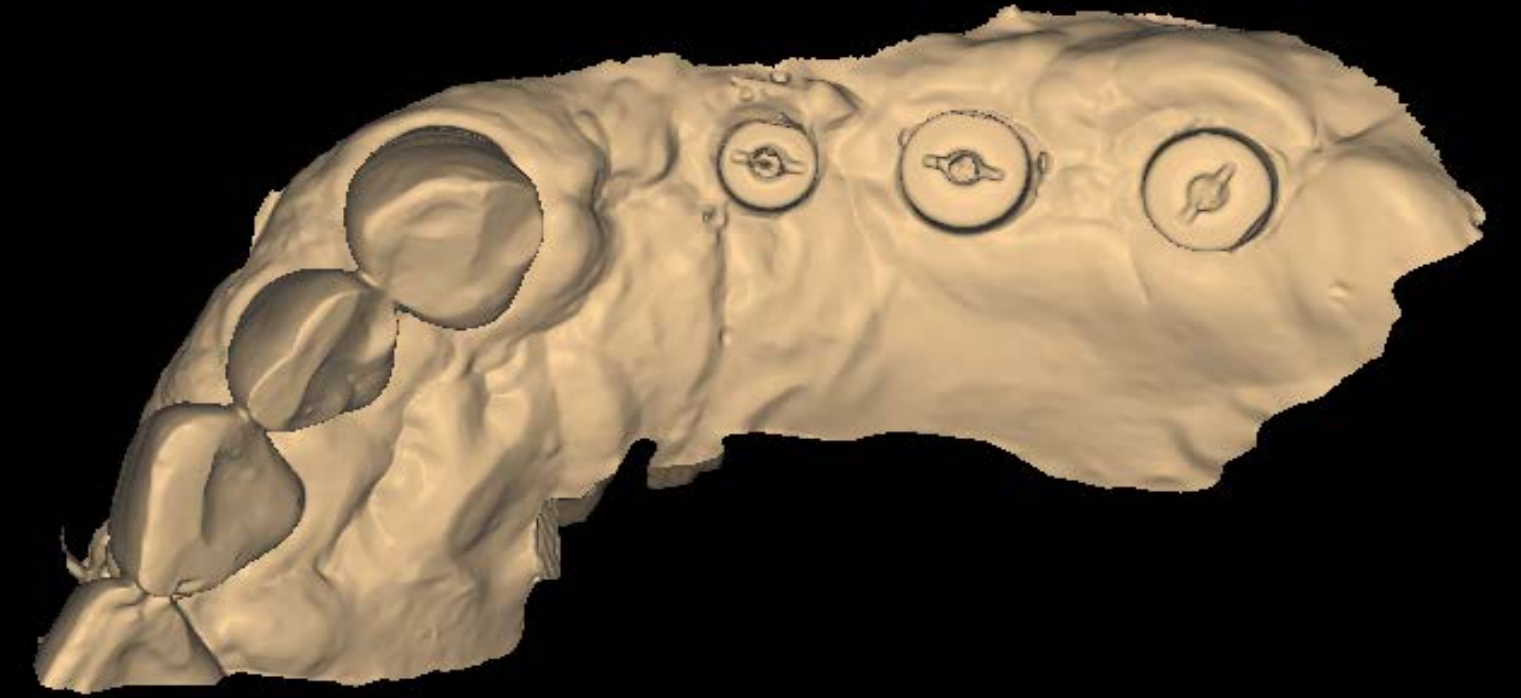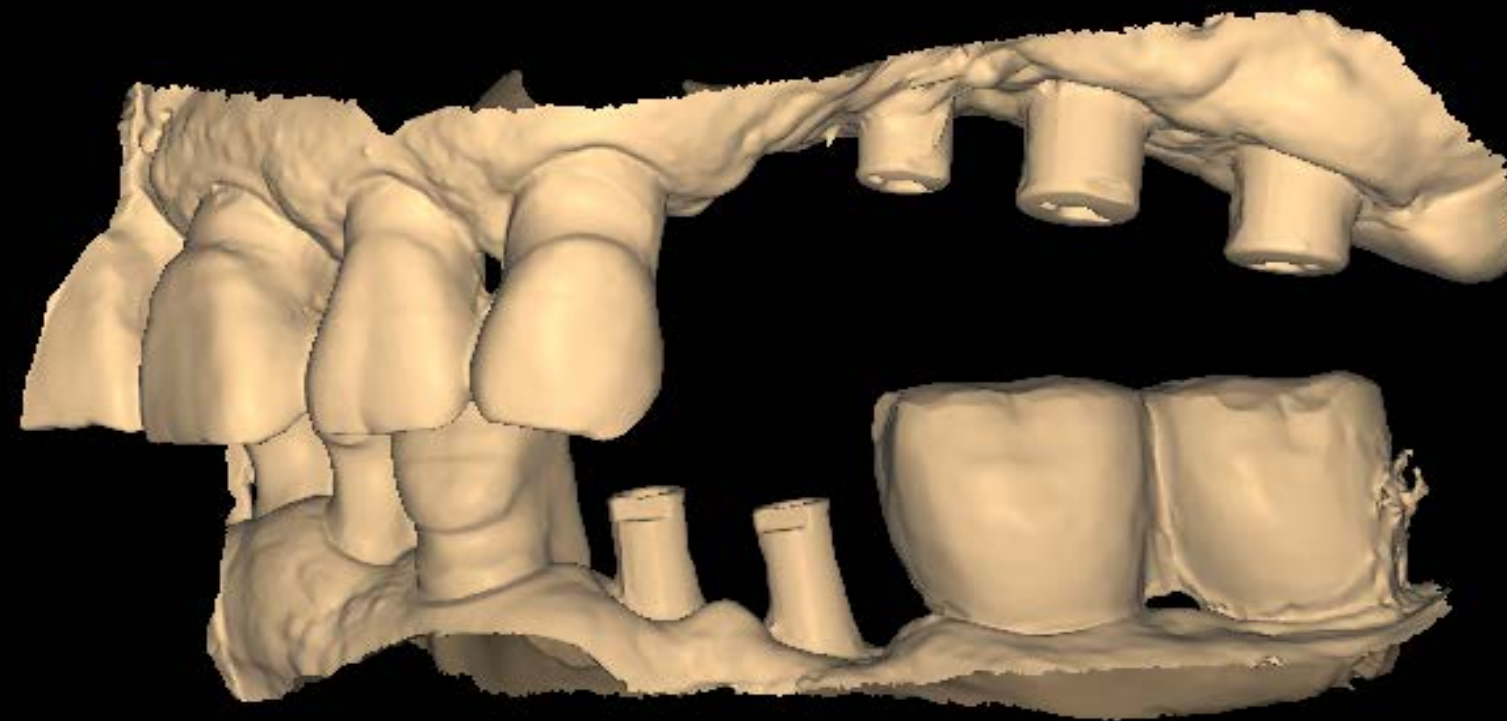



CAD

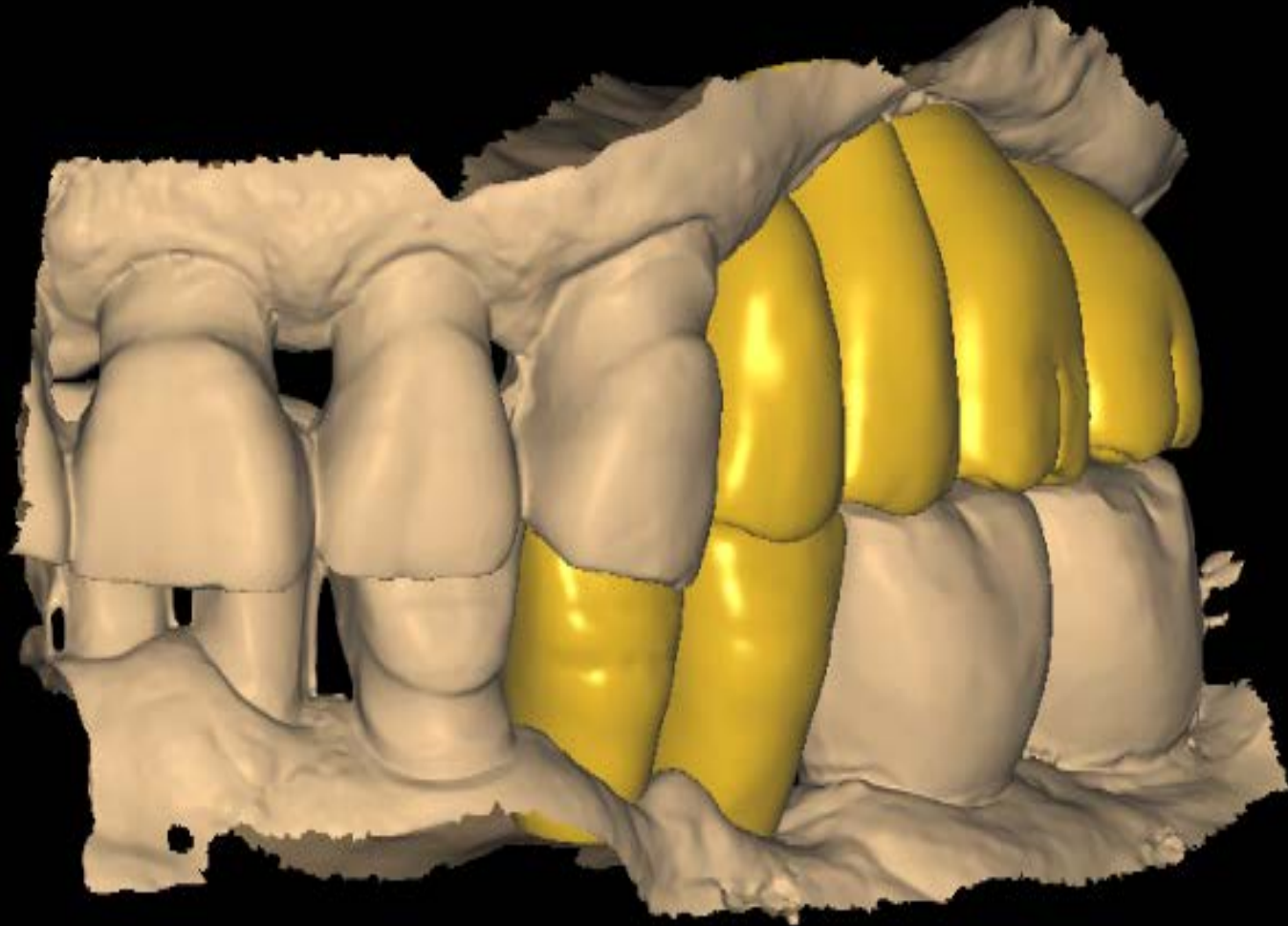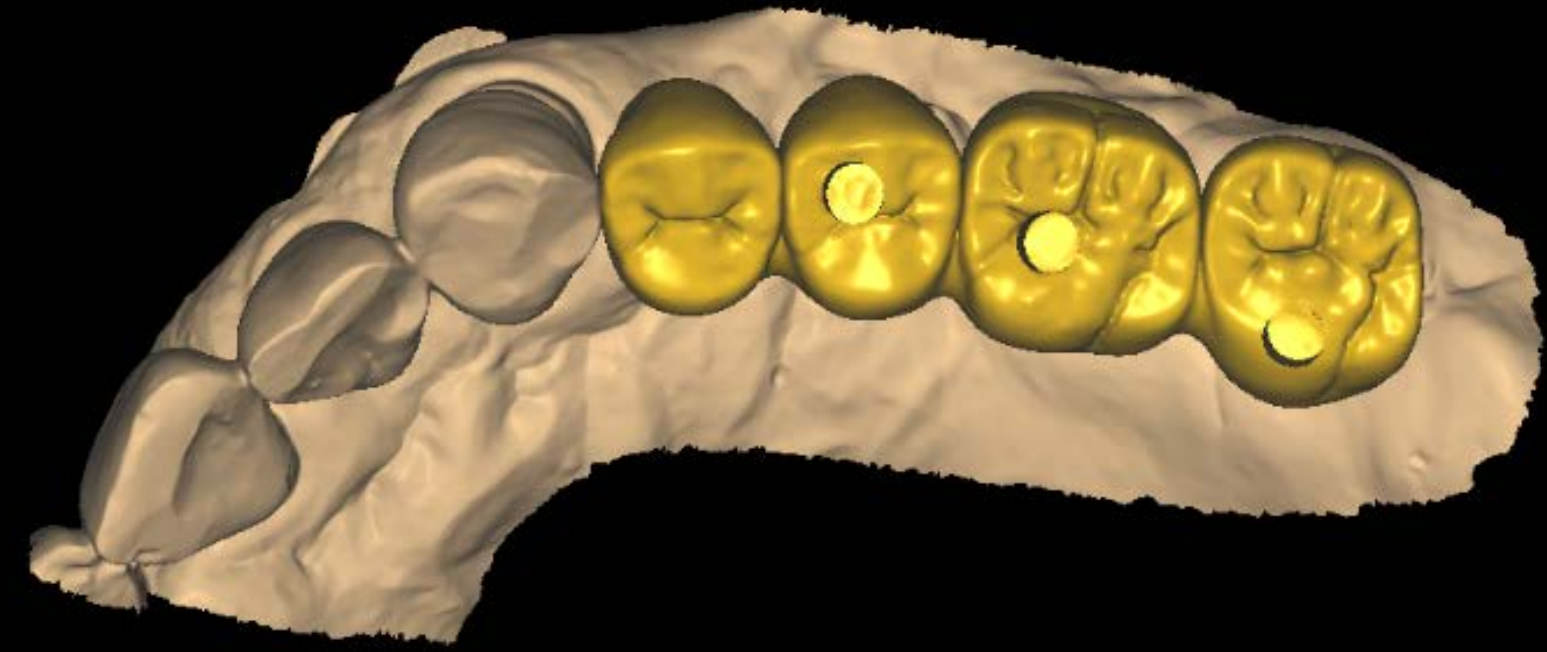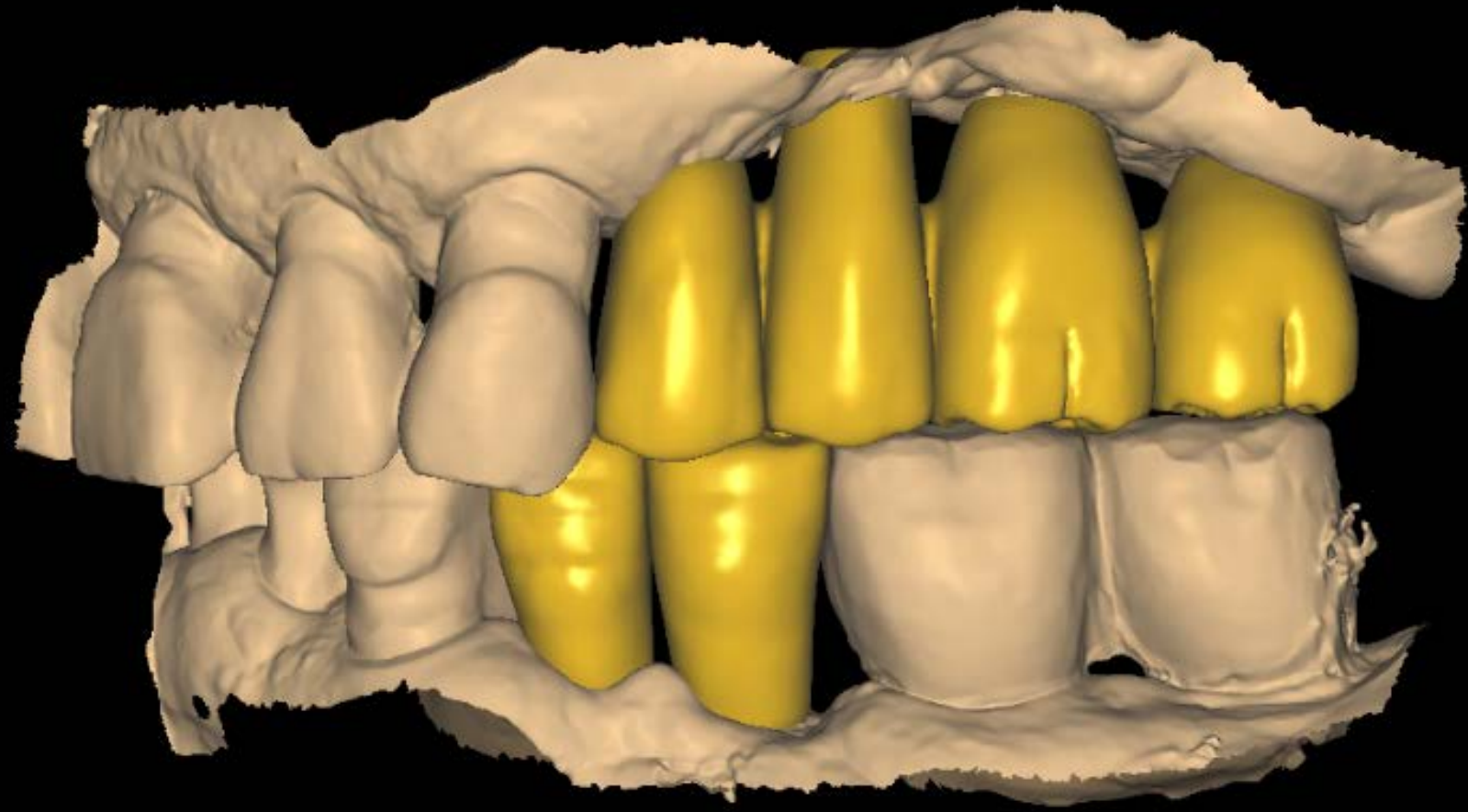



Final prosthesis





BTS


Putty + light body  
Easy and precise

기존의 analogue impression

Putty + regular body  
Difficult and occlusal error



Facial scanner  
IOS  
Virtual set up

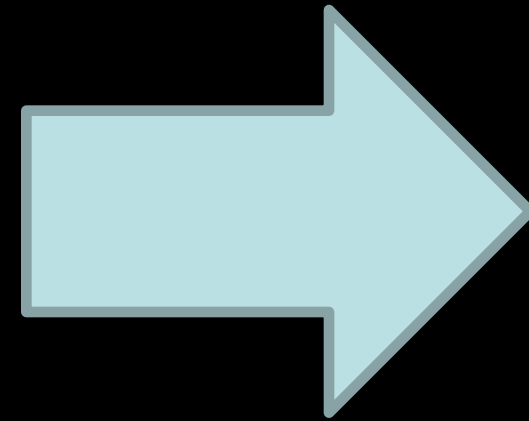

Very good with pts consulting and systemic tx

Dentium





Intra Oral Scan



Dentium





# Impression scan for working model

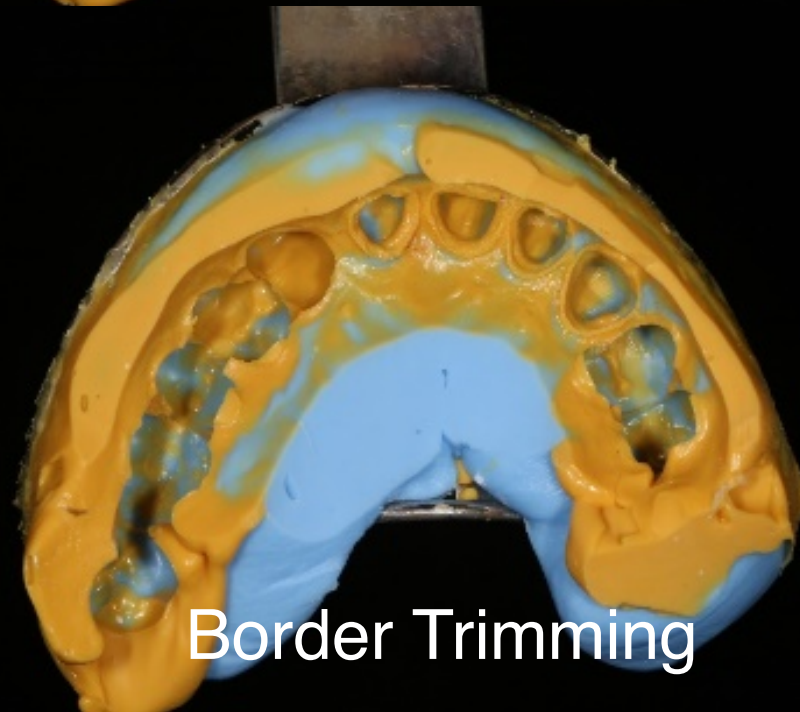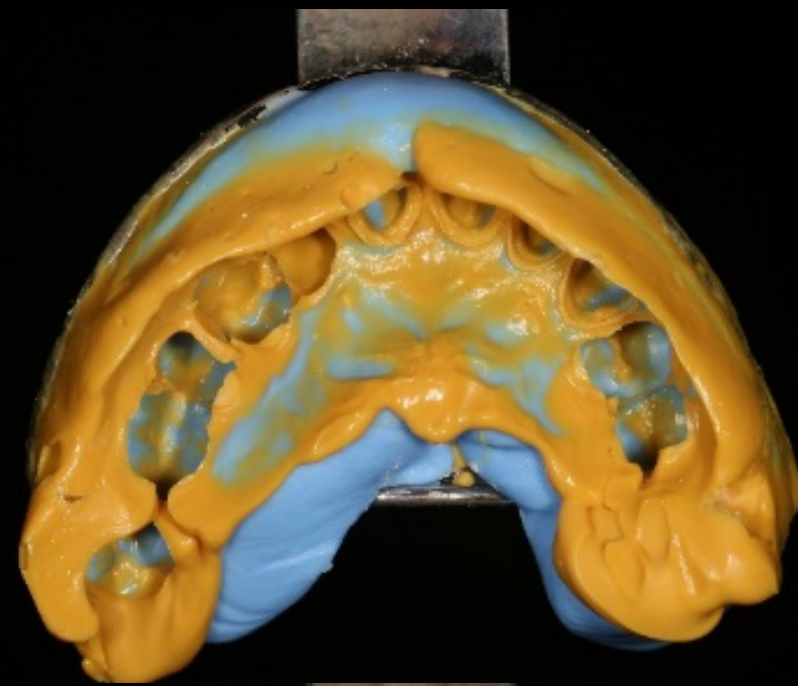

Border Trimming






Scan spray

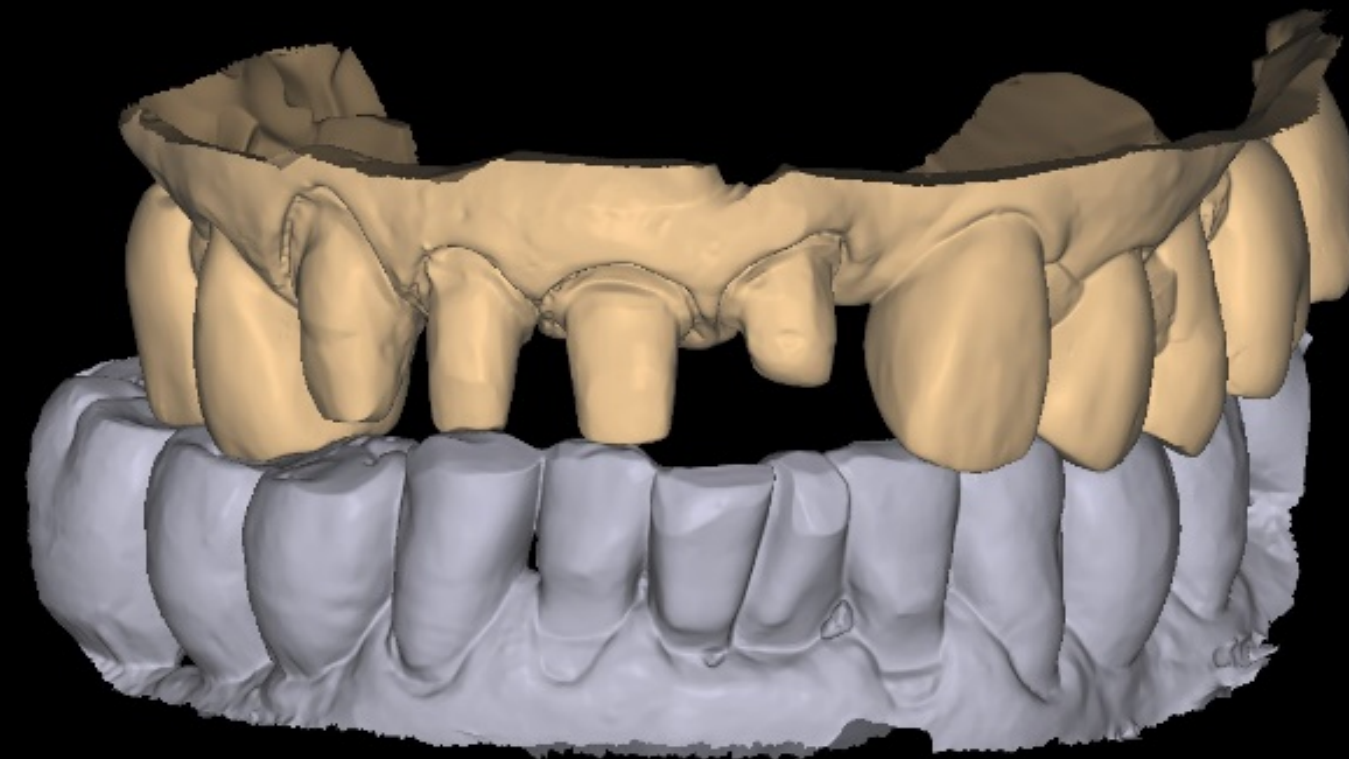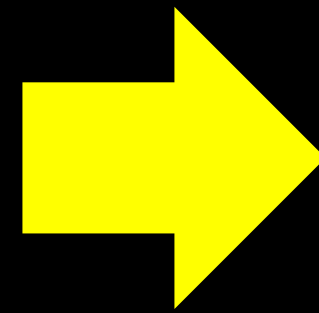

Dentium



Prosthesis design

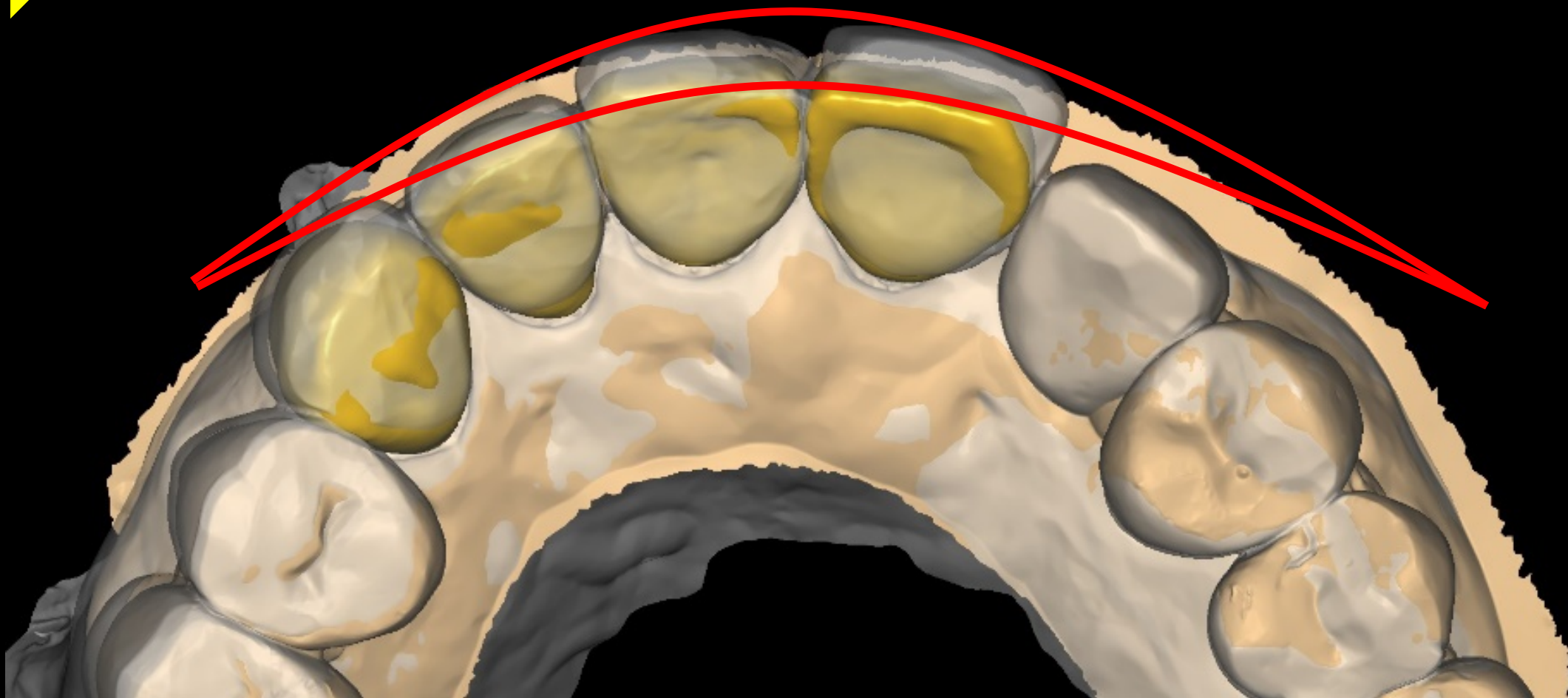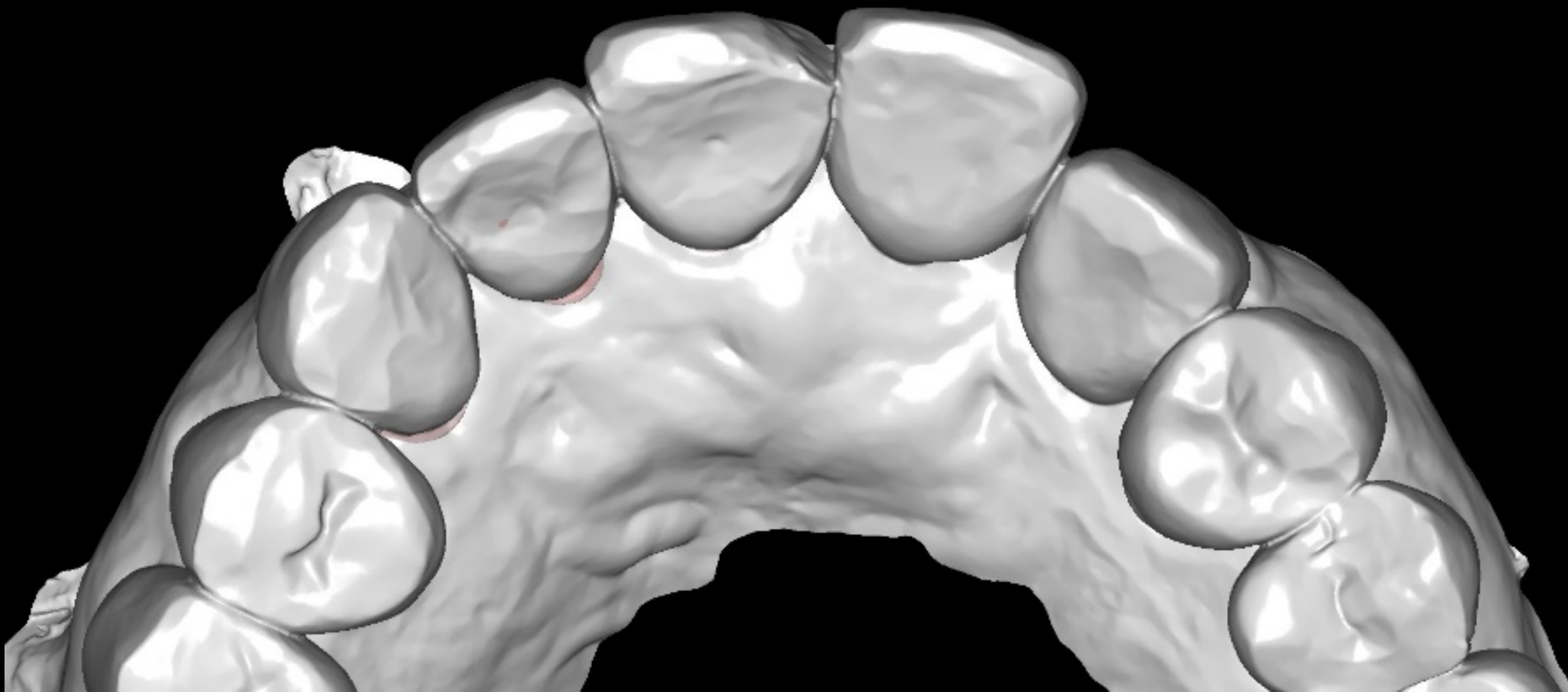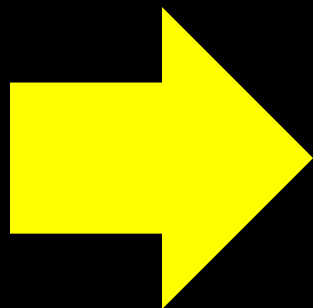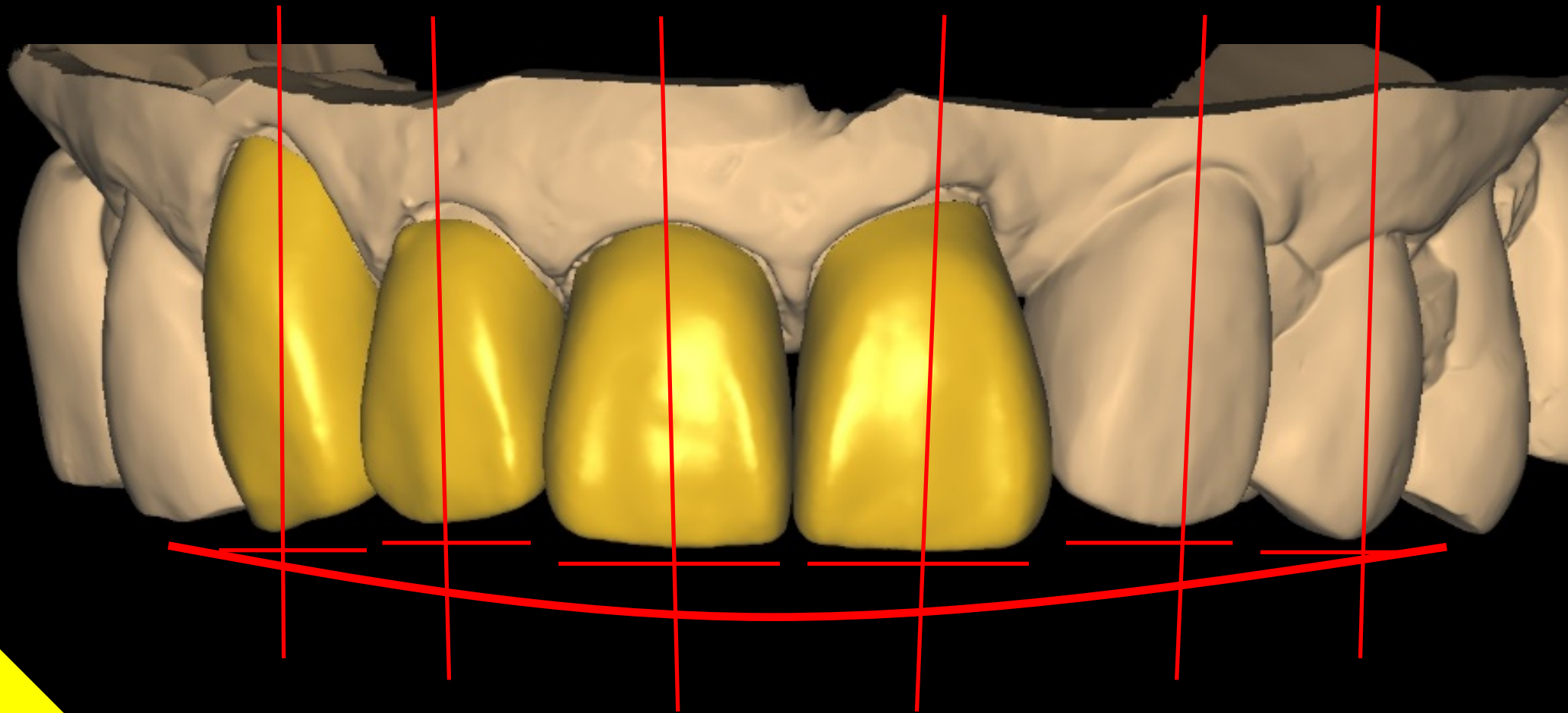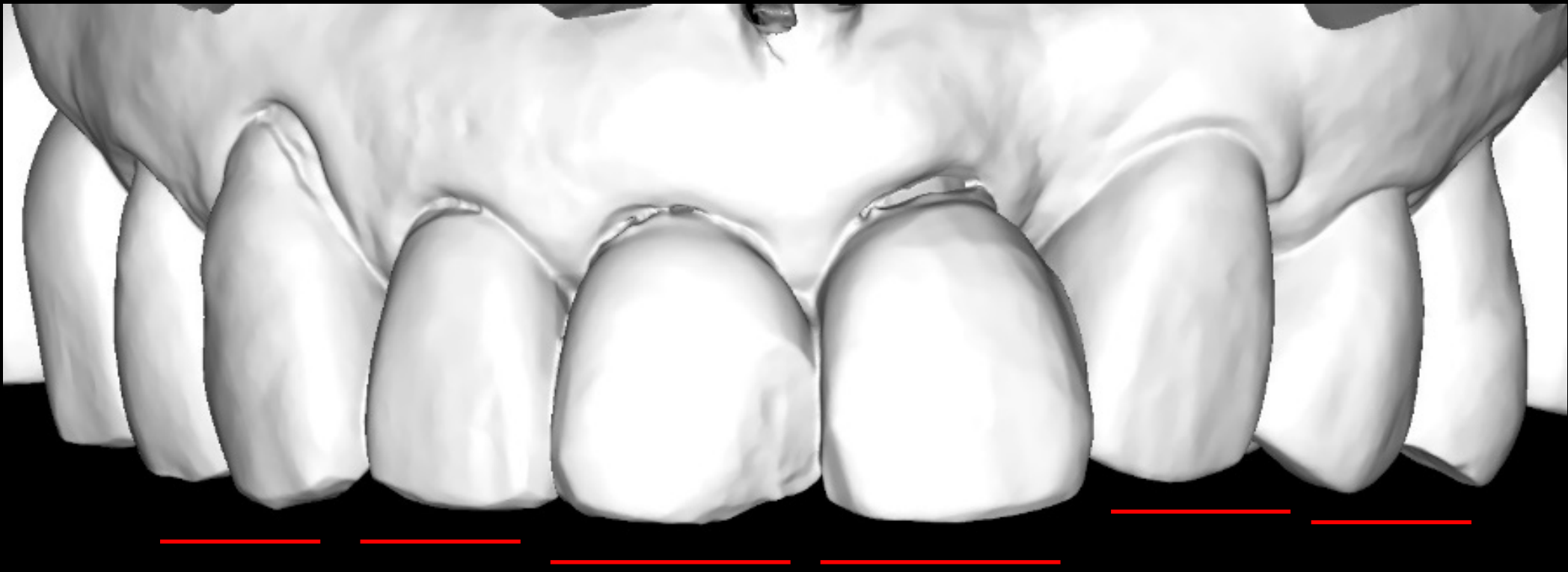

Dentium




Prosthesis design





Facial scanner

Good for pts consulting

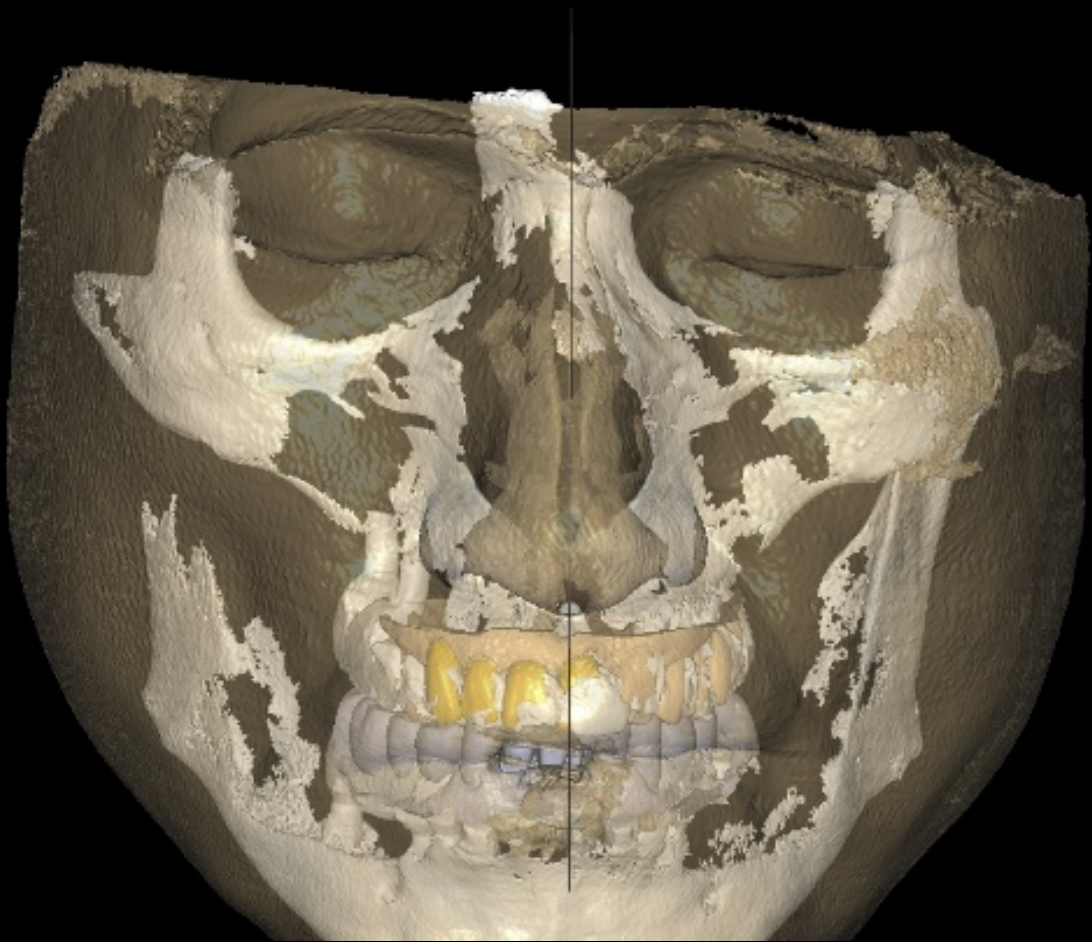

CT STL



More precise for lab

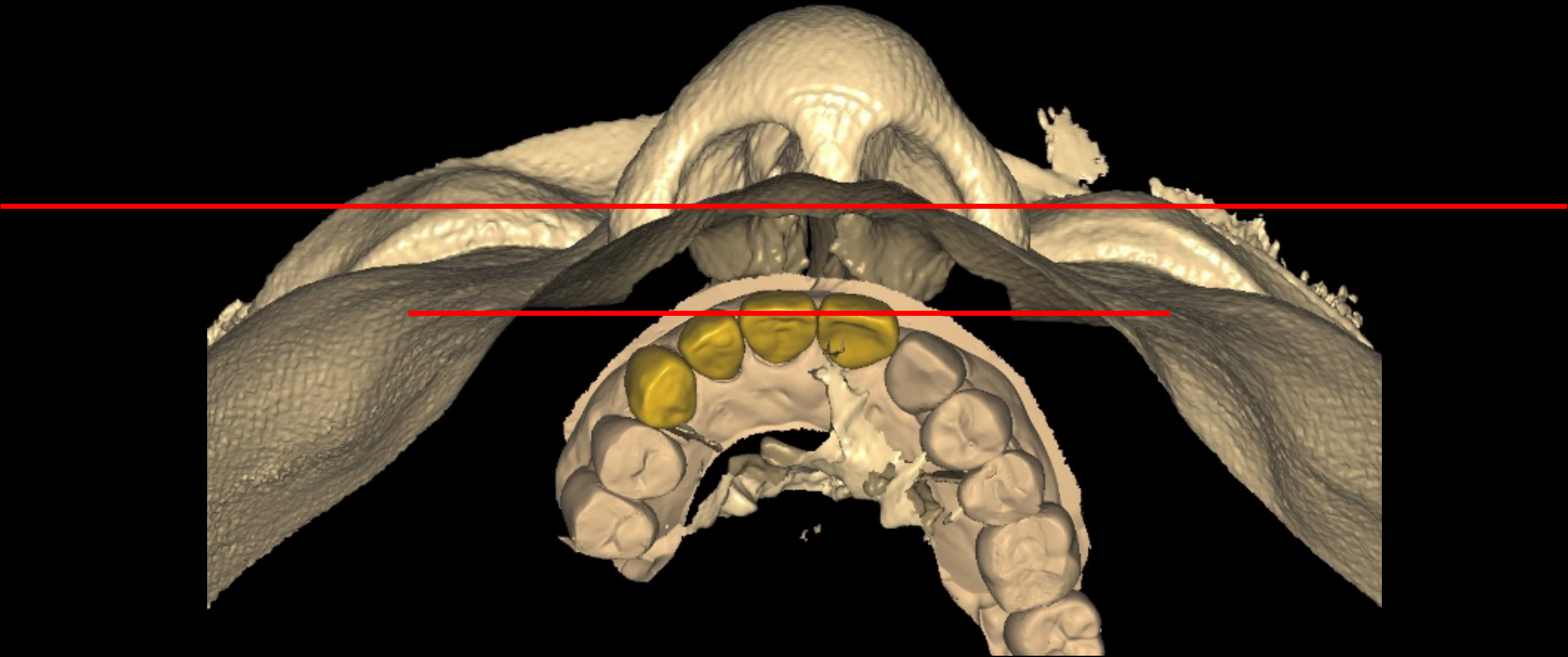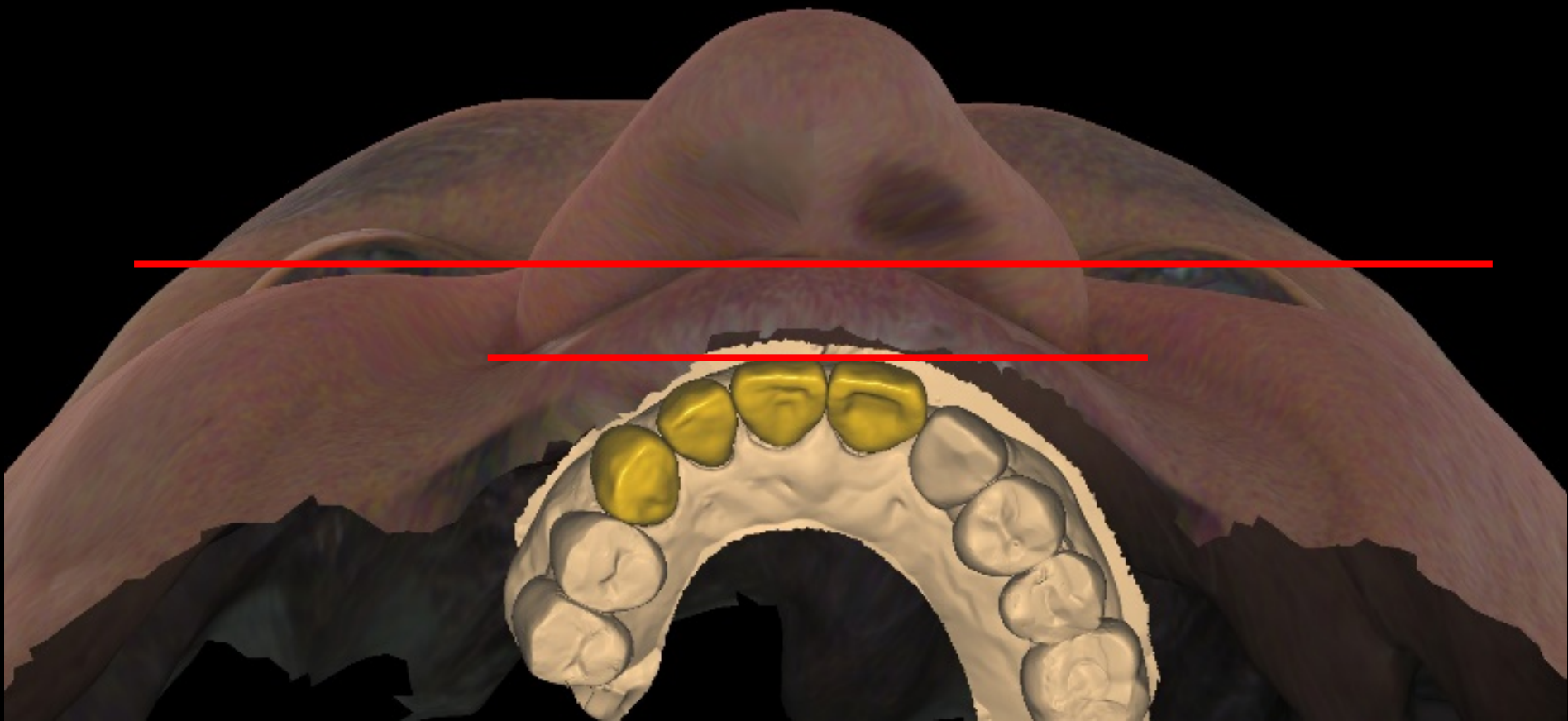

Dentium





# Prosthesis design

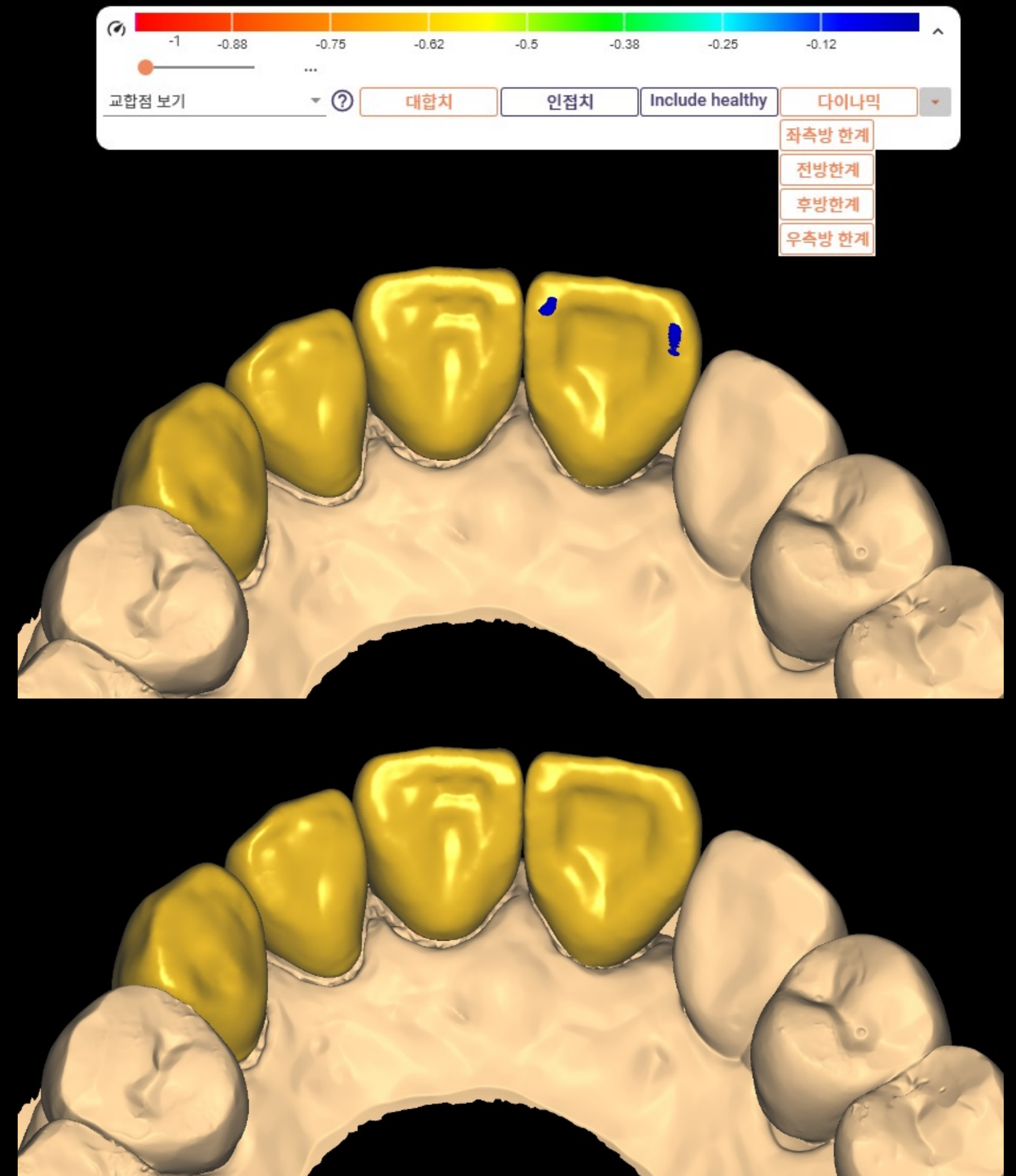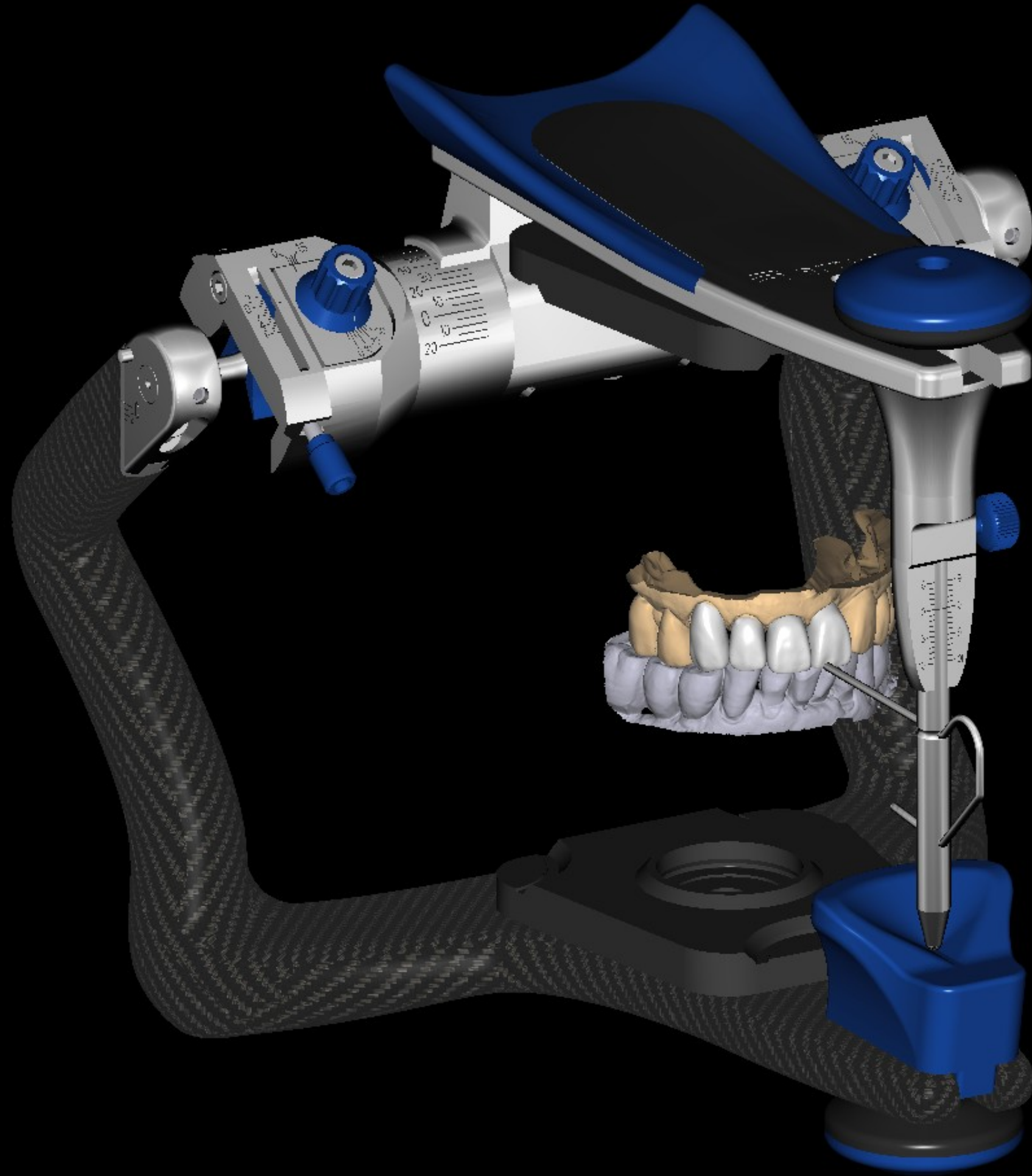

# Dentium



Prosthesis design

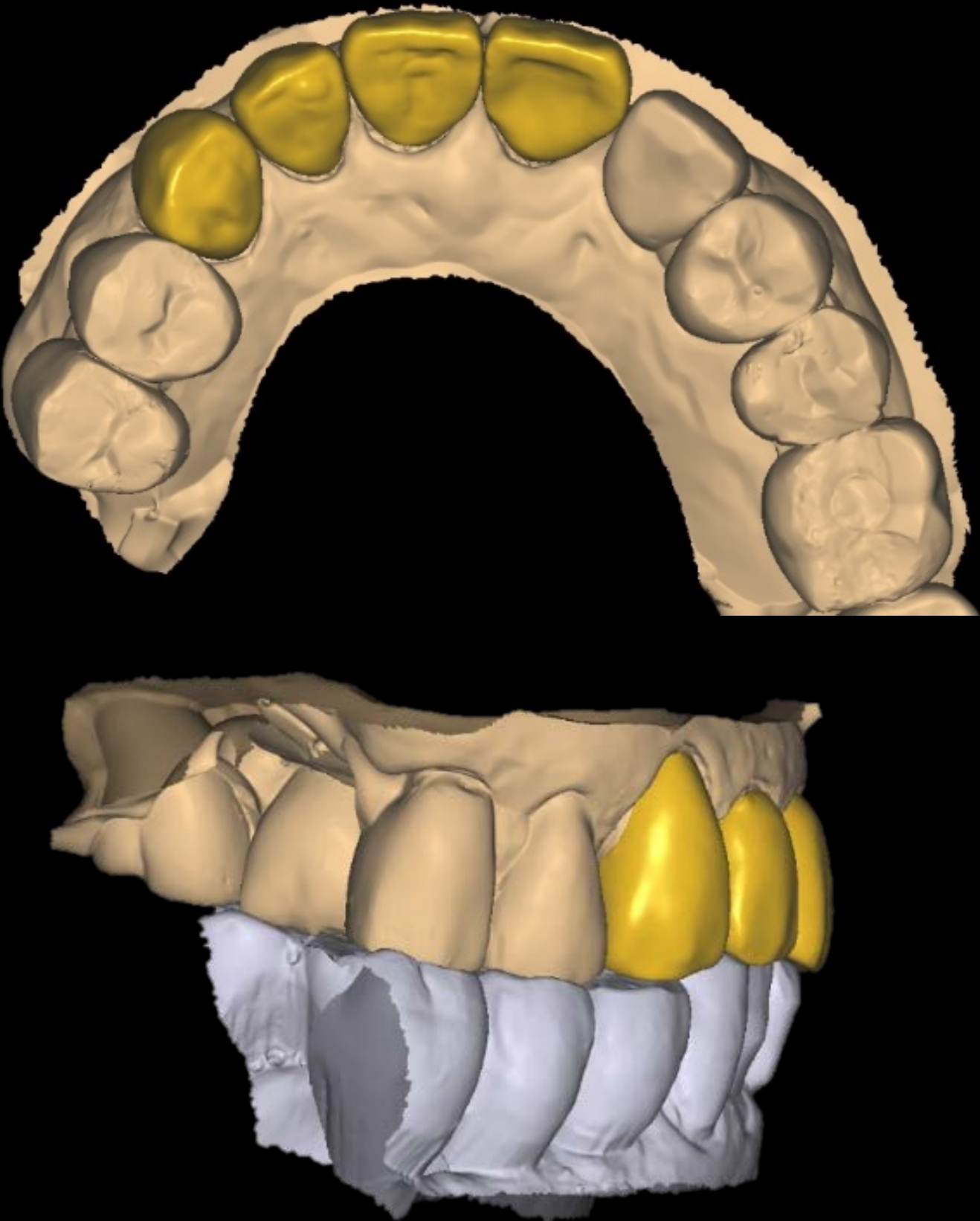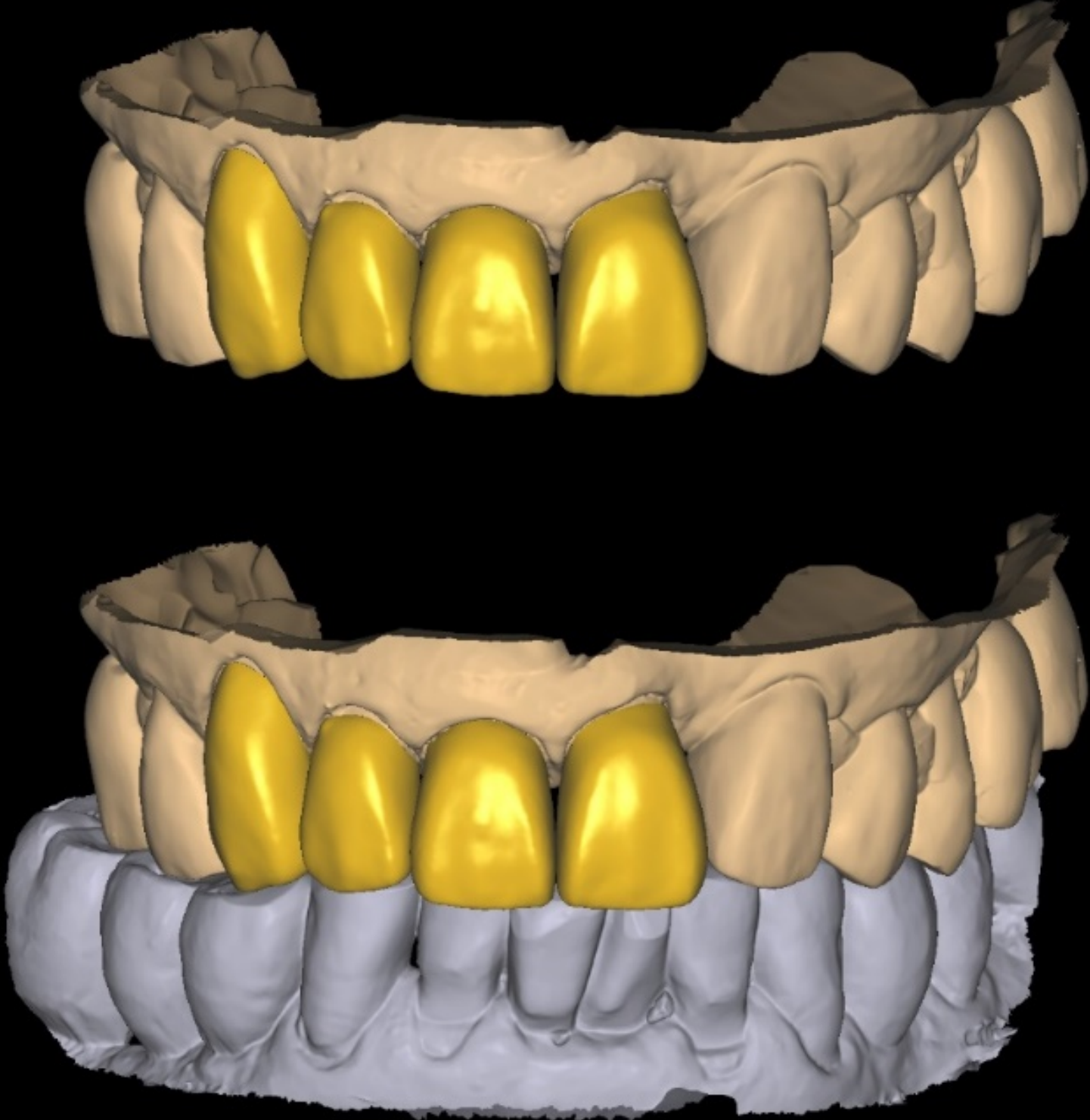

Dentium



## Shade taking (A3)



Dentium



Final prosthesis (2024-02-24)



Dentium



Final prosthesis (2024-02-24)





Before



After




1 Month Later






# Nose

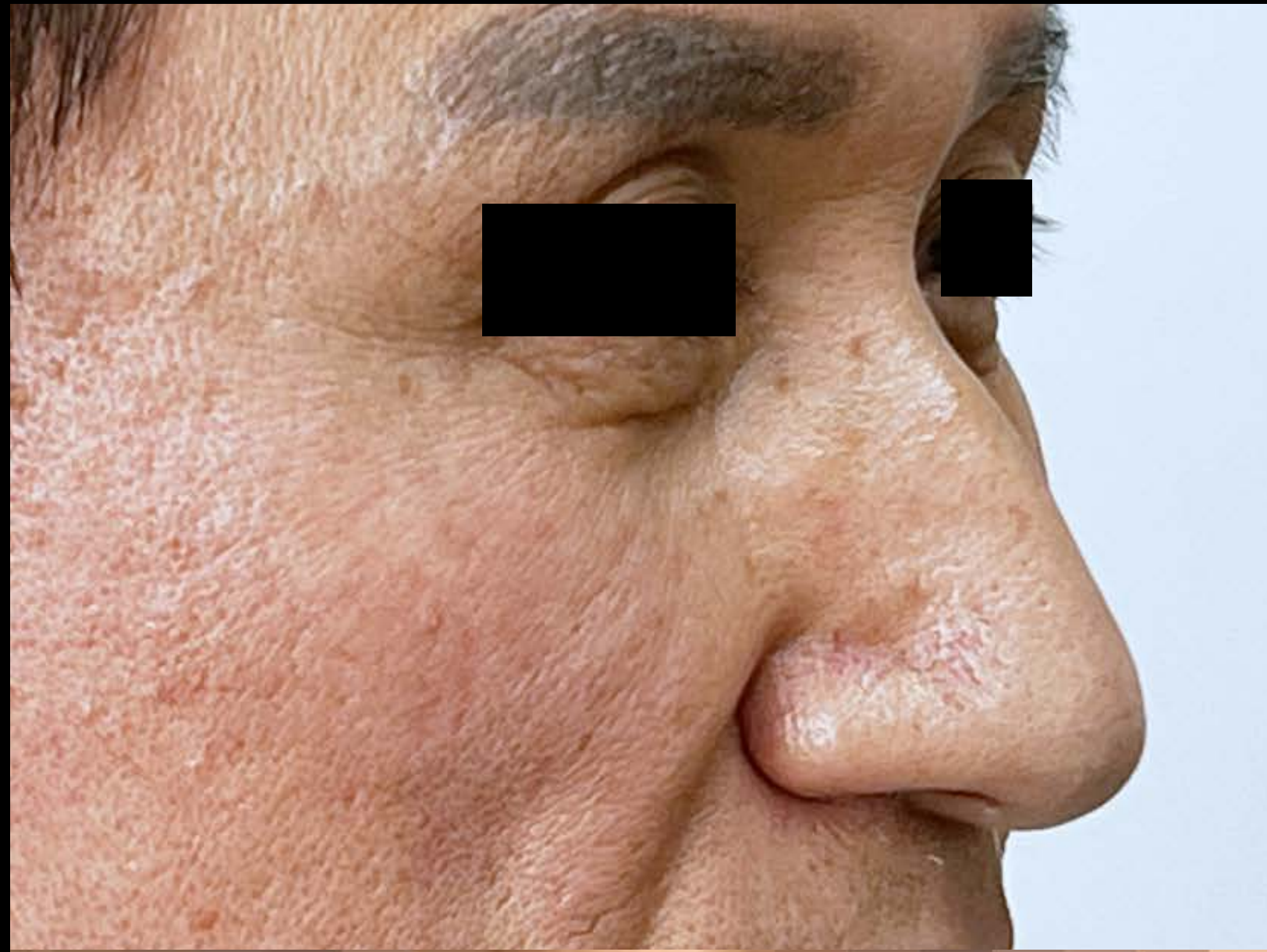

Before




After

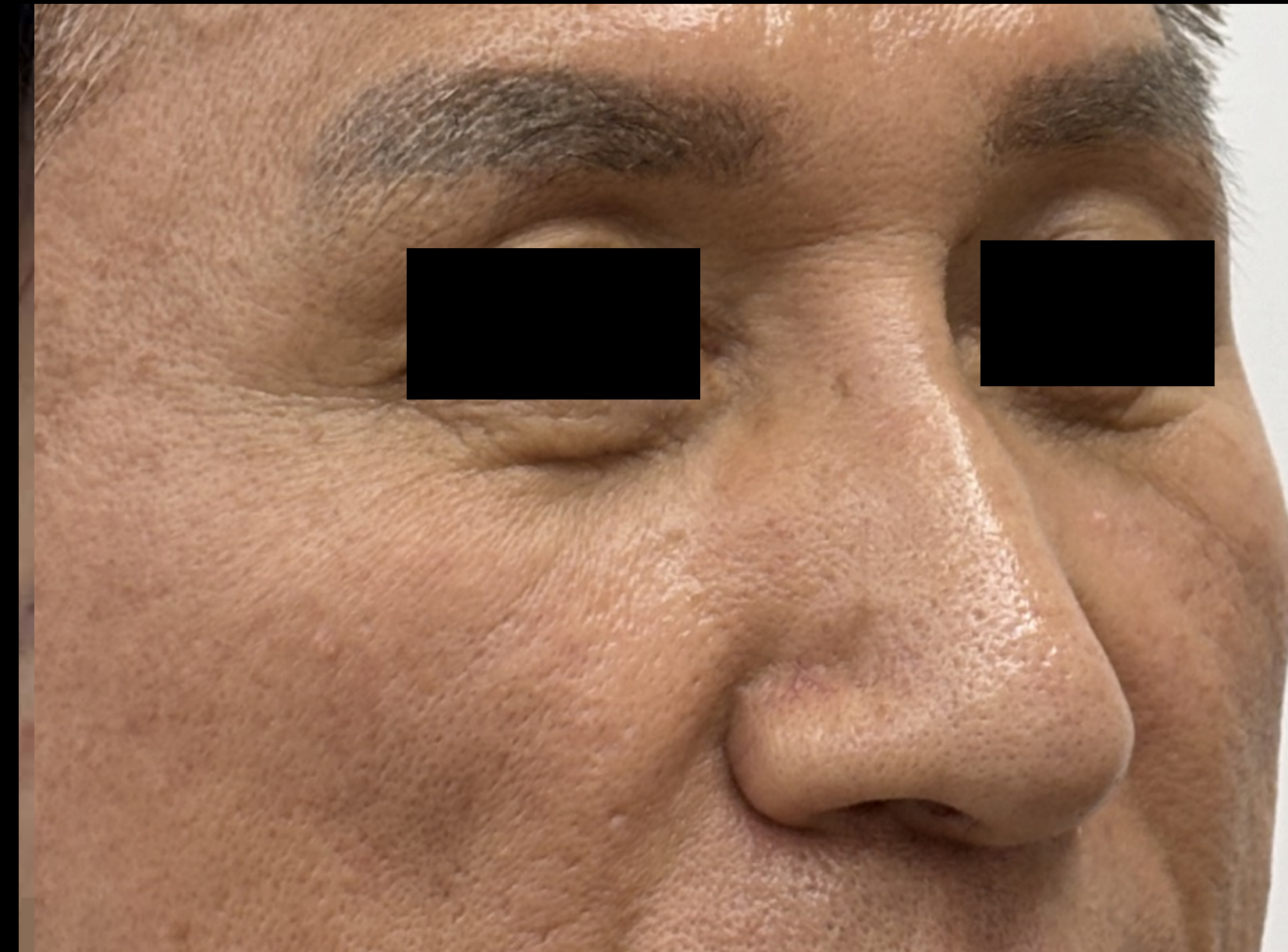

1 Month Later



# Nasolabial folds, Chin

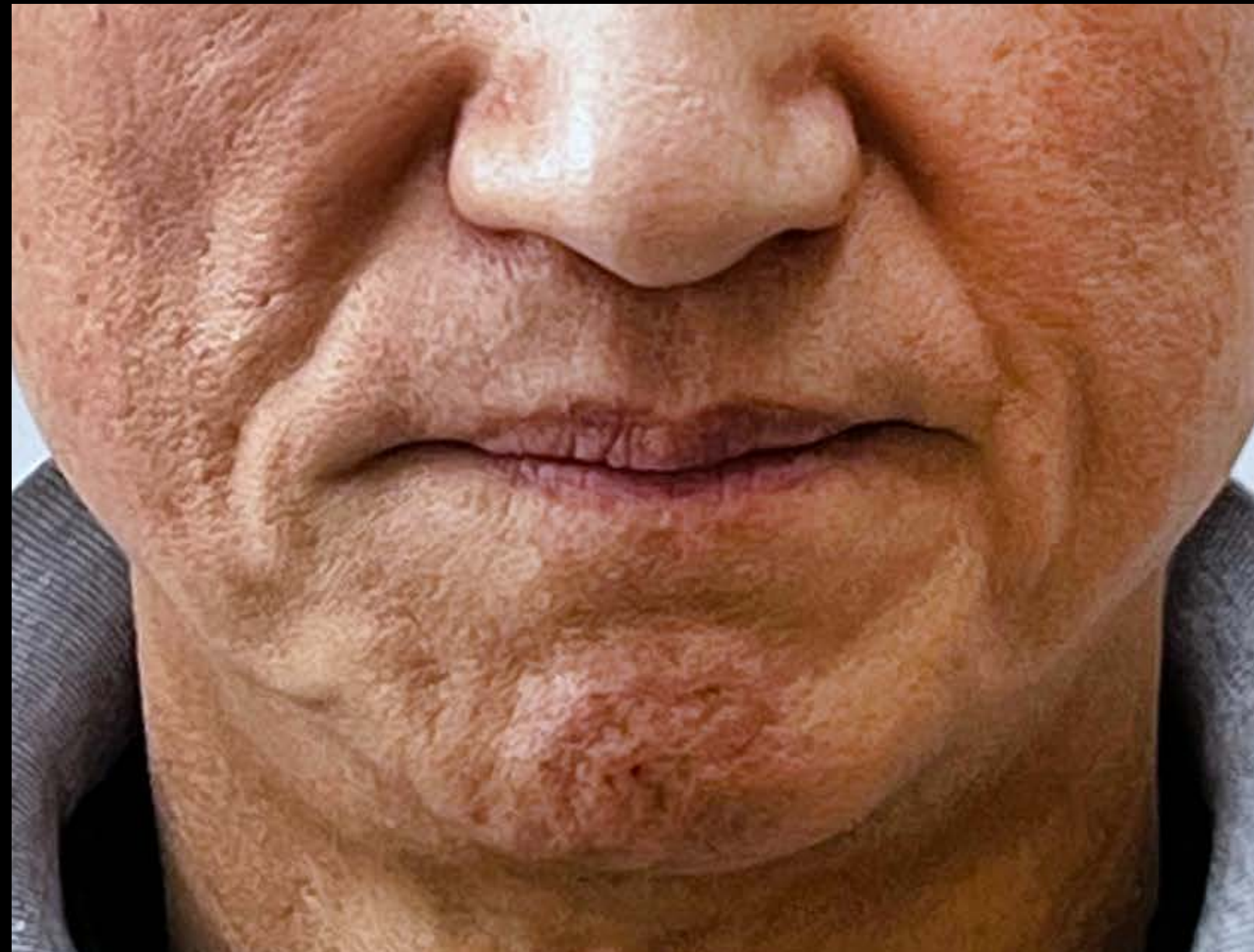

Before



After



1 Month Later



## A Change in appearance with aged person

- 1 Reduced soft tissue volume - collagen
- 2 A slight decrease in bone tissue
- 3 Reduced elasticity of the skin

-> Appears in the form of **sagging, facial depression**, etc



# The concept of fillers for elderly patients

## 1 Injection area

- 1) Nasolabial fold
- 2) Marionette line

## 2 Increase the volume


- 1) **Mainly injection** needle instead of cannula
- 2) **Mainly Bone contact**
- 3) Inject only fine wrinkles of the skin directly under the skin
- 4) Increased effectiveness when used with Botox



# MONALISA Lidocaine

Hyaluronic Acid Dermal Filler with Lidocaine

# MONALISA

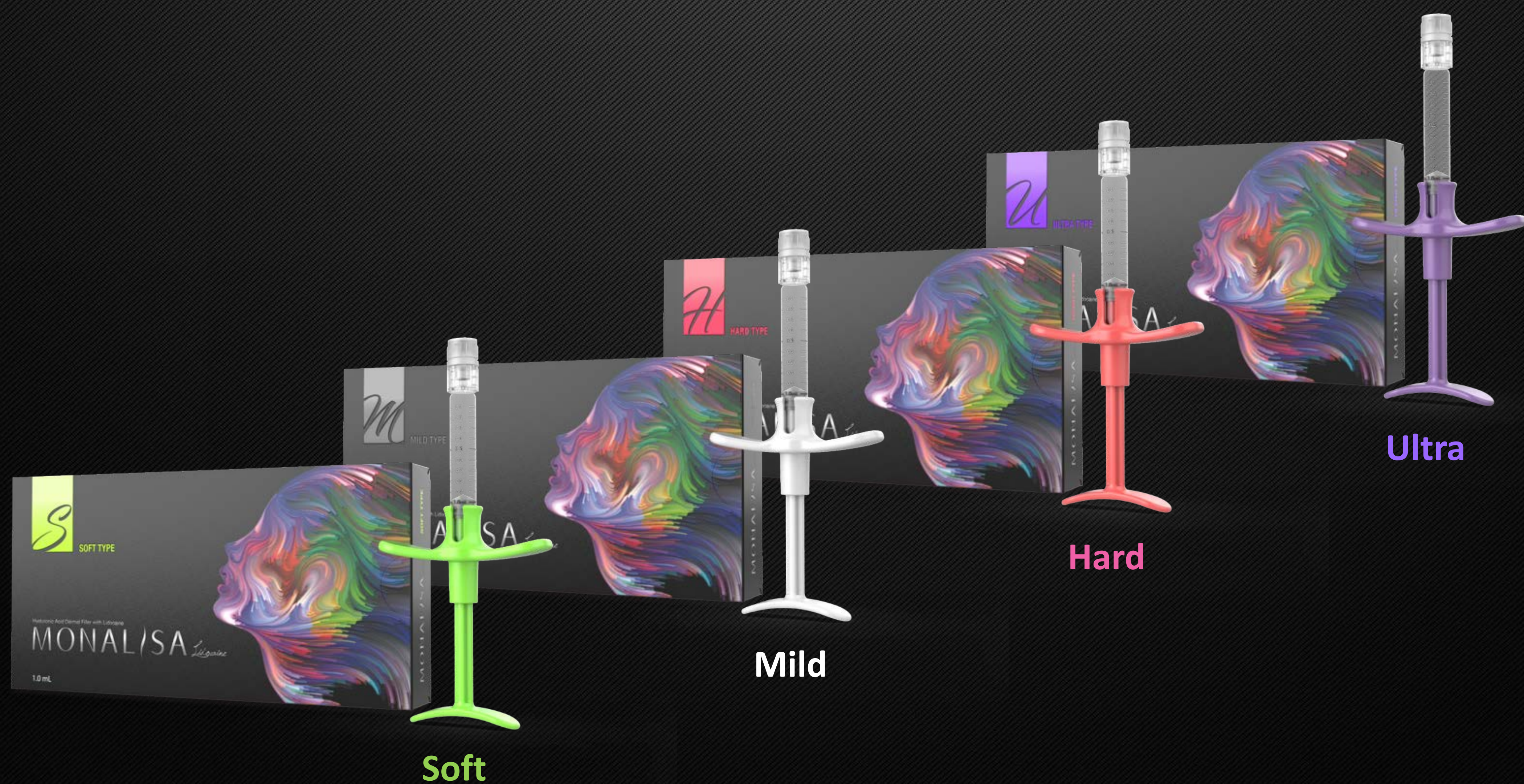

Soft

Mild

Hard

Ultra




# Physical Properties

Hyaluronic Acid Dermal Filler with Lidocaine


# MONALISA

## Particle size & shape


x 40

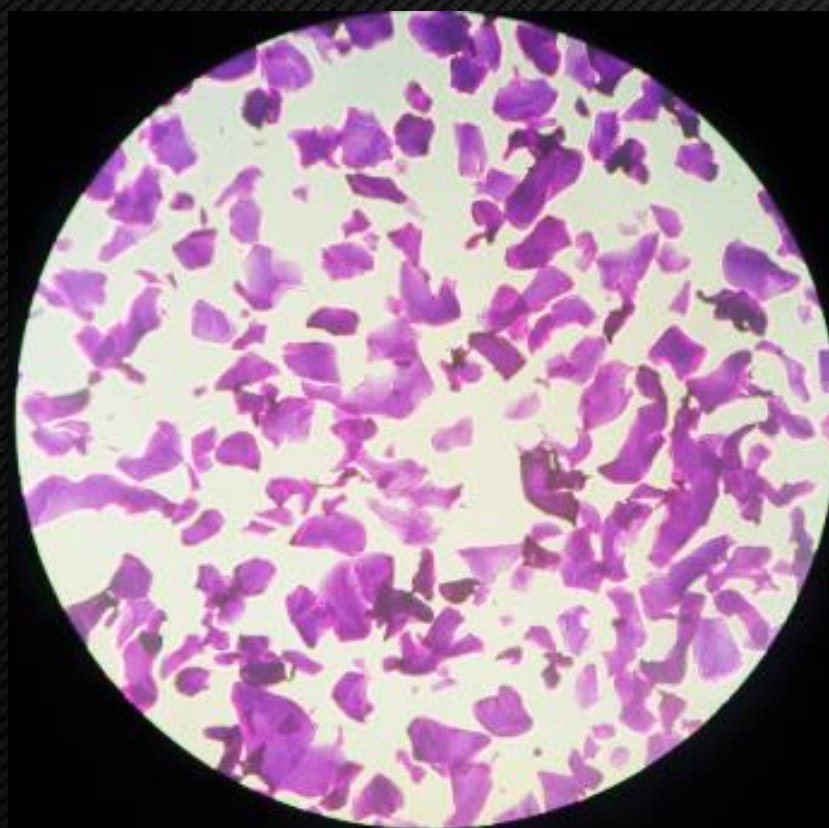


Soft

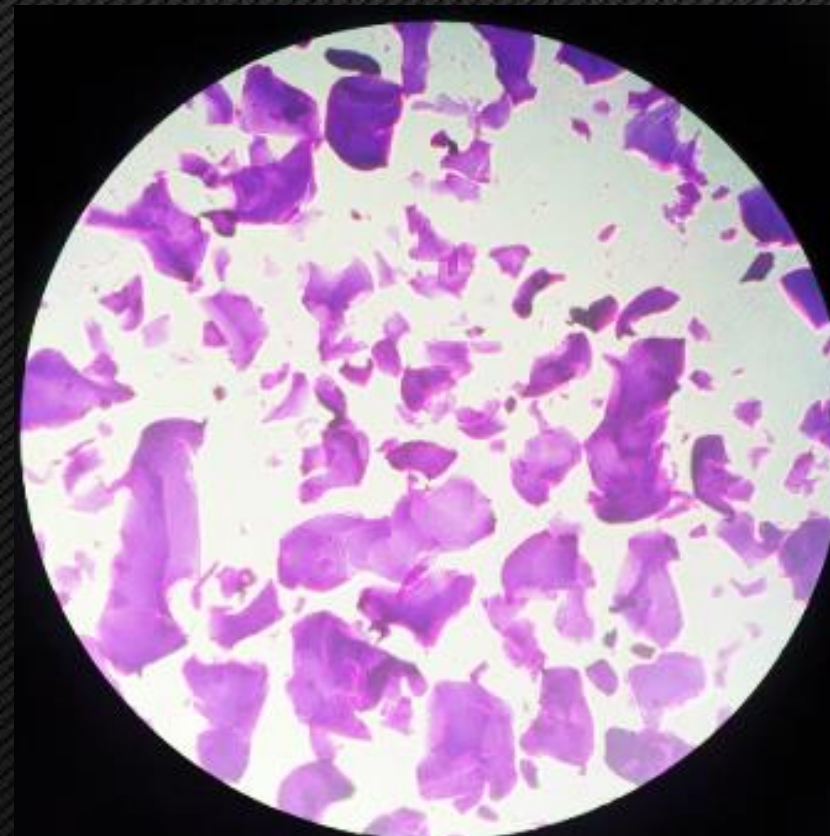

Mild / Smile

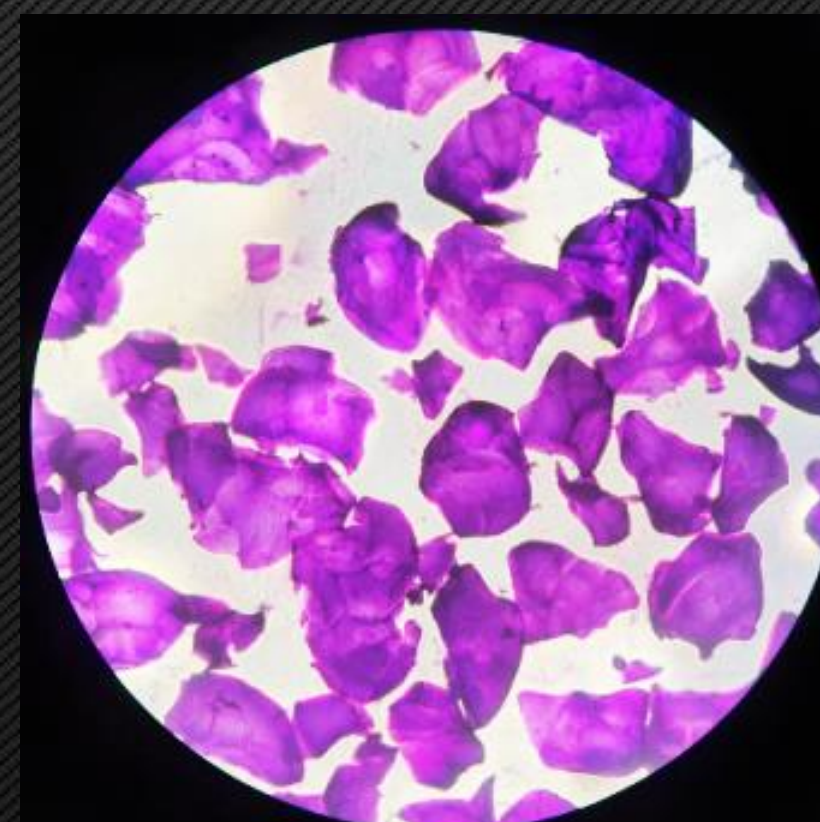

Hard

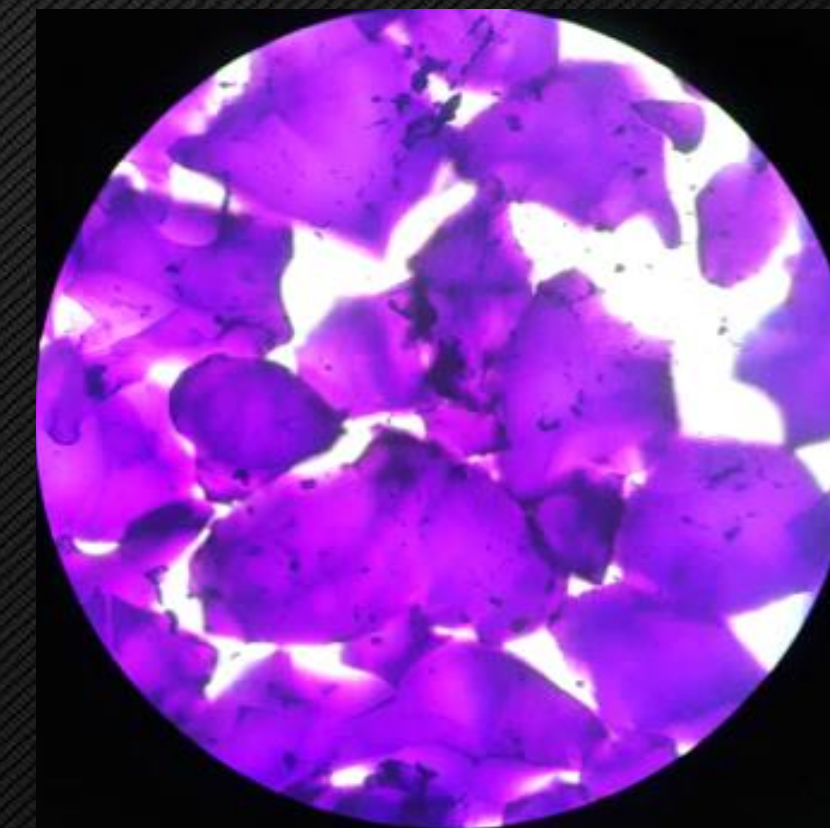

Ultra / High Elastic



# MONALISA Lidocaine

Hyaluronic Acid Dermal Filler with Lidocaine


# MONALISA




| Product             | SMILE                                               |  | HIGH ELASTIC                           |
|---------------------|-----------------------------------------------------|--|----------------------------------------|
| Form                | No impurities, transparent and colorless gel        |  |                                        |
| Components          | Cross-linked hyaluronic acid, 0.3% Lidocaine        |  |                                        |
| HA Concentration    | 24 mg/mL                                            |  |                                        |
| Particle Size       | 400 μm                                              |  | 900 μm                                 |
| Volume              | 1.0 mL                                              |  |                                        |
| Application         | Superficial dermis /<br>Middle layer of<br>subcutis |  | Deep to very deep<br>layer of subcutis |
| Needle size         | 27 GTW*                                             |  | 25 GTW*                                |
| Storage Temperature | 2 ~ 25 °C                                           |  |                                        |

\* TW : Thin wall




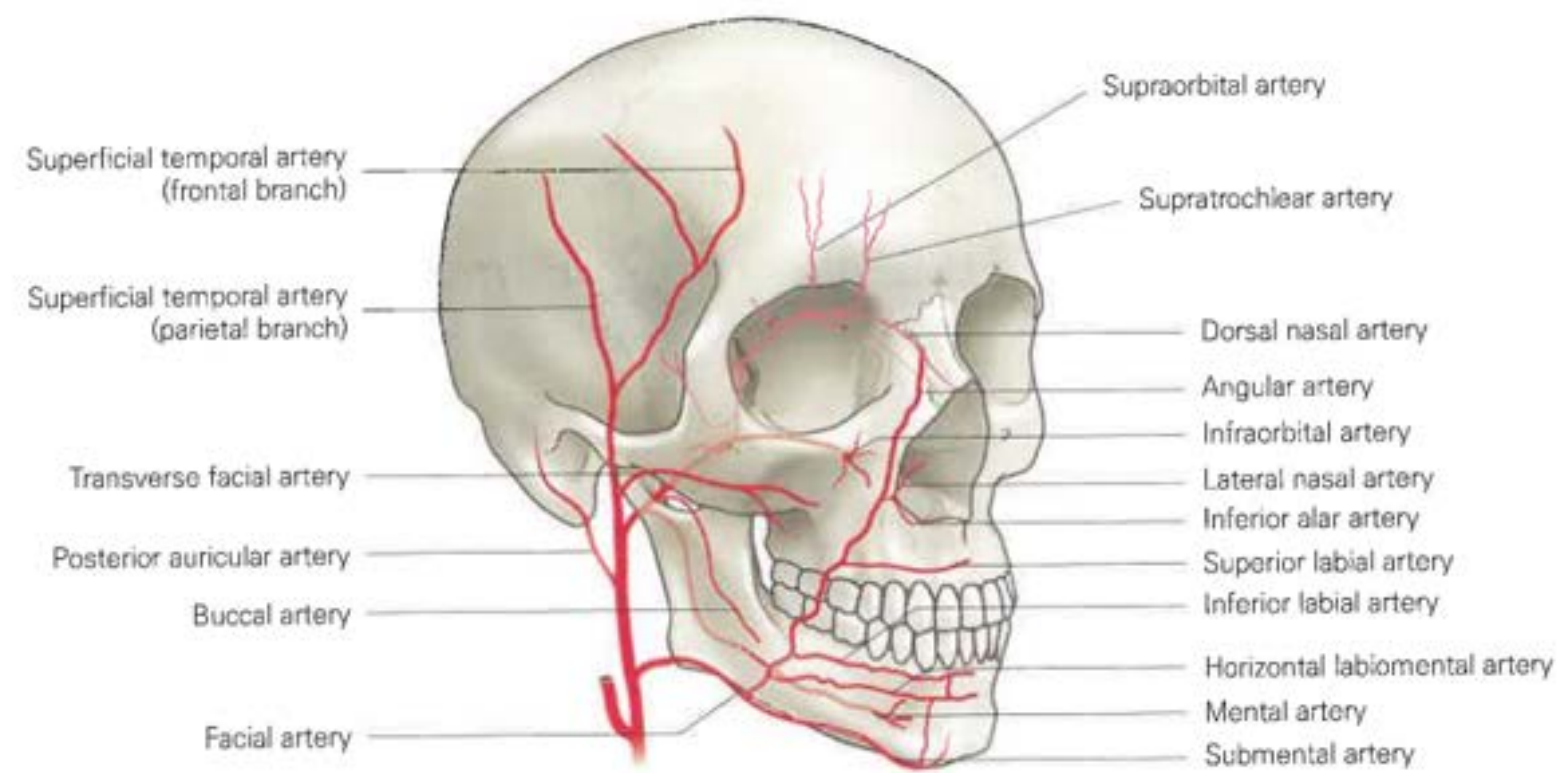


▲ Arterial branches of face

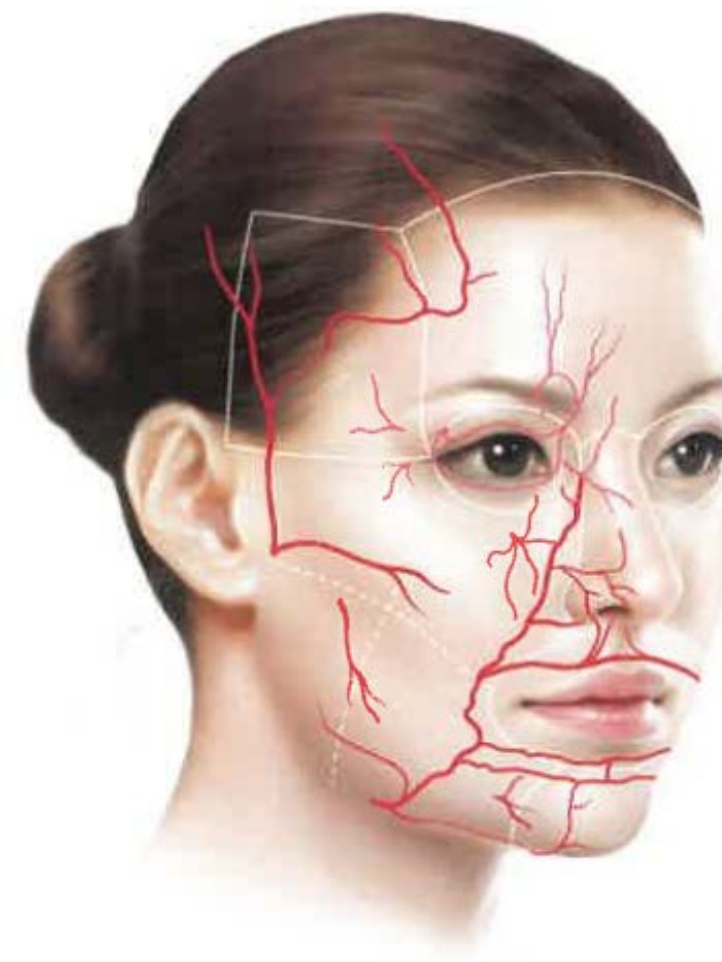


▲ Arteries supplying facial regions

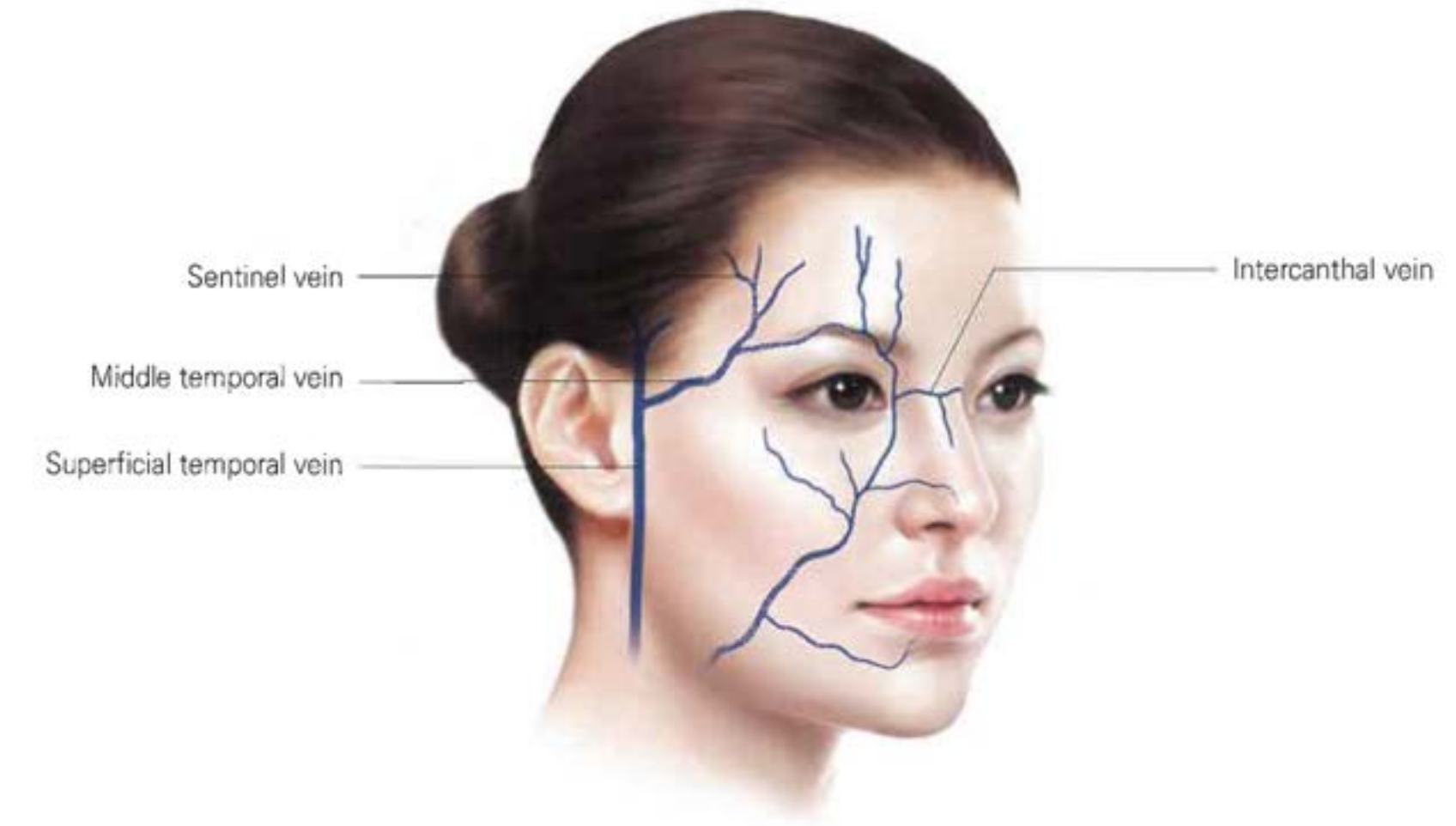

▲ Veins of face

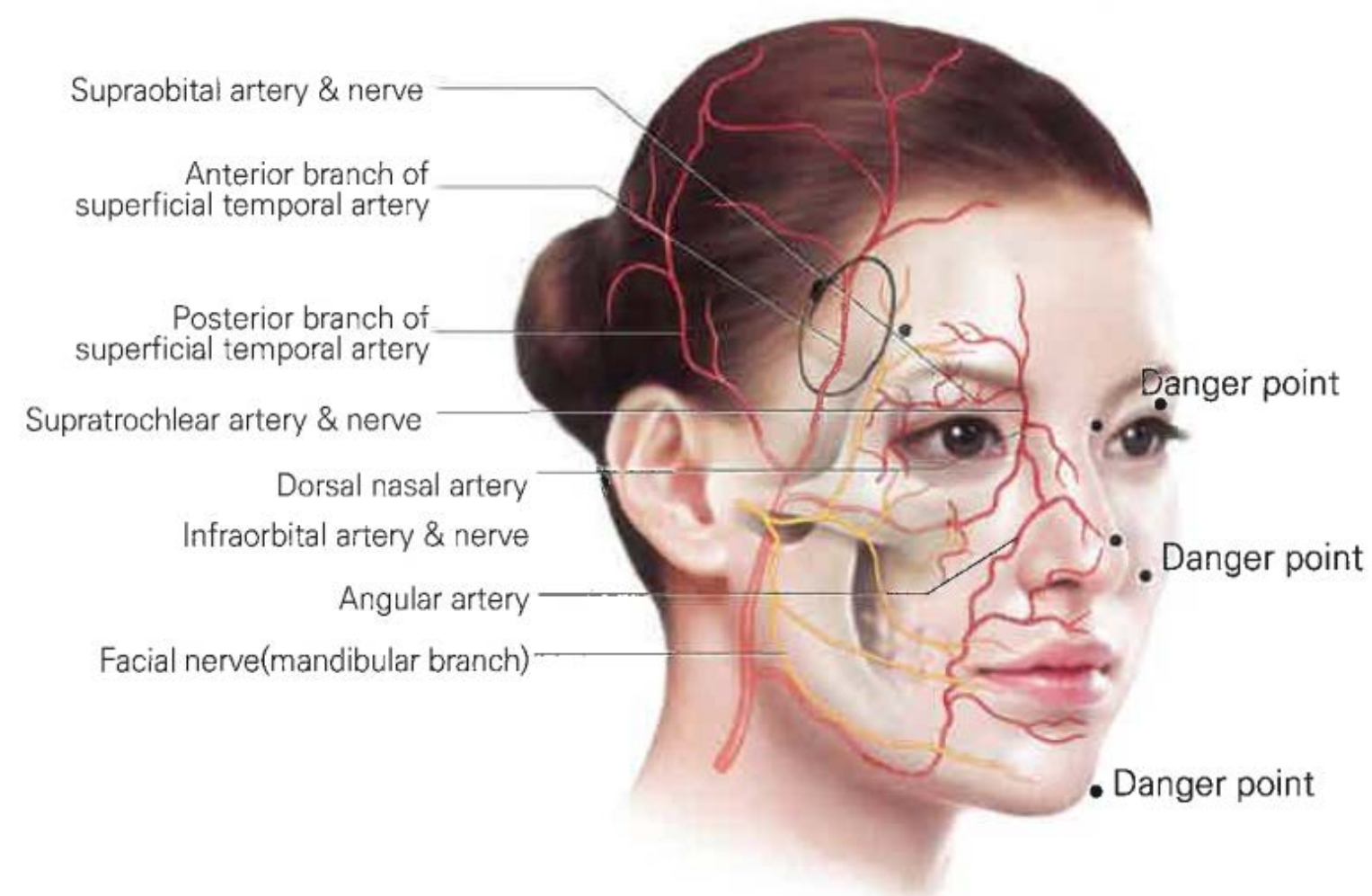

▲ The blood supply of lateral forehead area

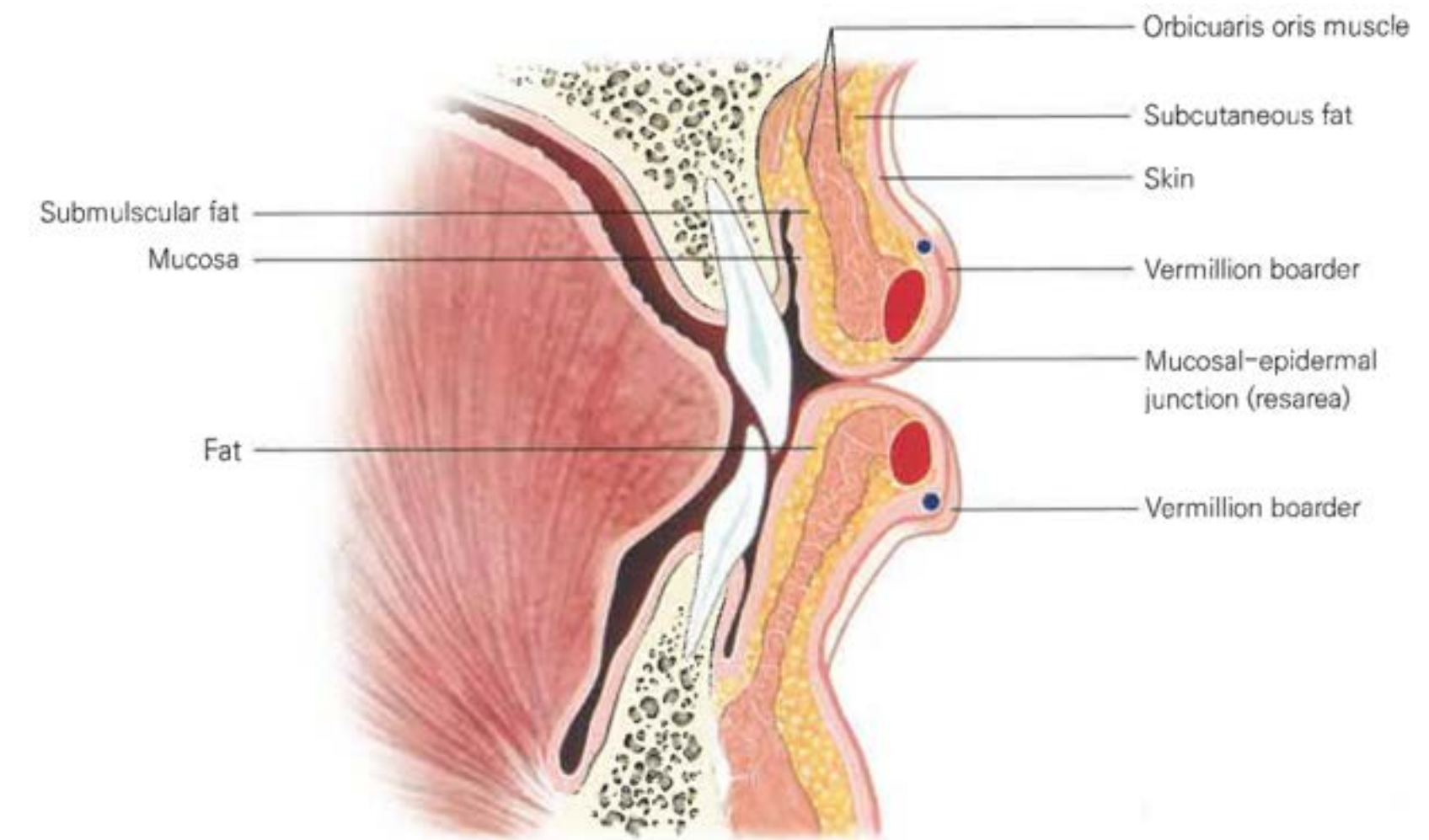

▲ Sagittal section of the lip



## The Side effects of fillers

1 Blindness

2 Necrosis of the tip of the nose

3 Edema

1)early edema

2)delayed edema



Safe filler

Natural image filler injection

Mutilayer

- bone contact
- just below the dermis



Hyaluronic Acid Dermal Filler with Lidocaine

# MONALISA



25G Needle



25G Needle / 23G Cannula



# Application site



## **Bone contact touch (Deep fatty layer)**

- Long-lasting
- Deep wrinkle improvement



## **Superficial skin (Superficial fatty layer)**

- maintain a natural and smooth appearance
- low invasiveness





# Clinical case : Nasolabial folds & Marionette Line & Lip

Before



After





Before & After





Before & After









# Bright MTA Sealer *Plus*

Root Canal Sealing Material






## Excellent sealing and regeneration

Ultra Antibacterial Effect

Calcium Silicate → Mineralization


High Flowability → Convenience Handling

Excellent Radiopaque

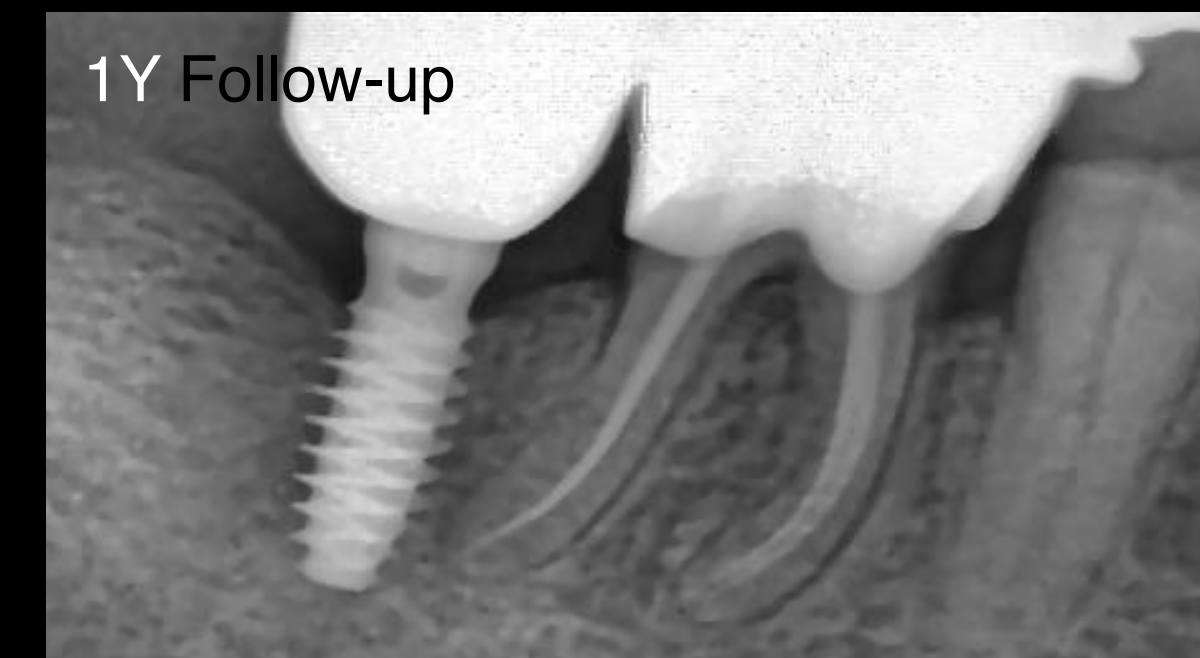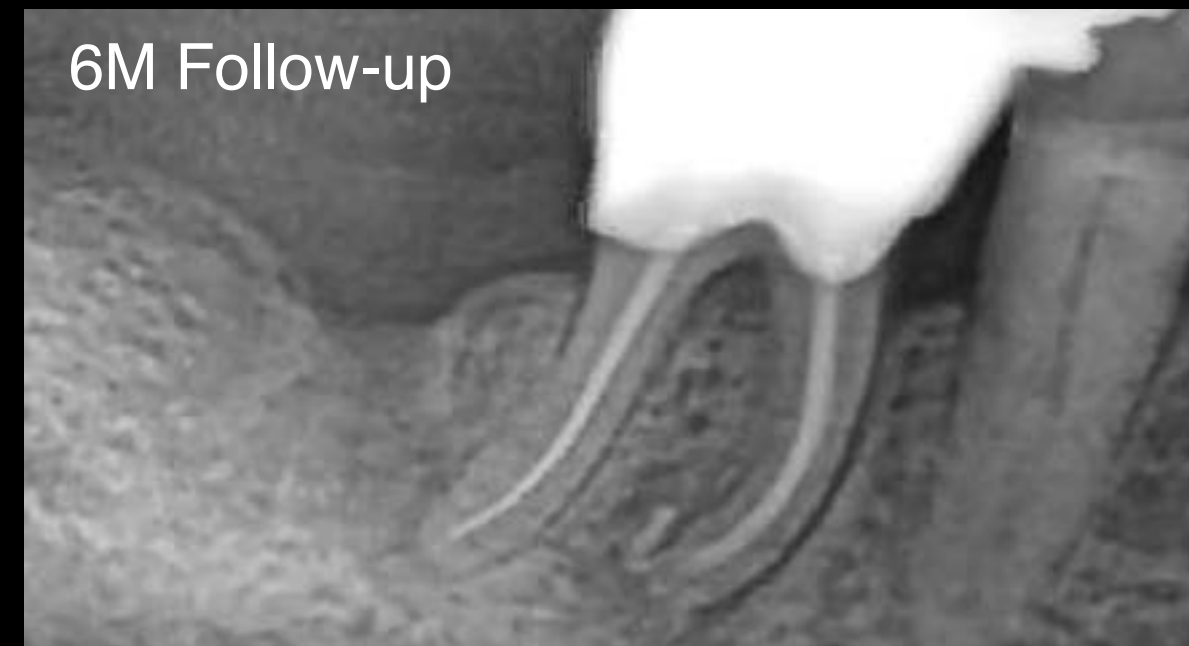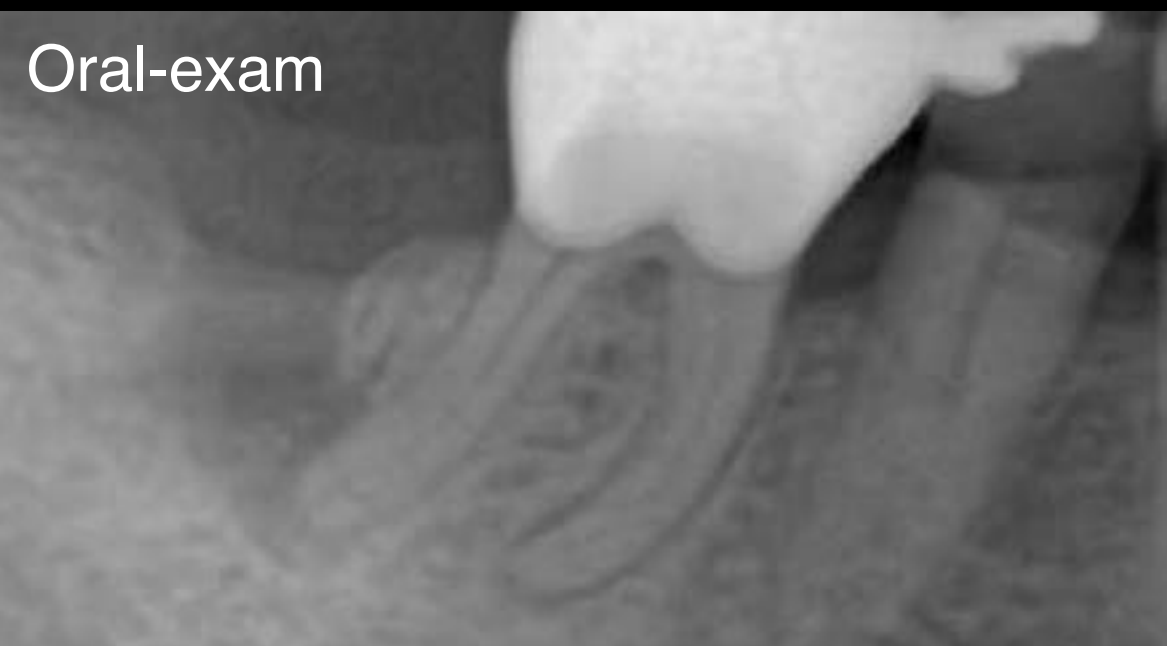




Oral exam (2023-11-18)

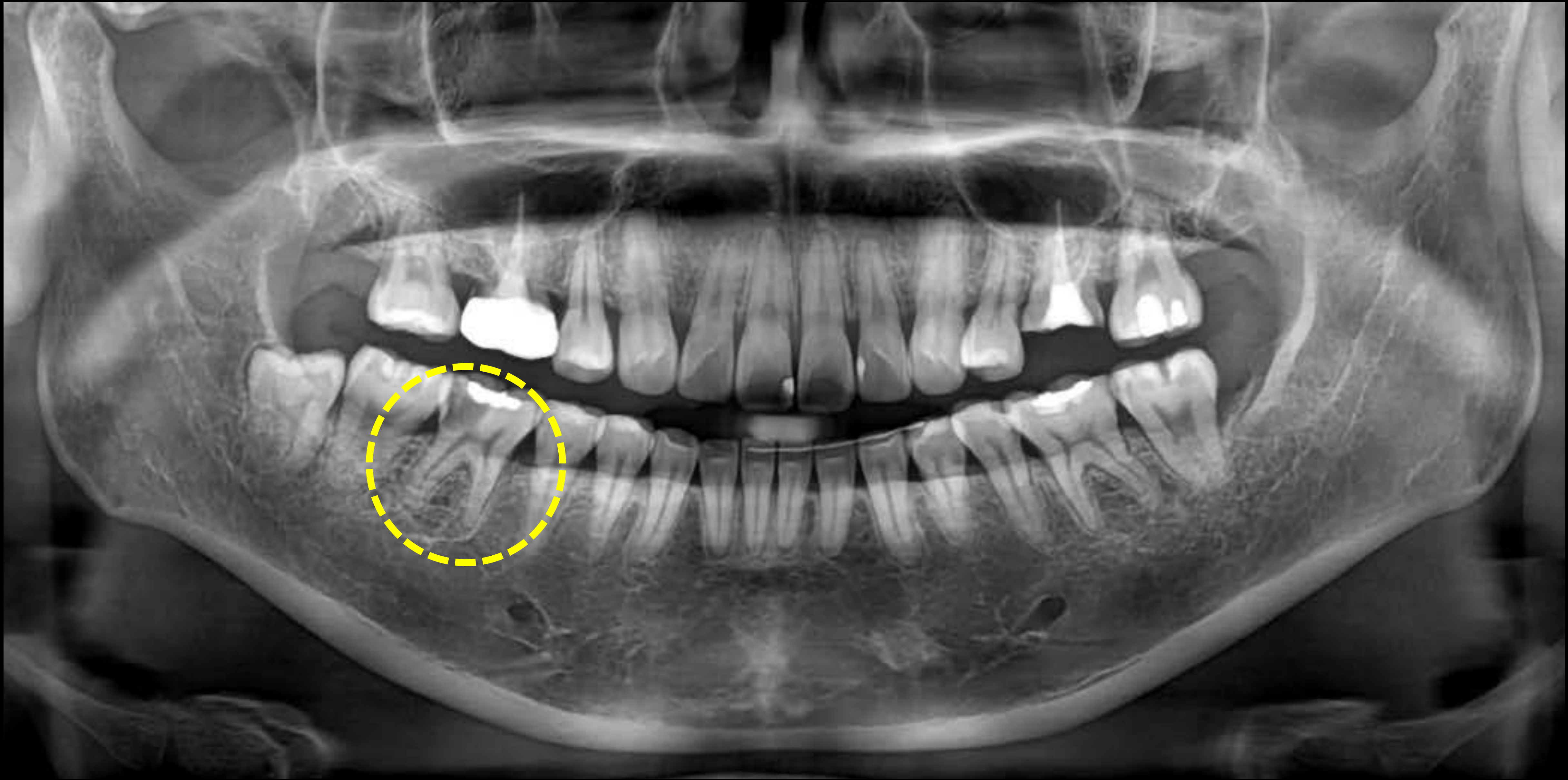




< Mesial canal >

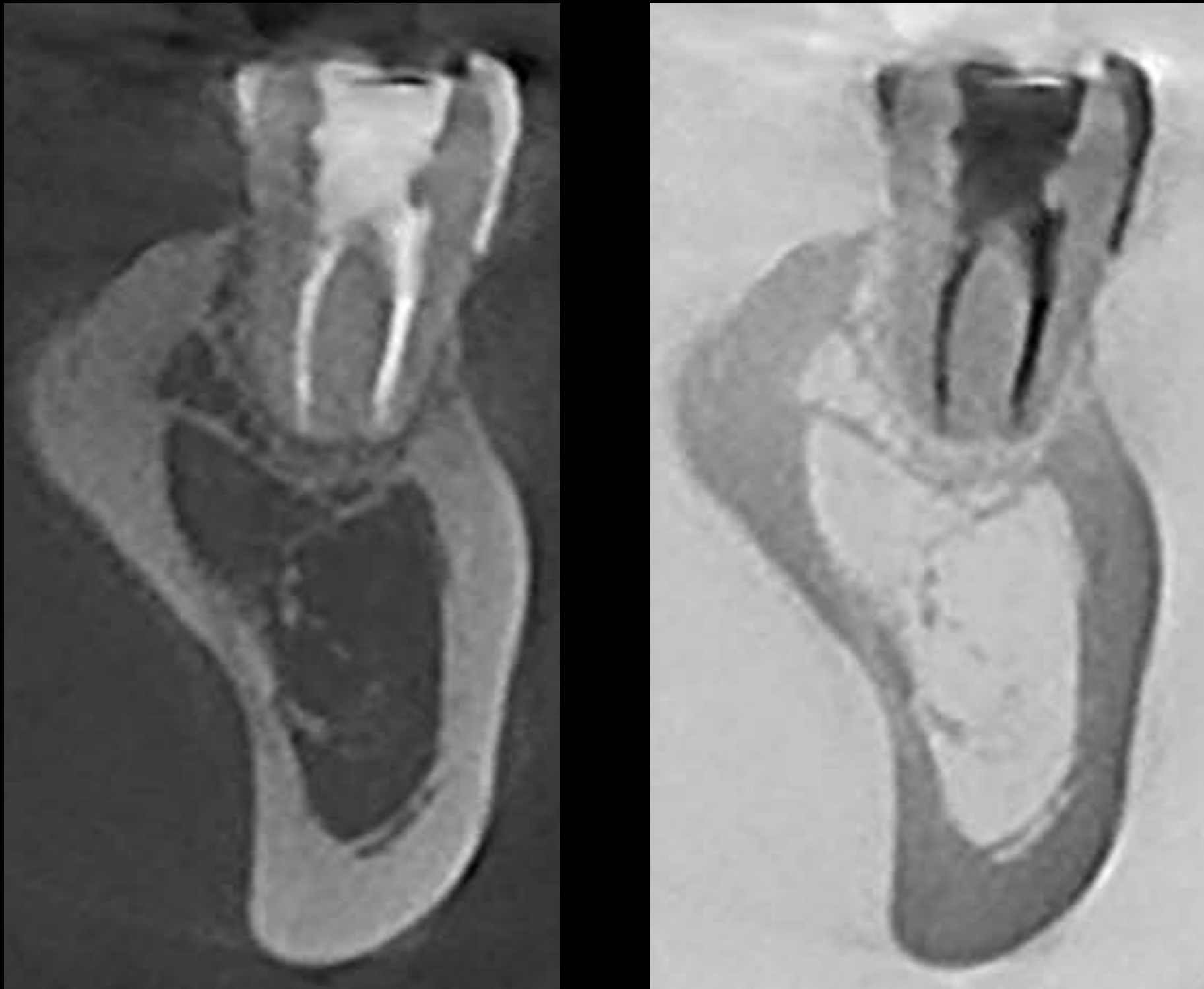


< Distal canal >

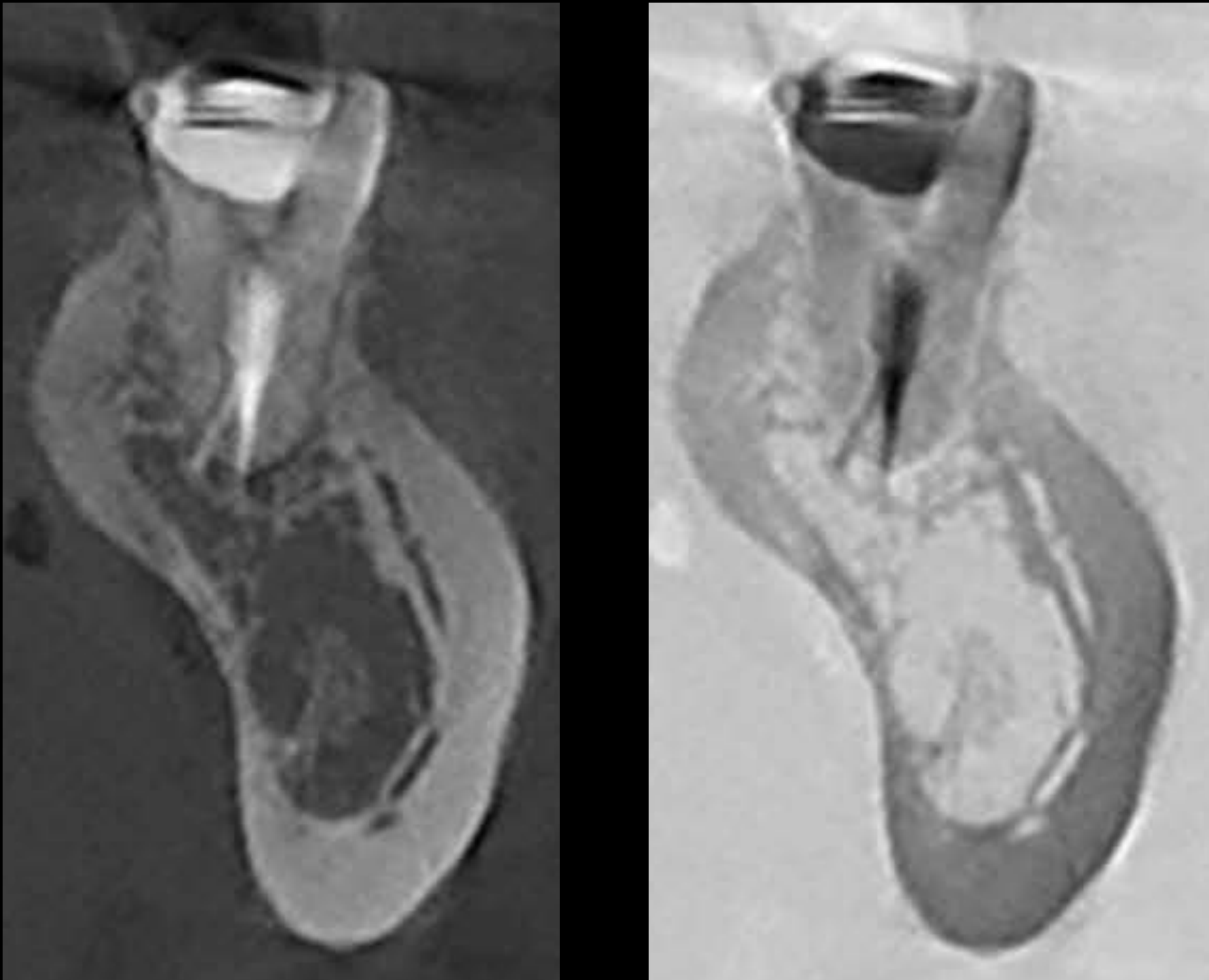



Follow up : 7 week (2024-03-27)

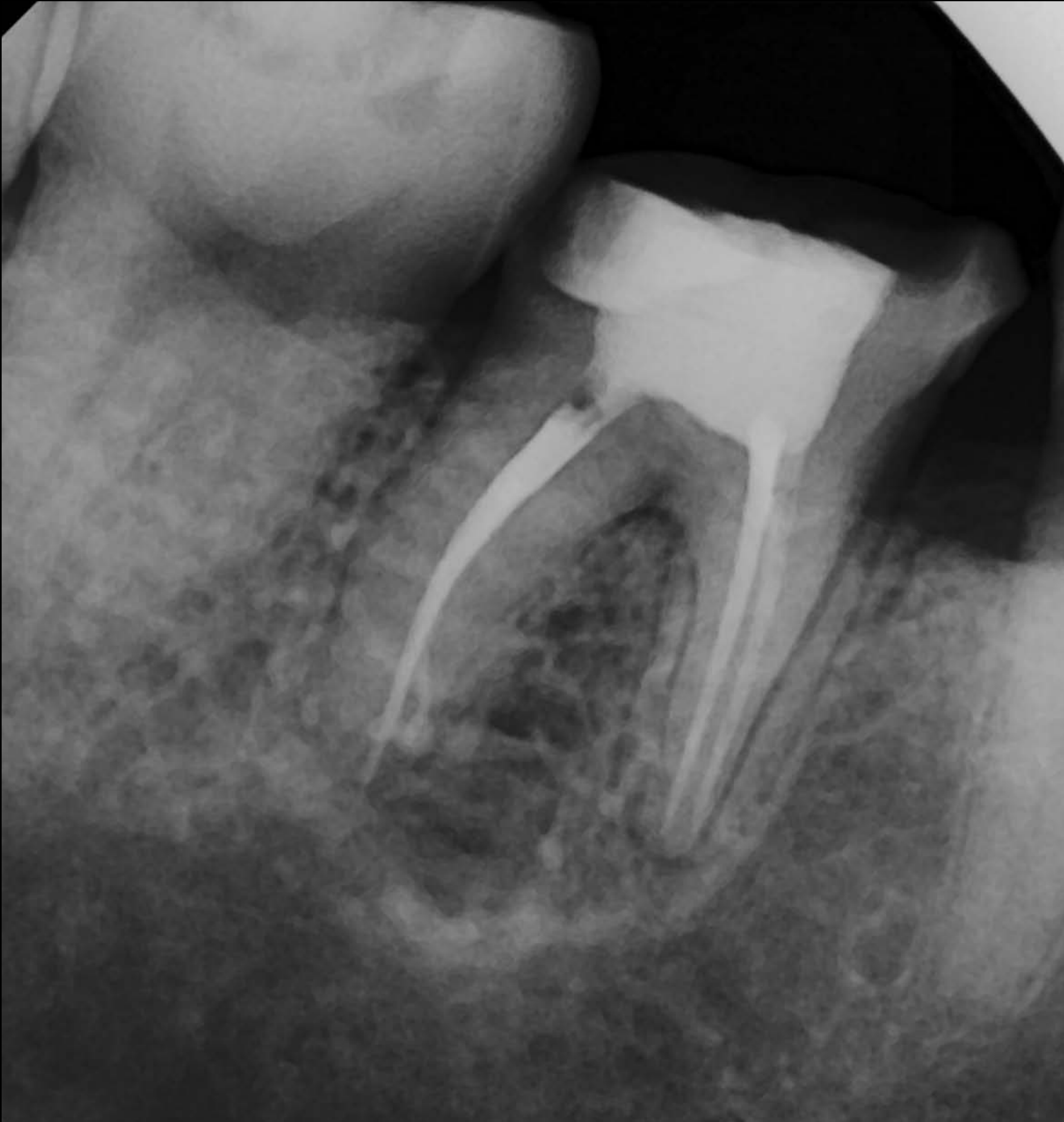



Main ingredient of MTA

CaO   SiO<sub>2</sub>



It rises to pH 12 and form Apatite -> Mineralization



- #46
- Natural Tooth
- Shade A3.5
- Shade Block







# Working model








# Working Model

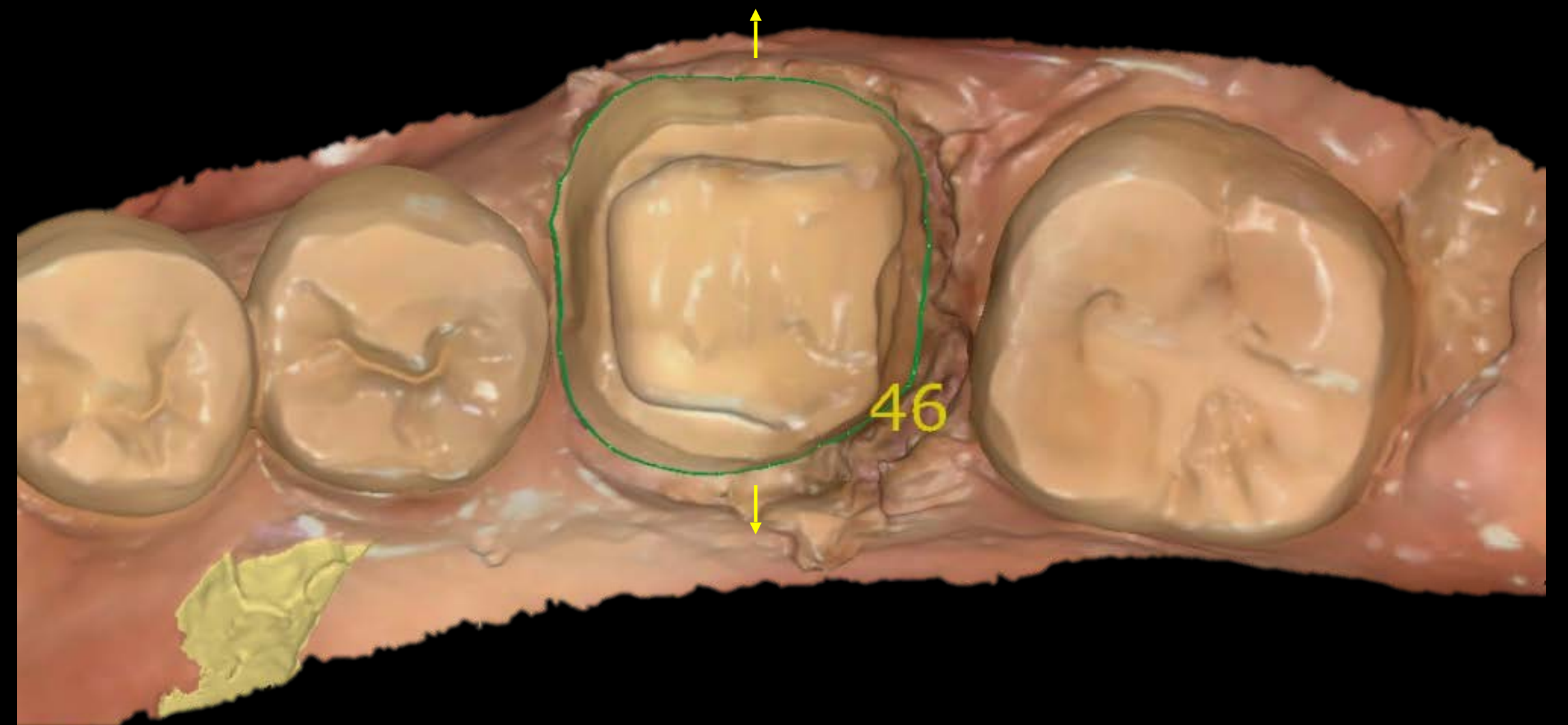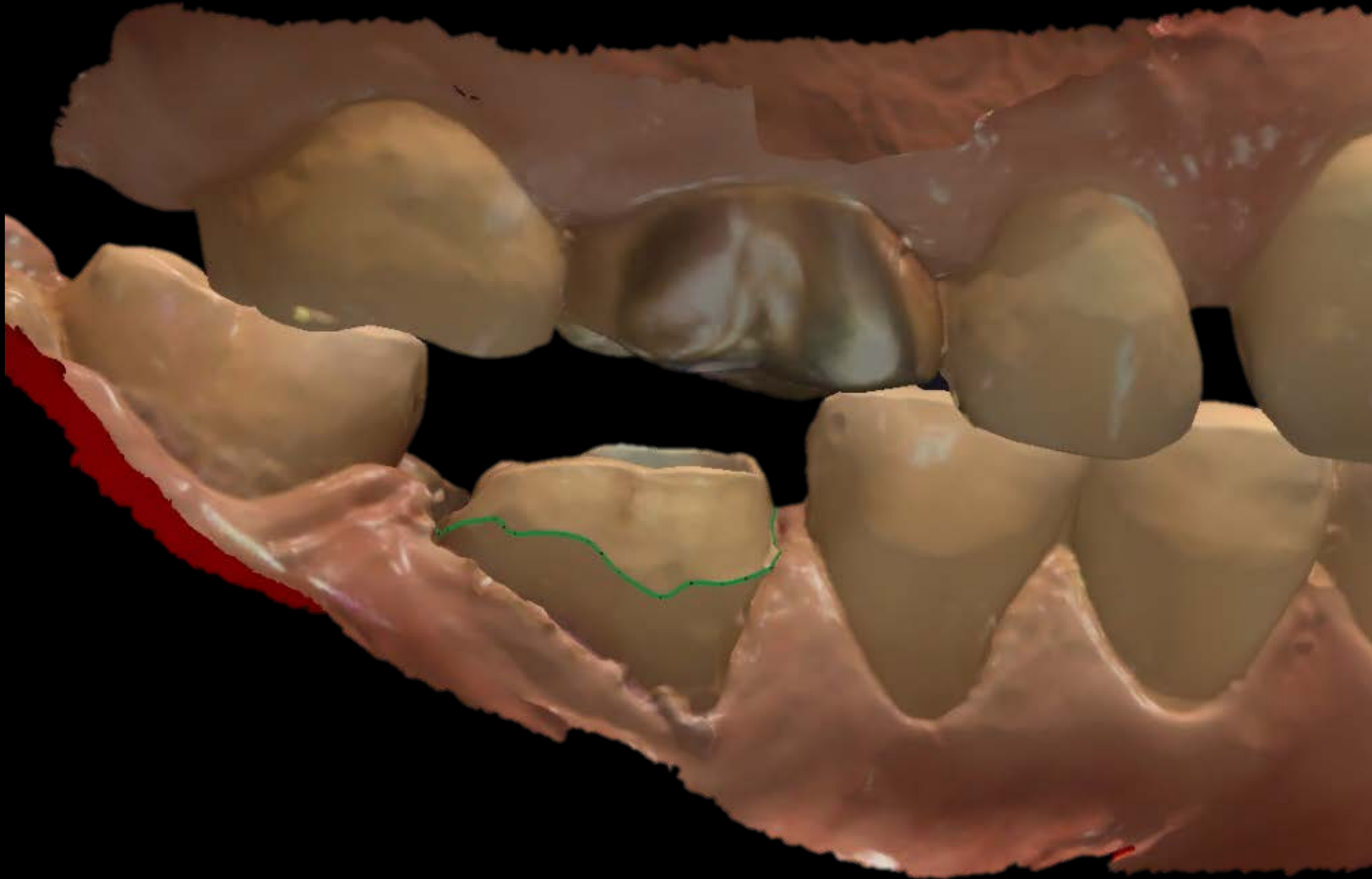



# Cad Design

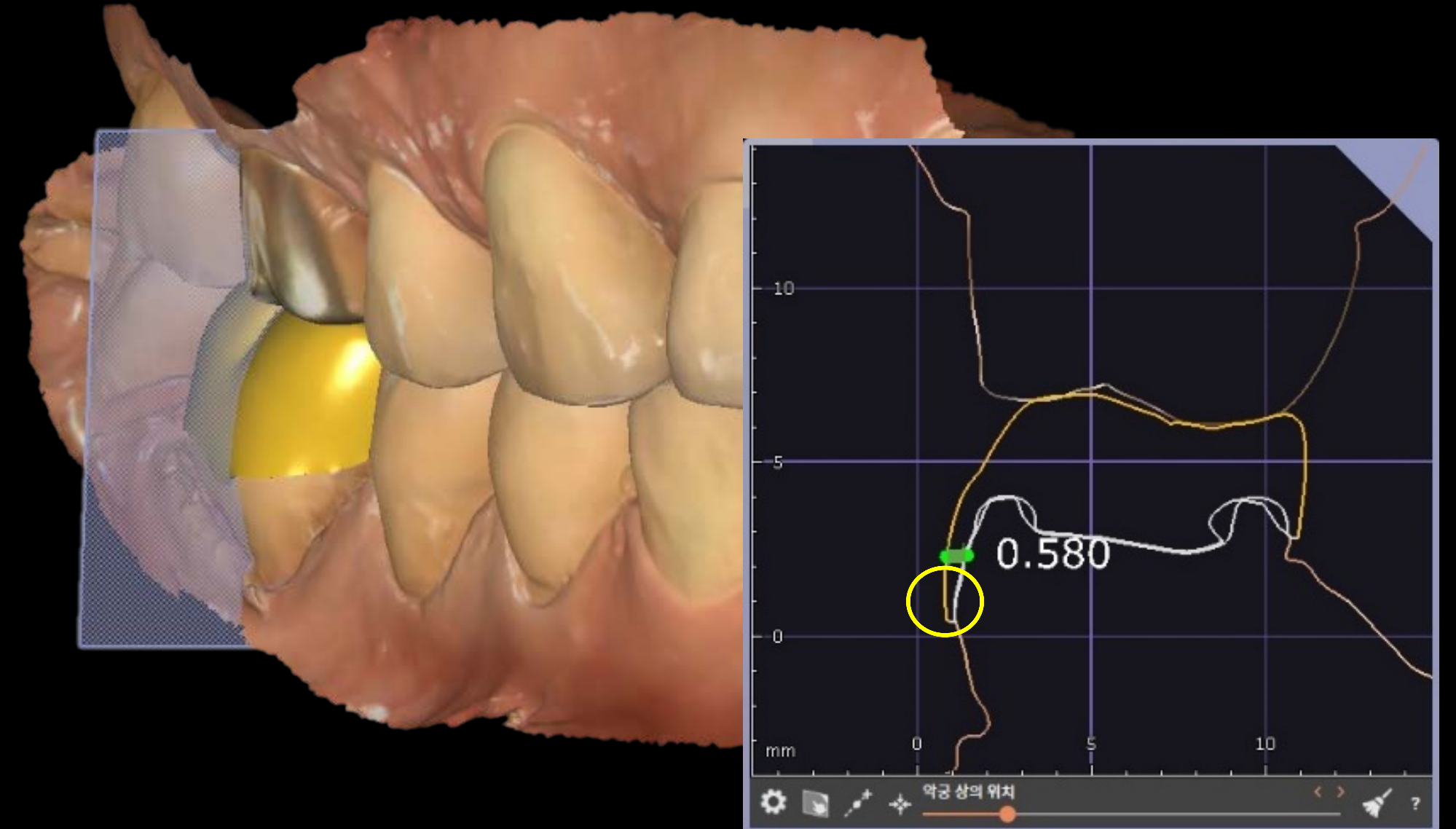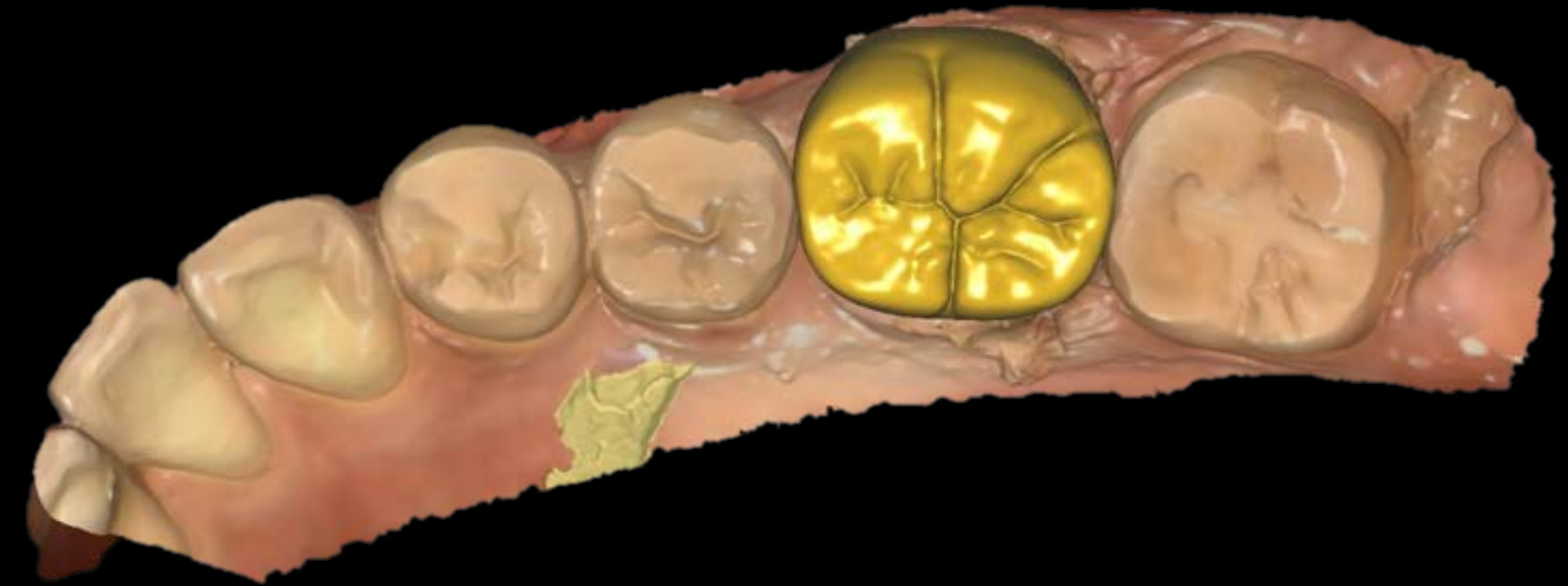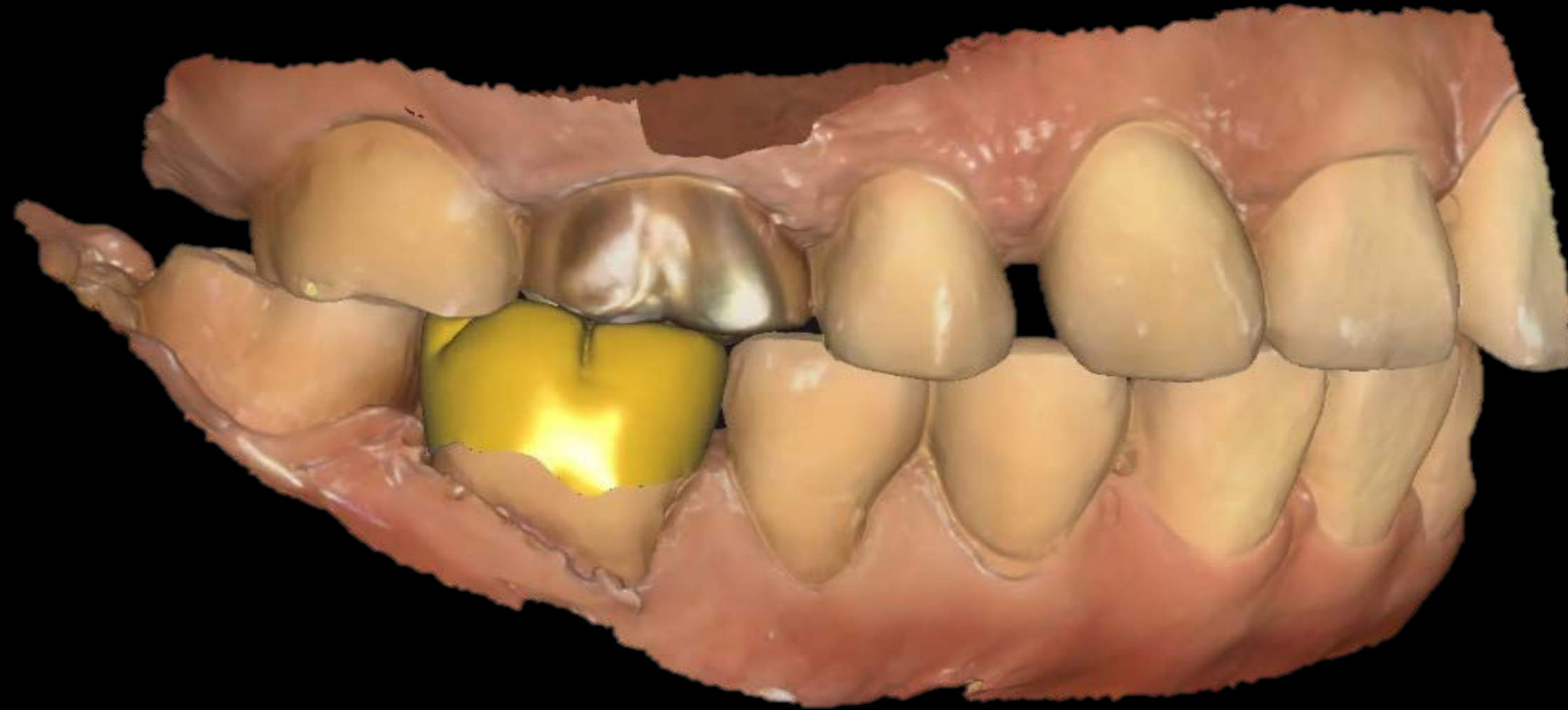




# Shade Taking





# Hyper Dent

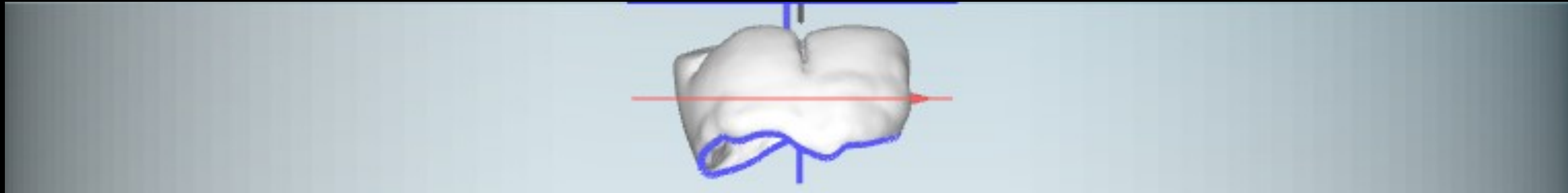

## Shade Block A2



# Coloring



Inner – White Opaque  
Cervical, Groove – A3.5  
Occlusal – Light Gray



# Glazing





# Final Prosthesis






# Bright Primer

- **Universal Primer containing 10-MDP**
- Improves Bond strength between cement and restorations
- Excellent bond strength to Zirconia & Metal & LS2




## Bond Strength


@ Zirconia

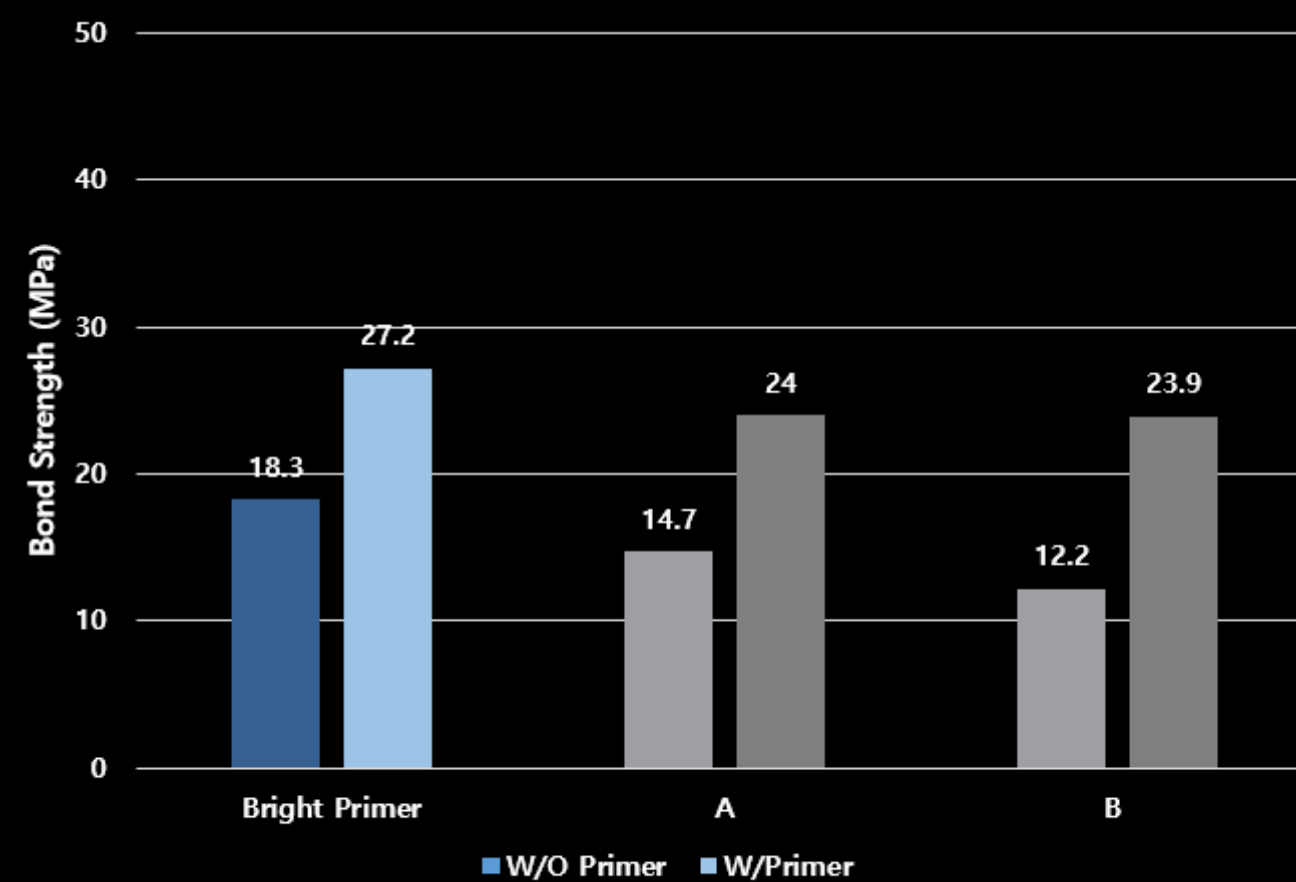

@ Metal

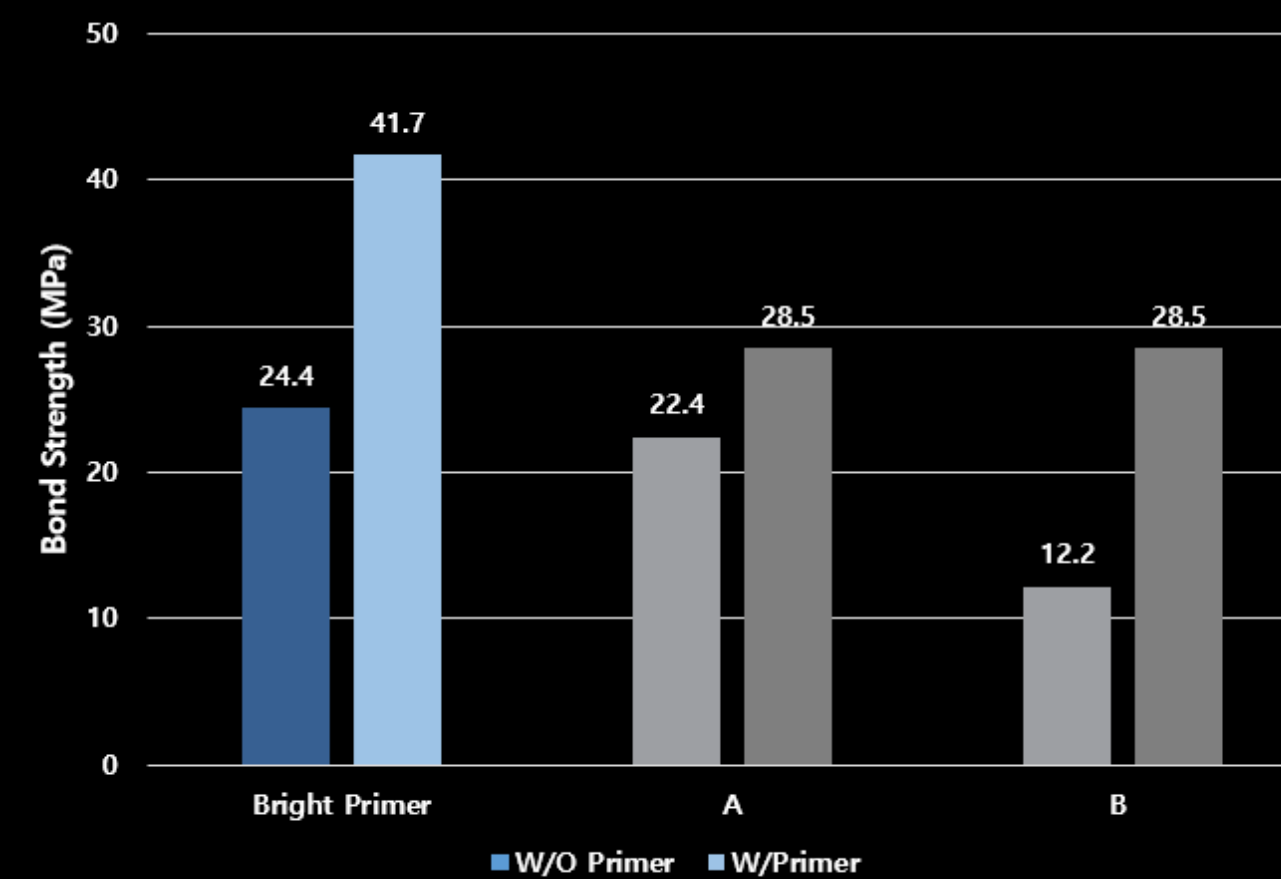

@ LS2

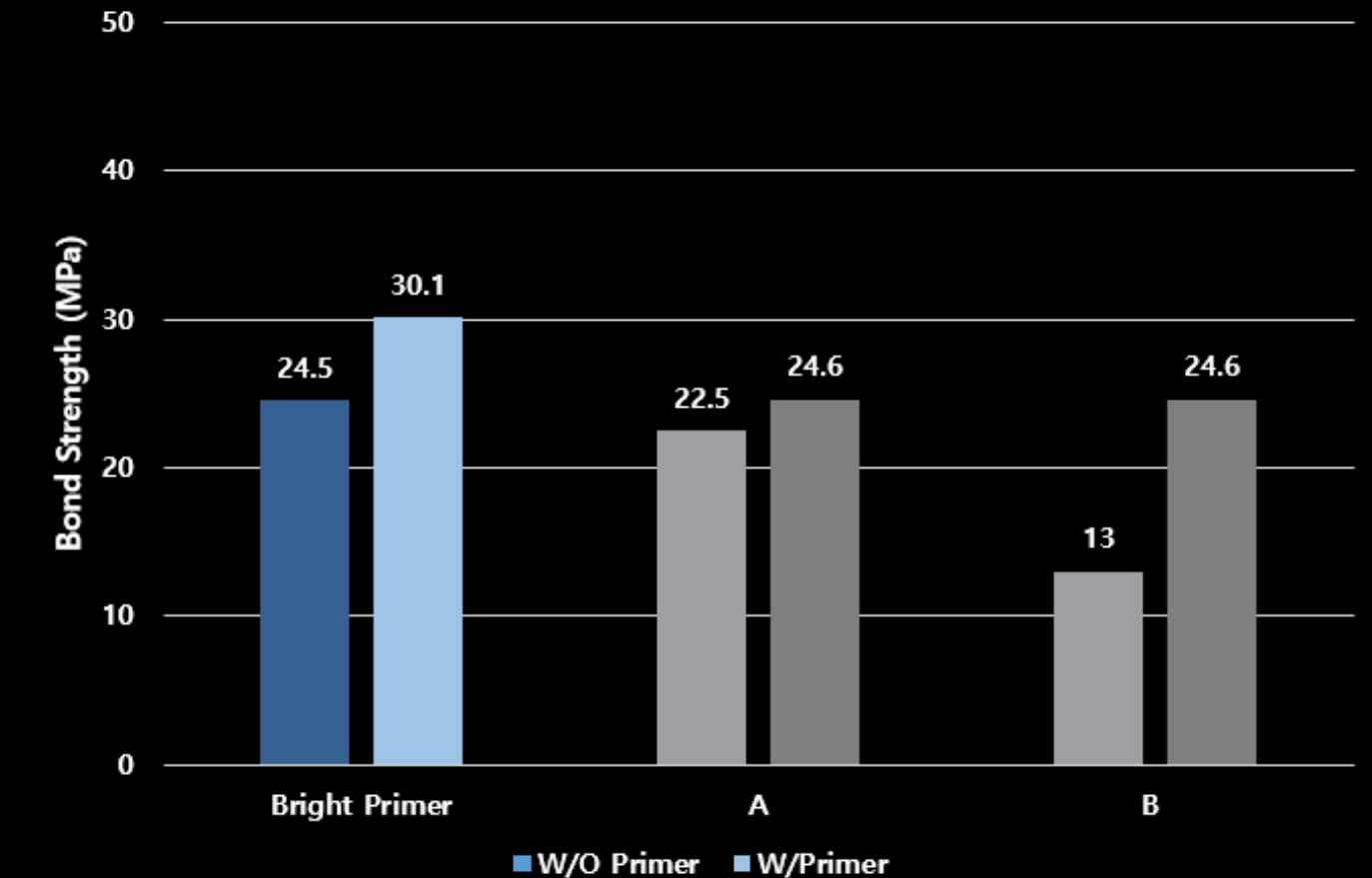




# Bright Resin Cement 2




- **Strong (10-MDP) → Strong bonding to indirect restoratives**
- Stable storage stability in both refrigerated and R.T storage
- Long Working time ( $\geq 2$  min) / Short Setting time ( $\leq 4$  min)



Bond Strength

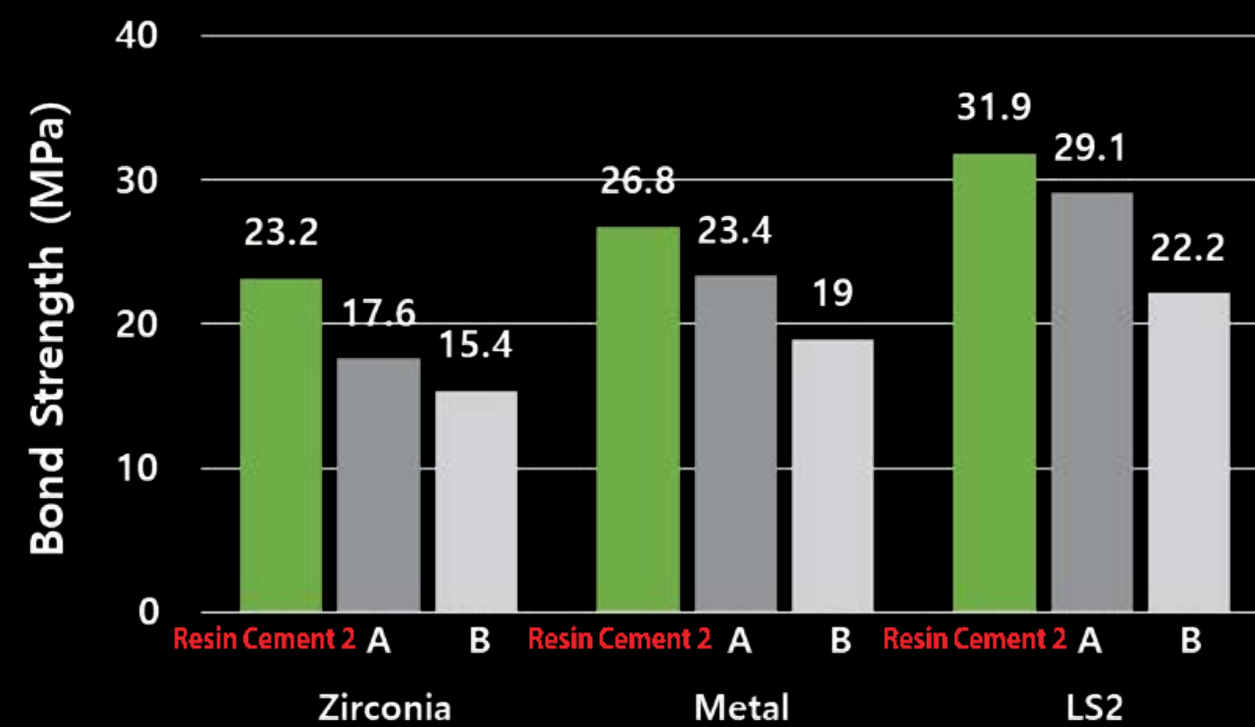


Adhesion durability (72h F/U)

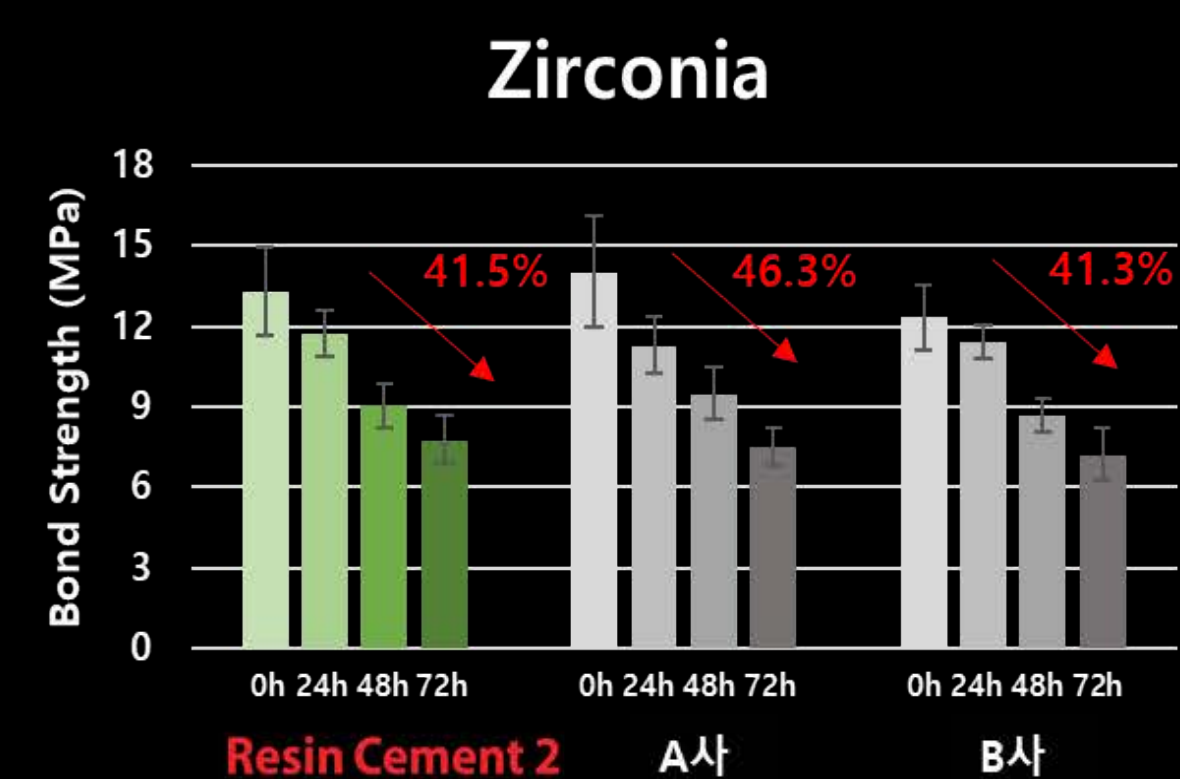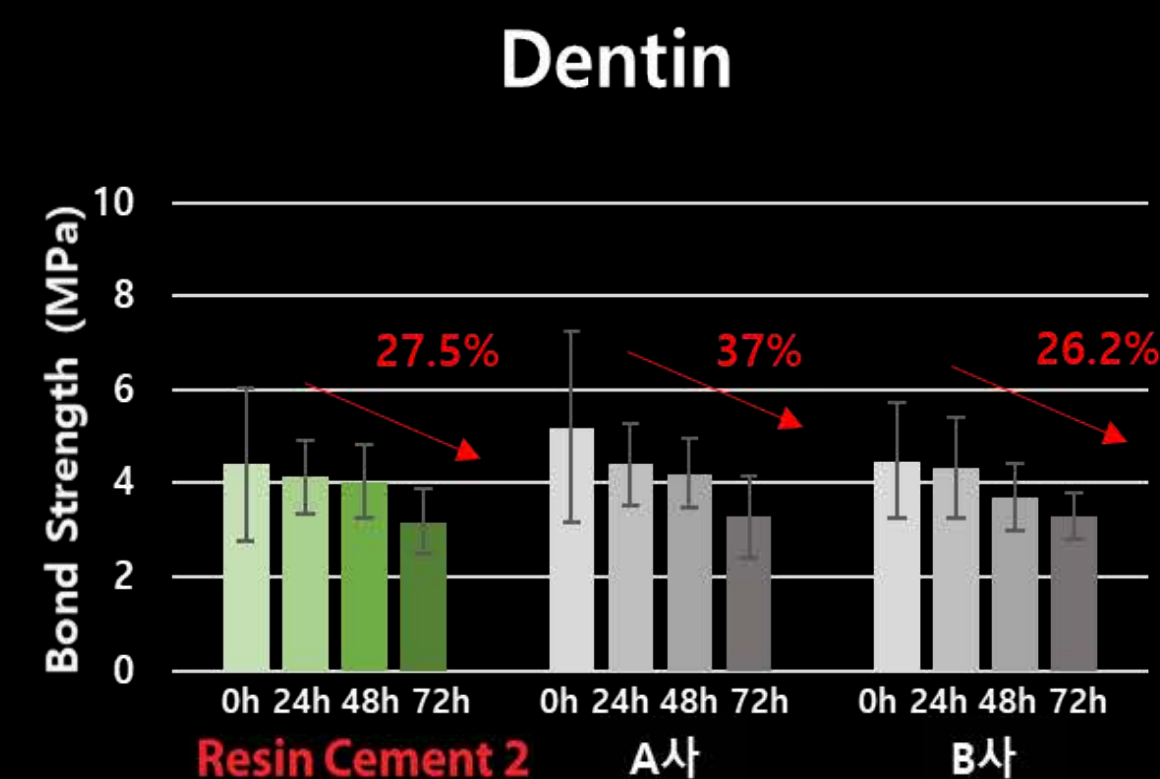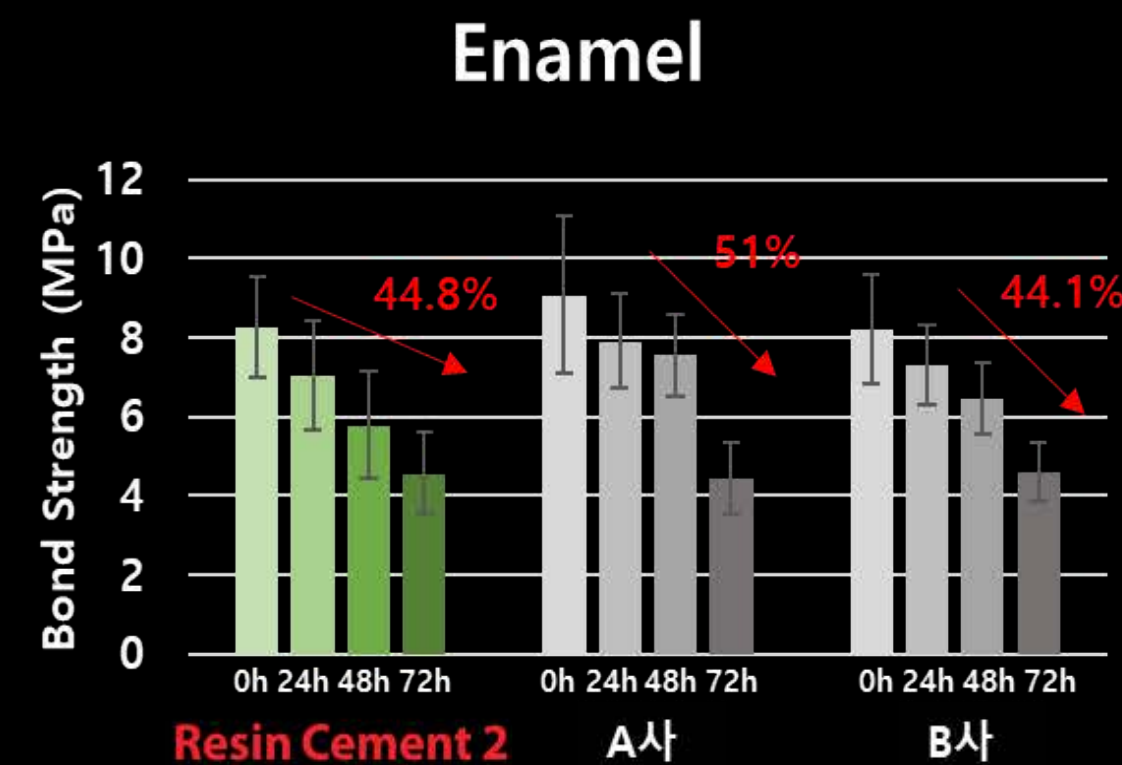




Oral exam (2021-09-16)

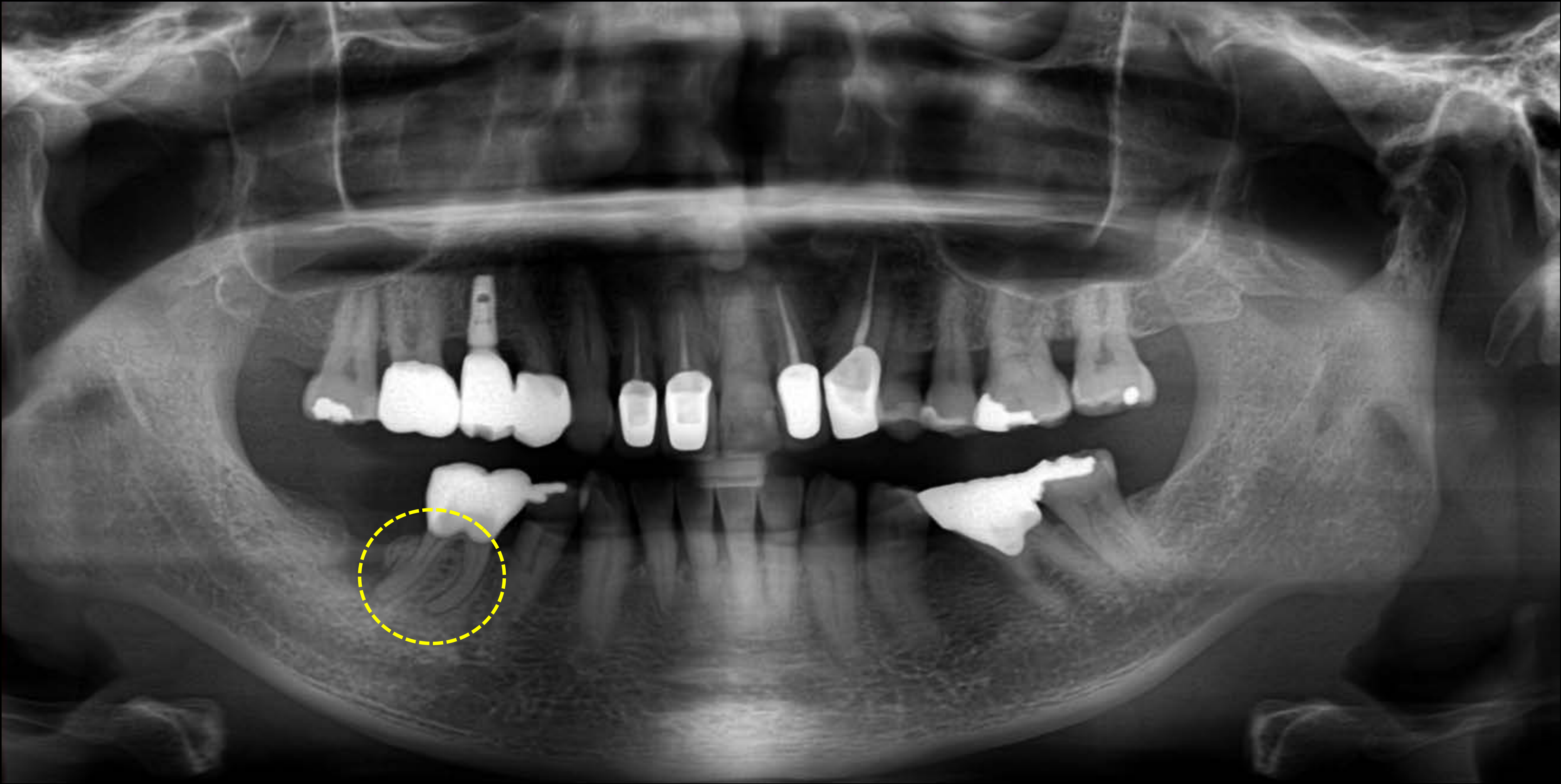




Oral exam (2022-03-10)

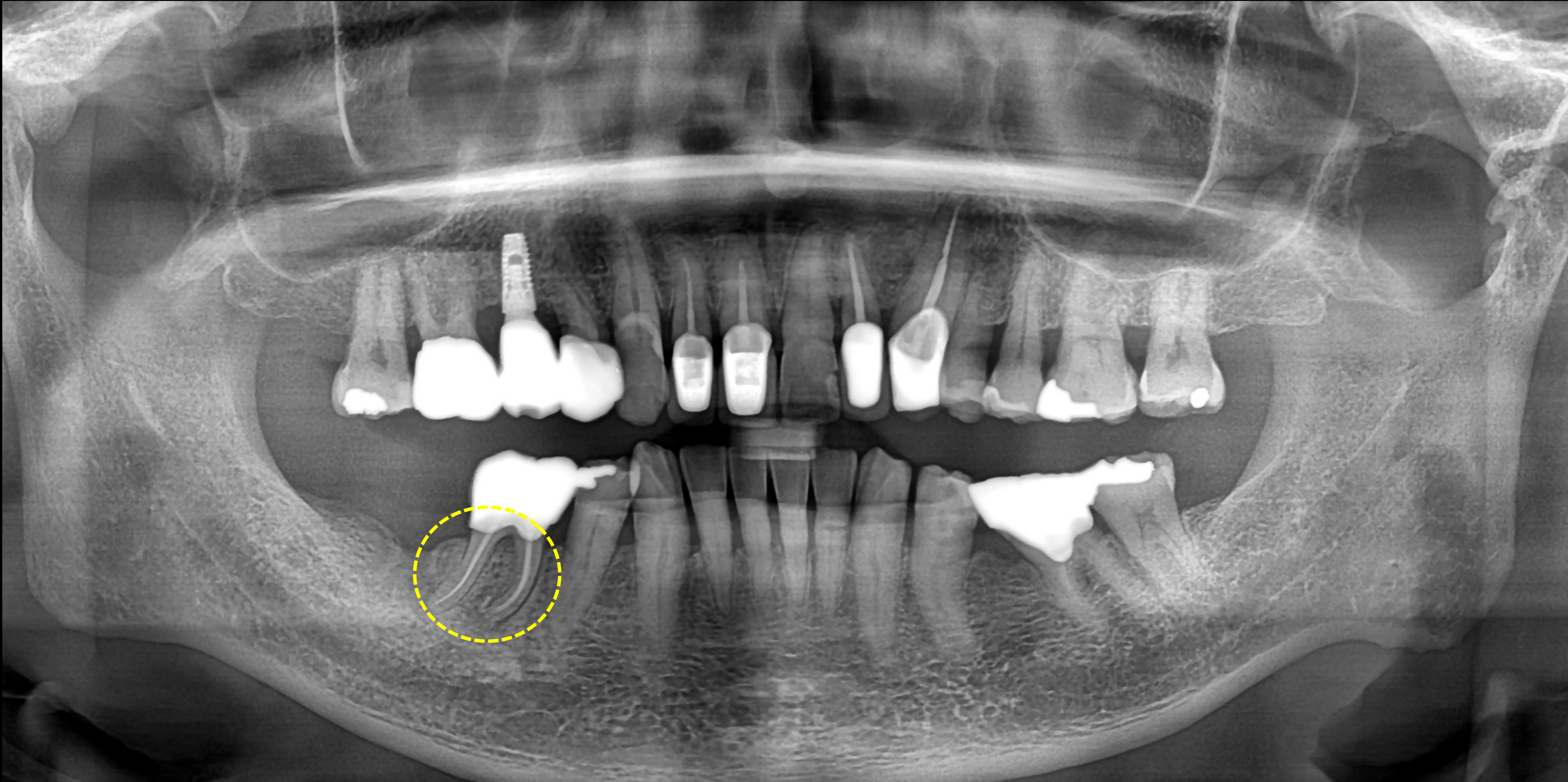




Oral exam (2022-11-30)

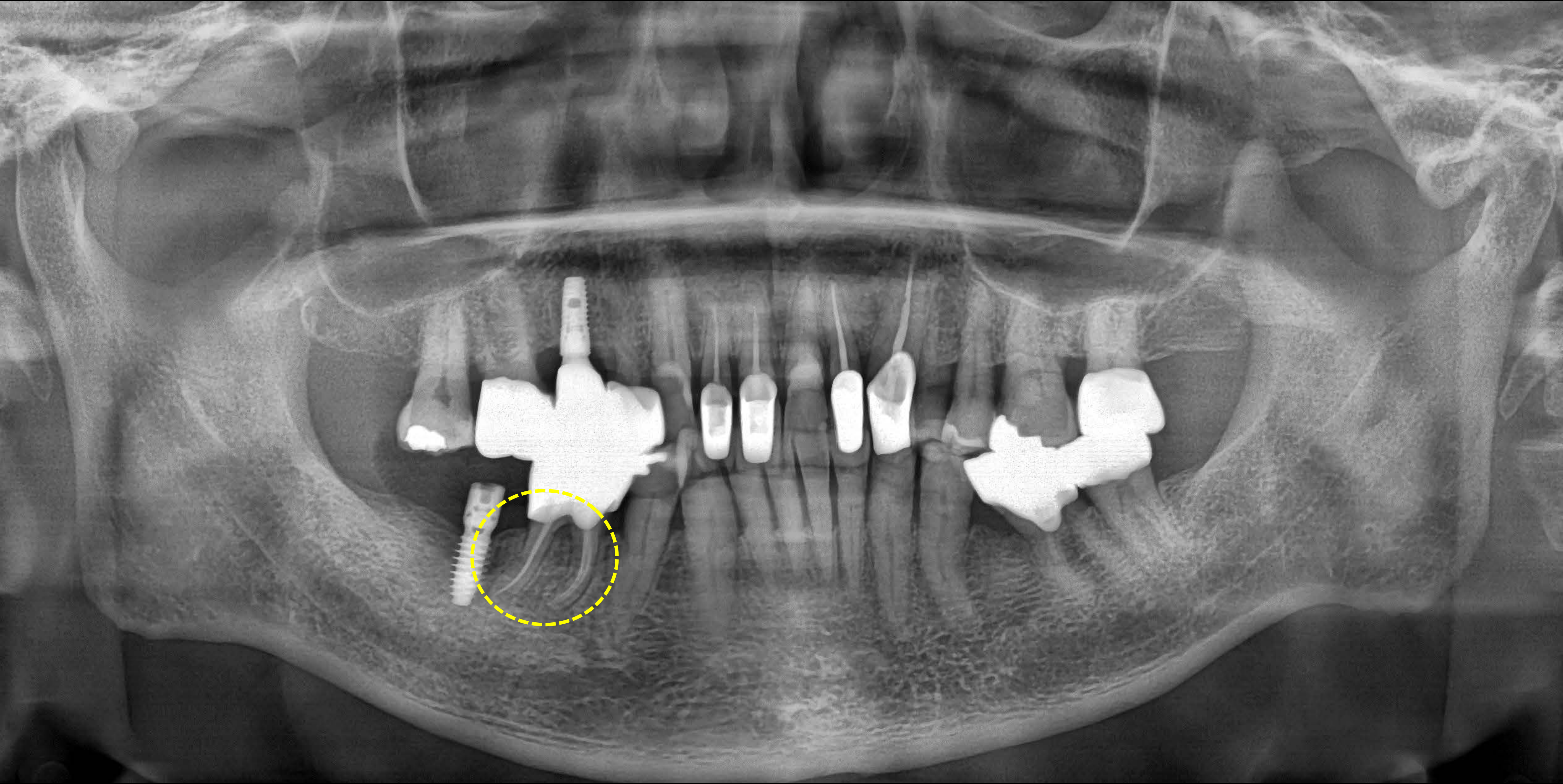





Oral exam (2023-03-22)

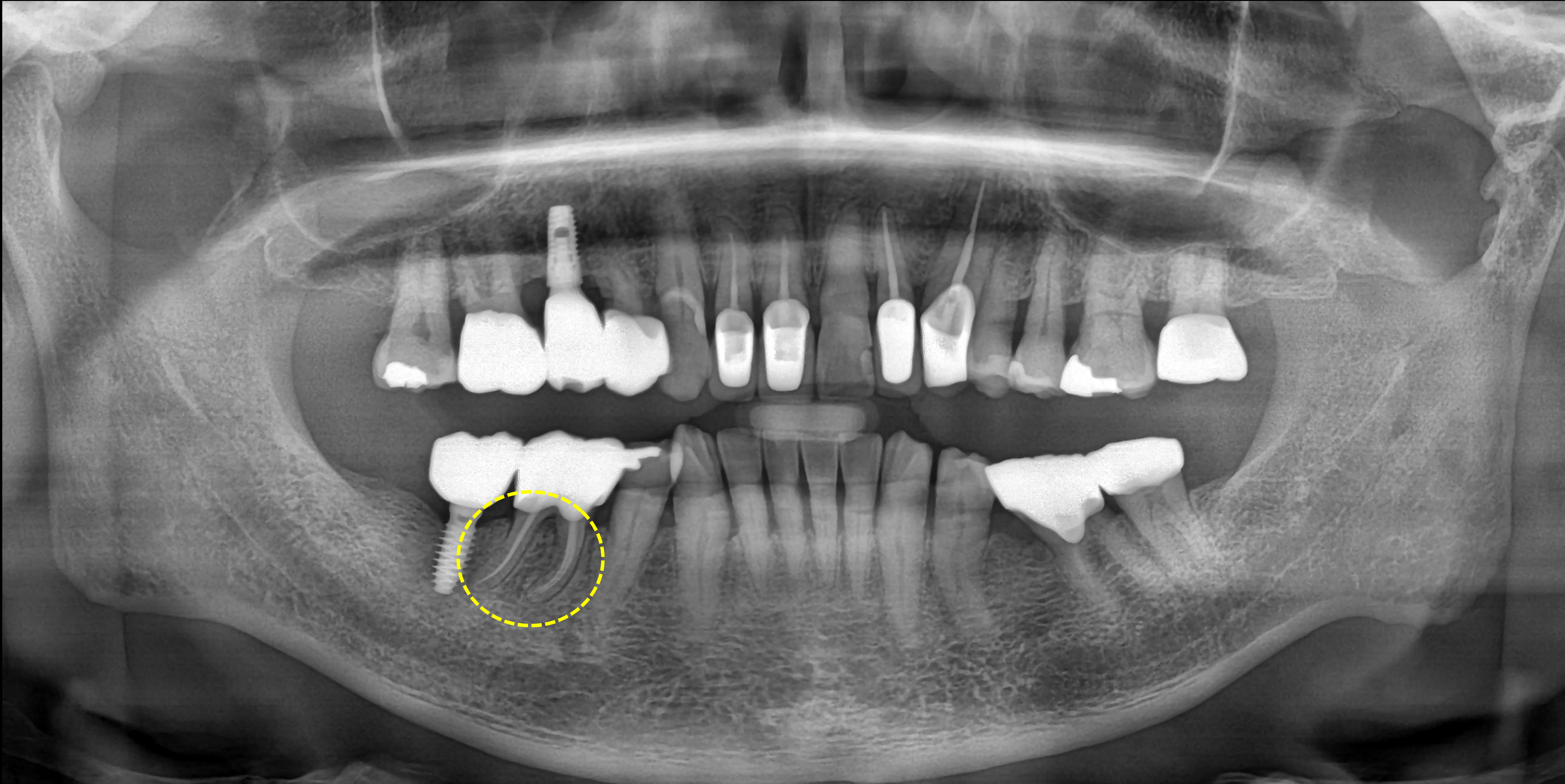




Pre-op (2021-09-16)

Follow-up : 1 year 4 month (2023-03-22)

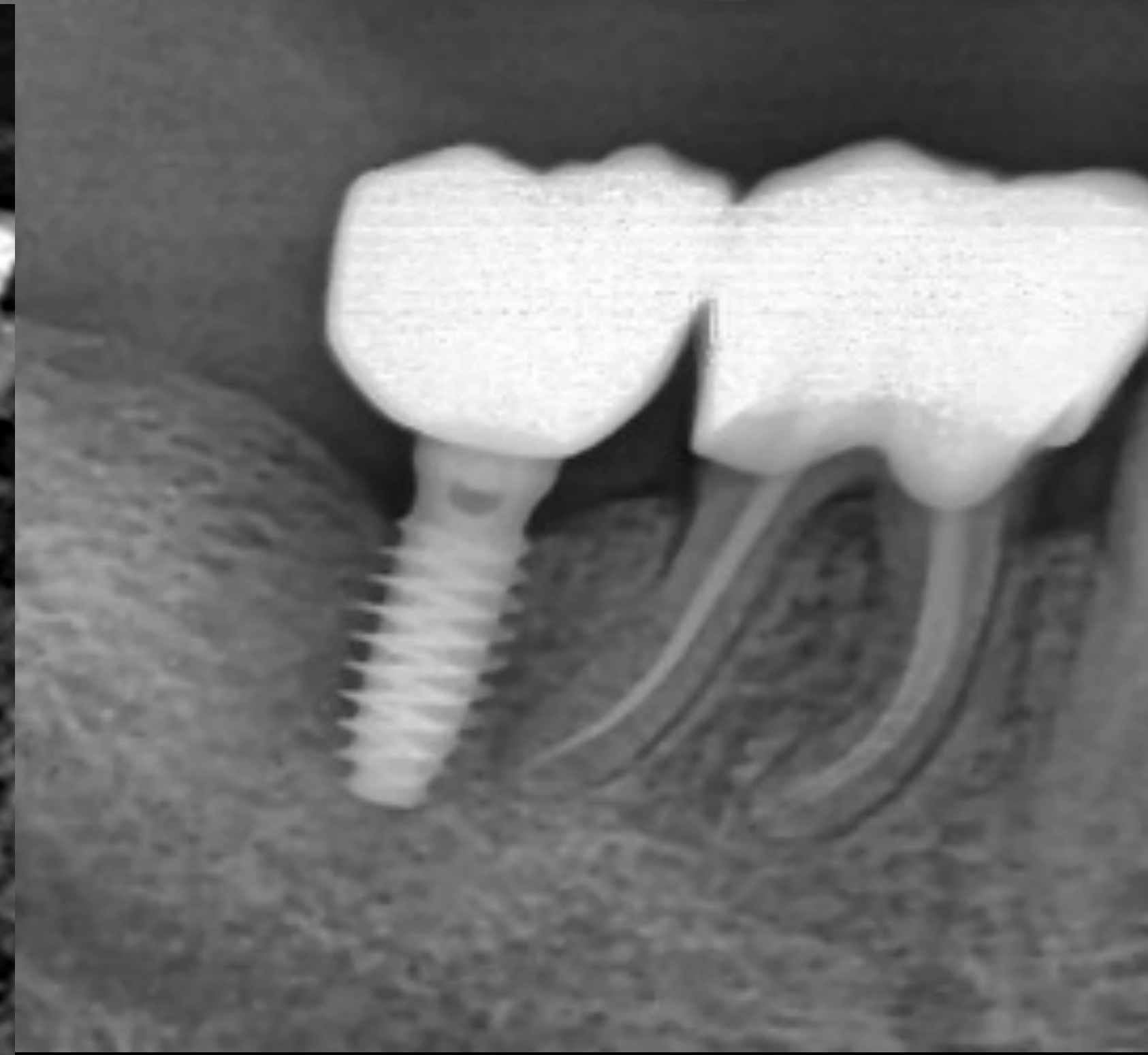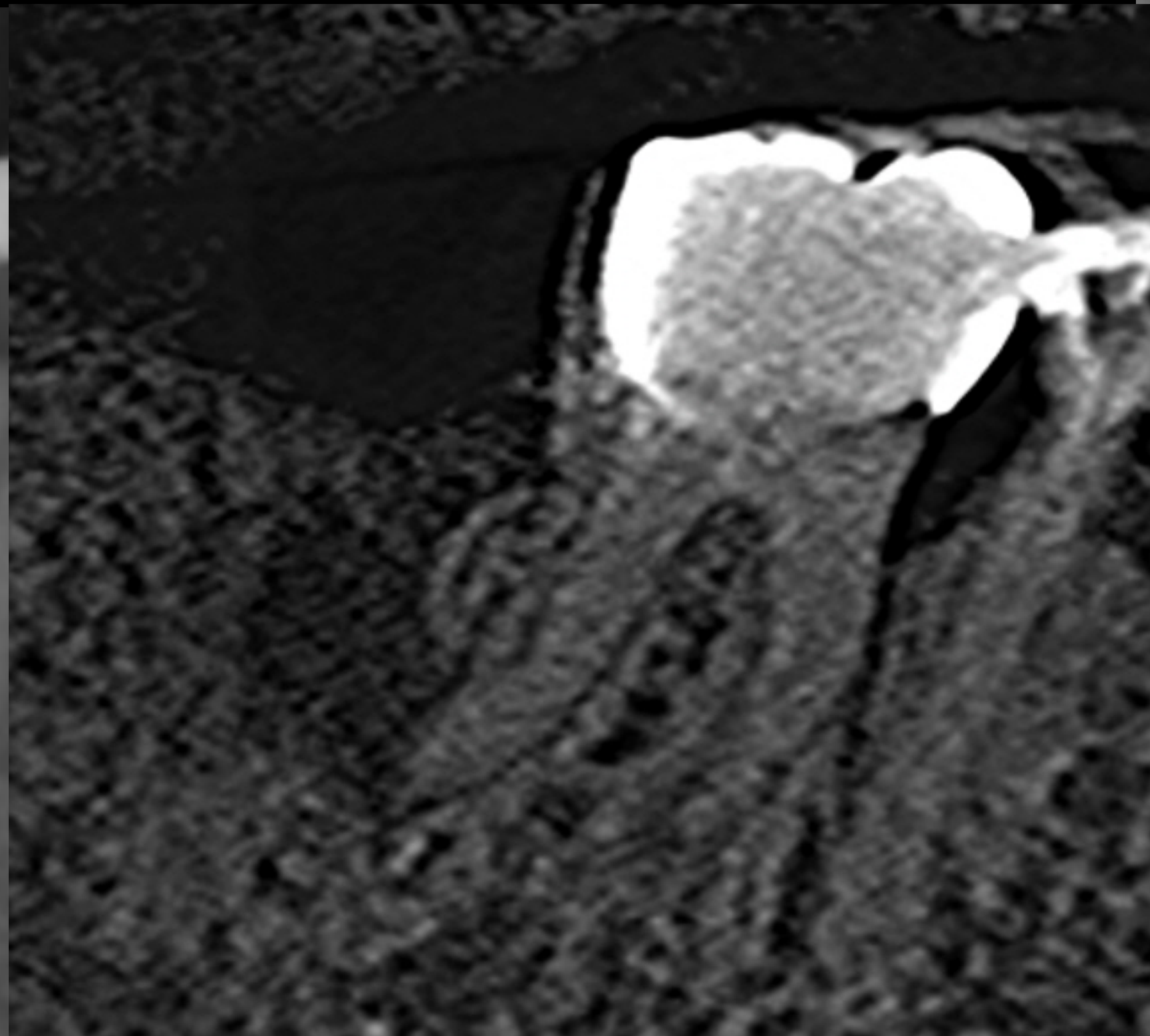




Follow-up : 4 month (2022-03-10)

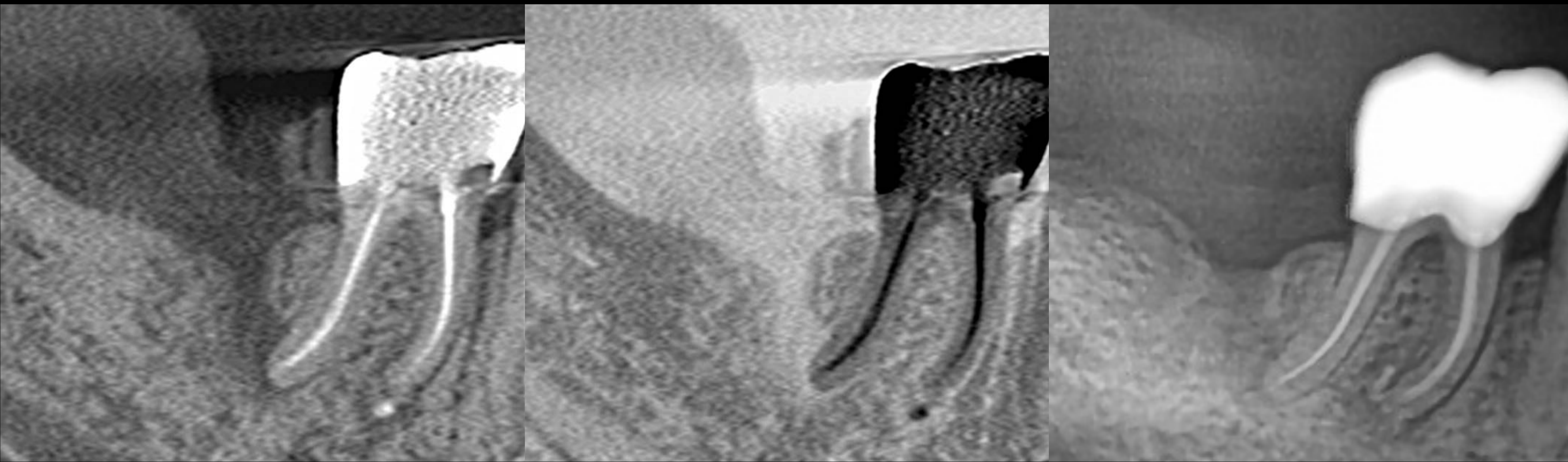





Follow-up : 1 year (2022-11-30)

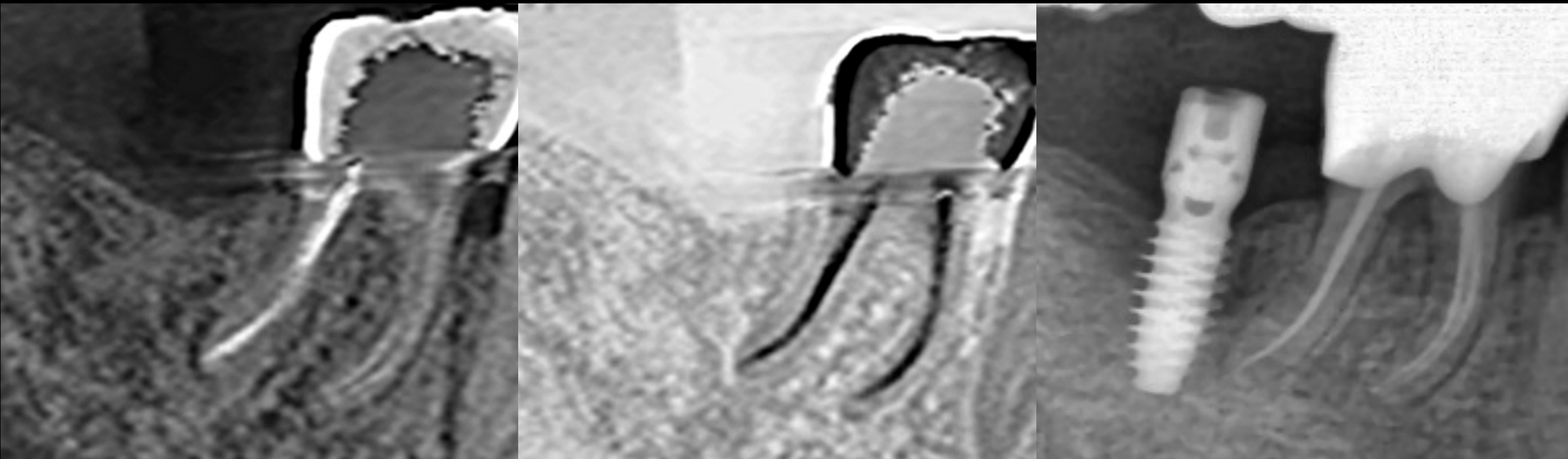




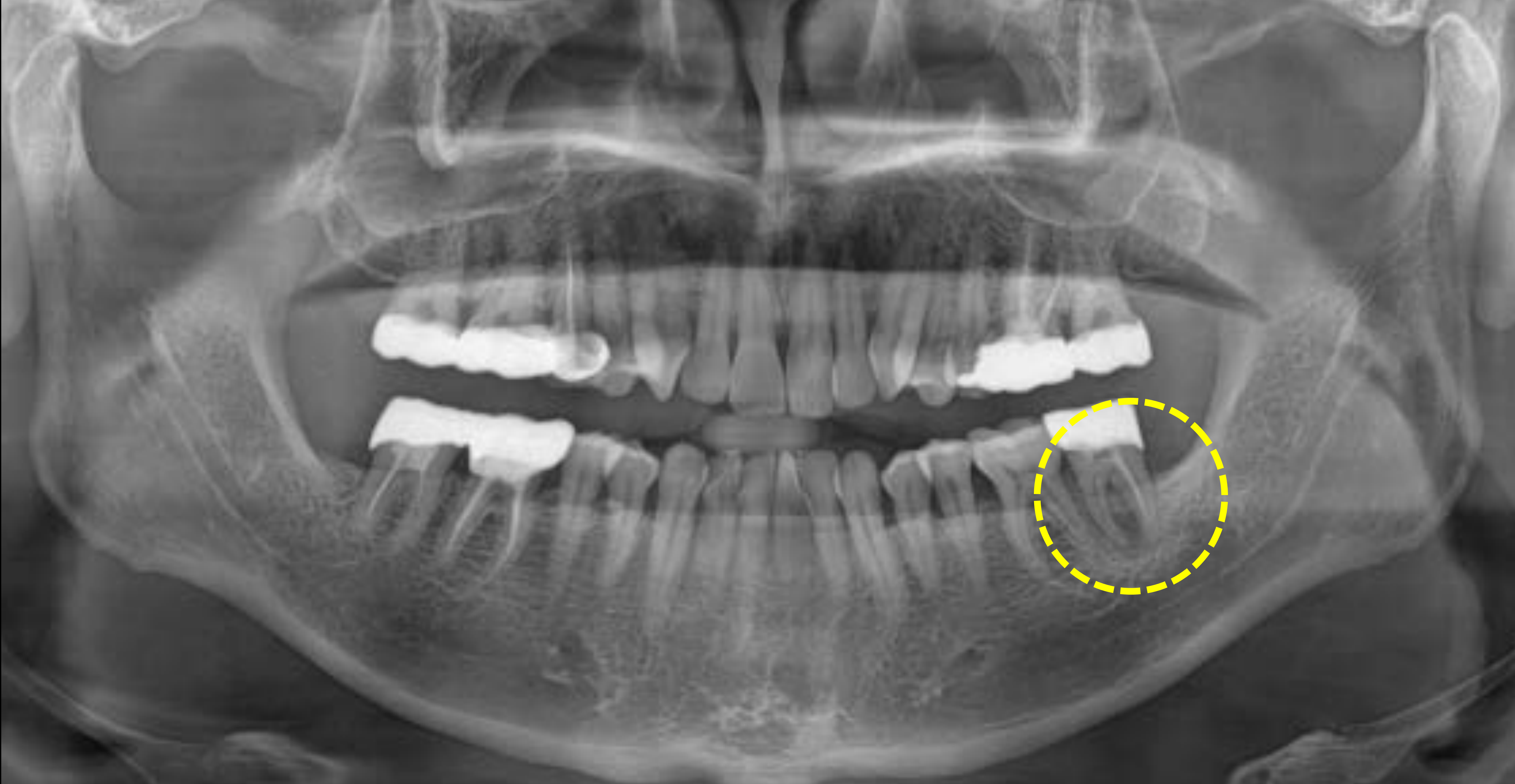



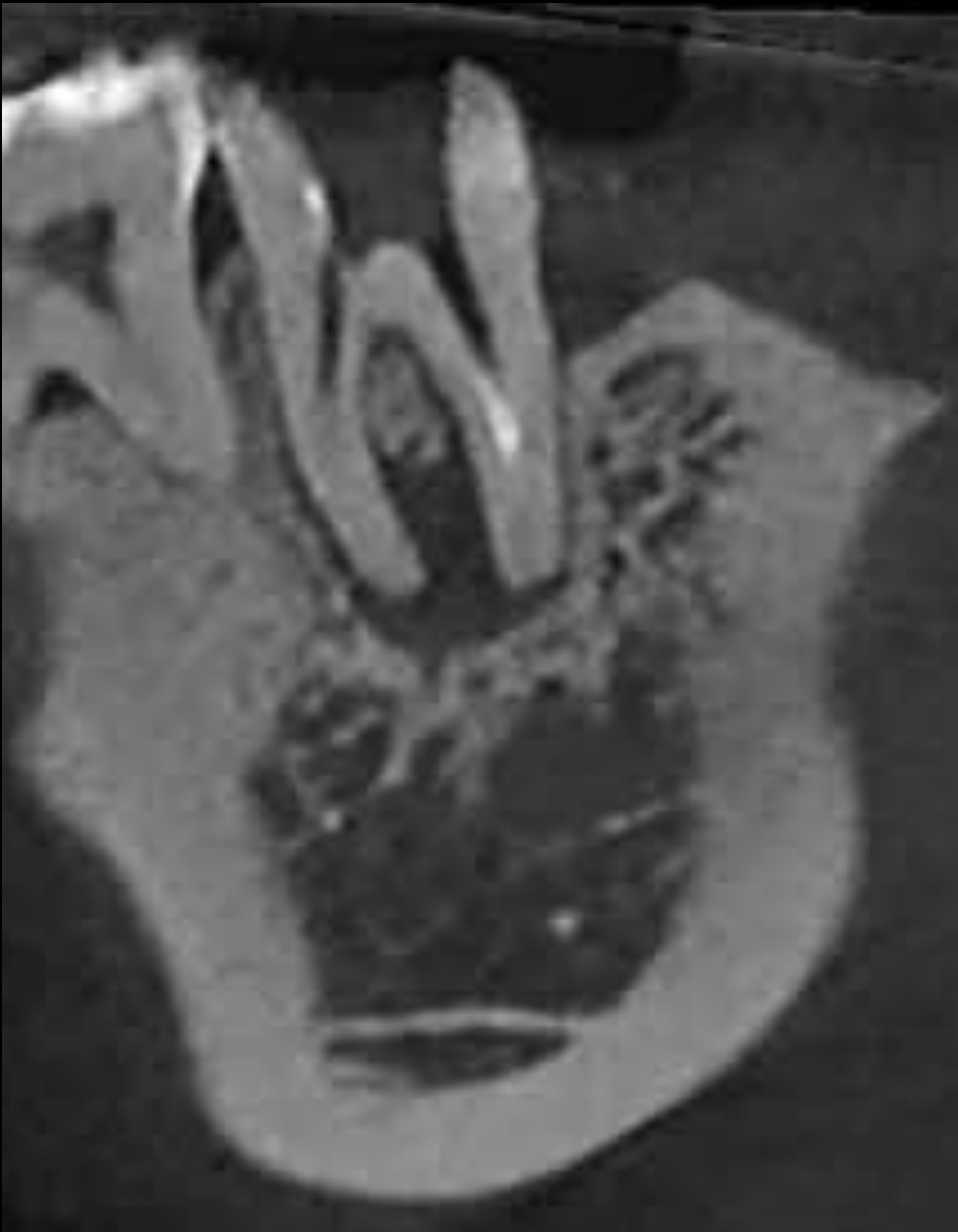




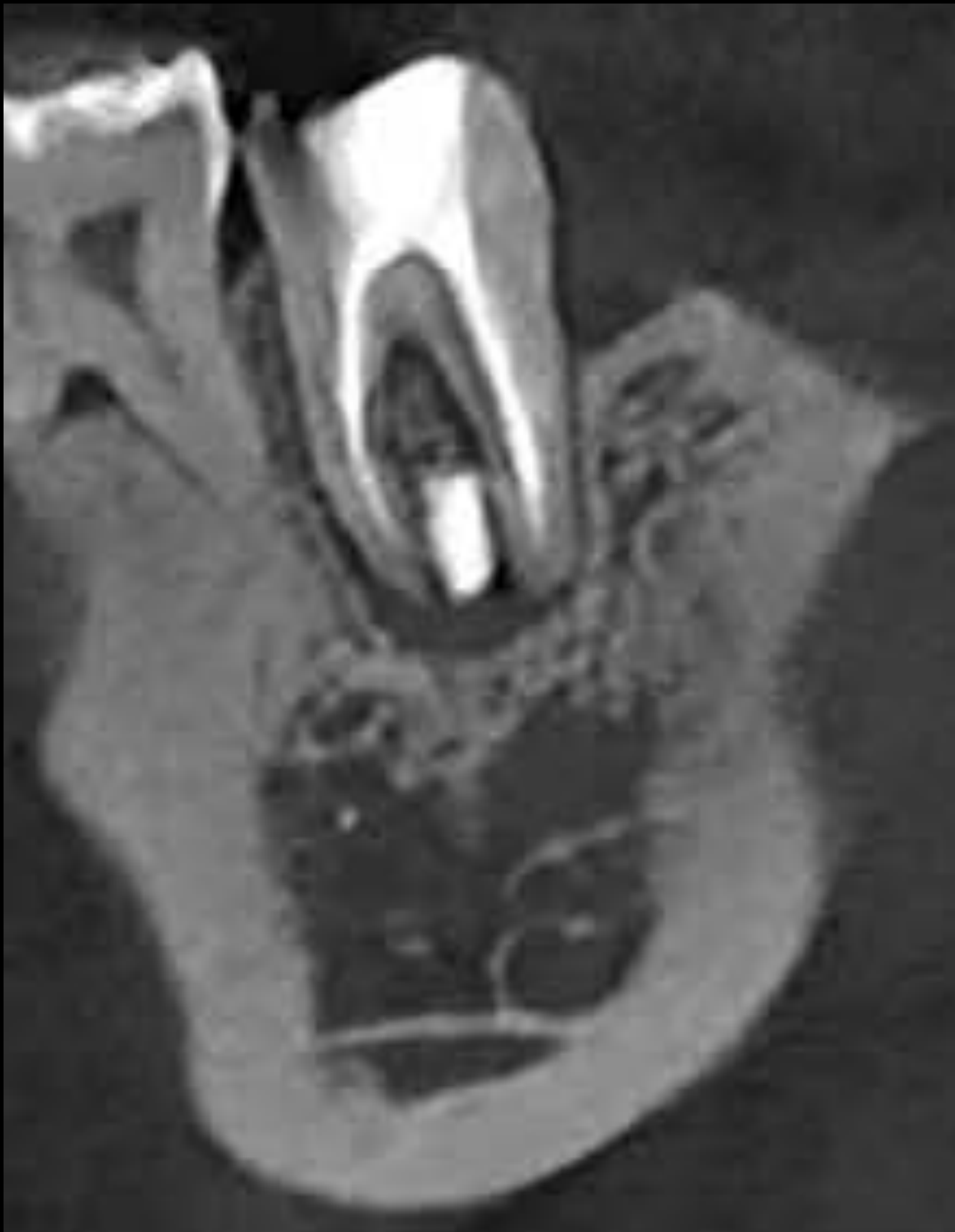




### #37 Mesial canal MTA Filling (2024-02-27)

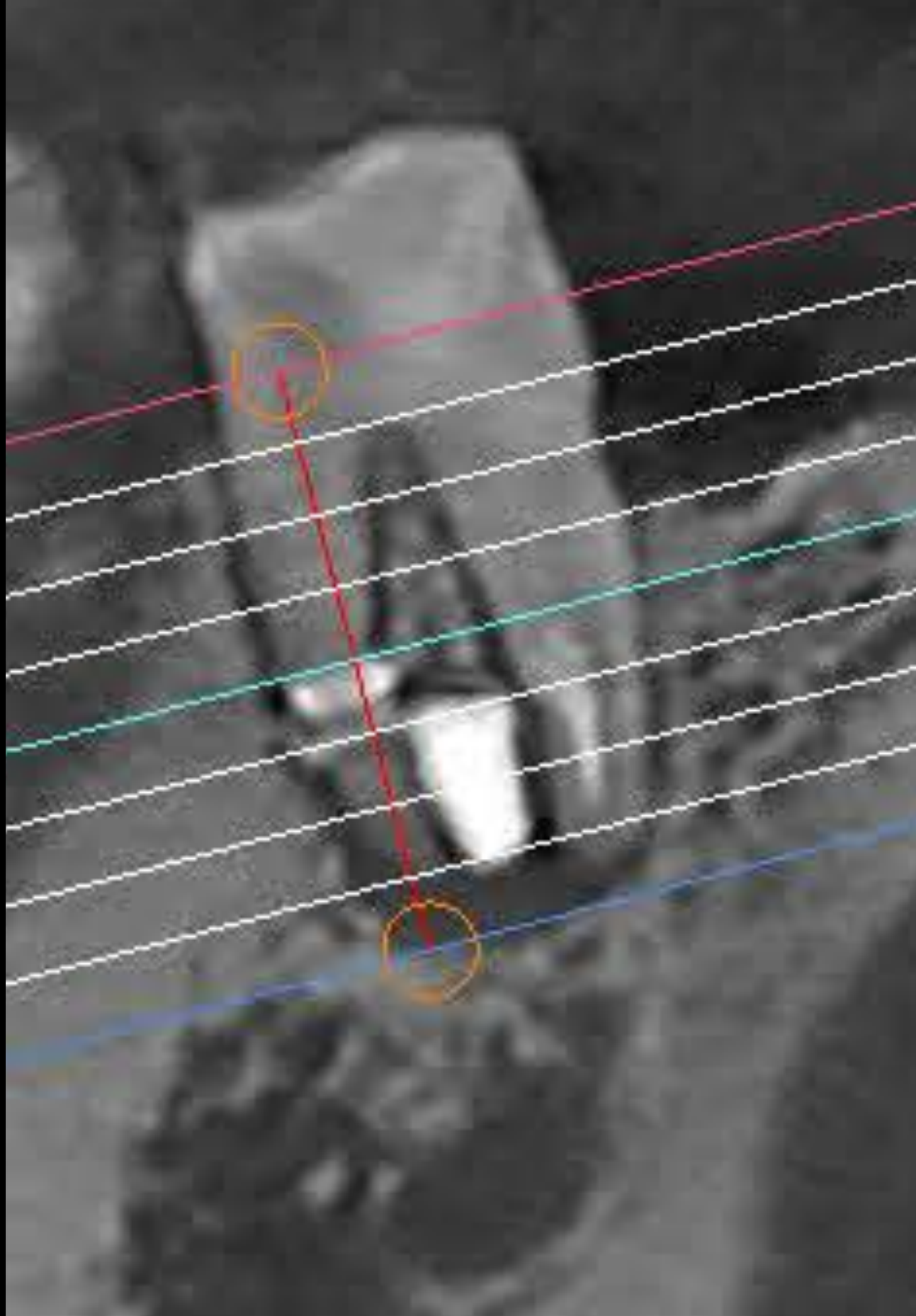




# #37 Digital canal GP-Cone + MTA Sealer filling (2024-02-27)

< Pre-op CT >

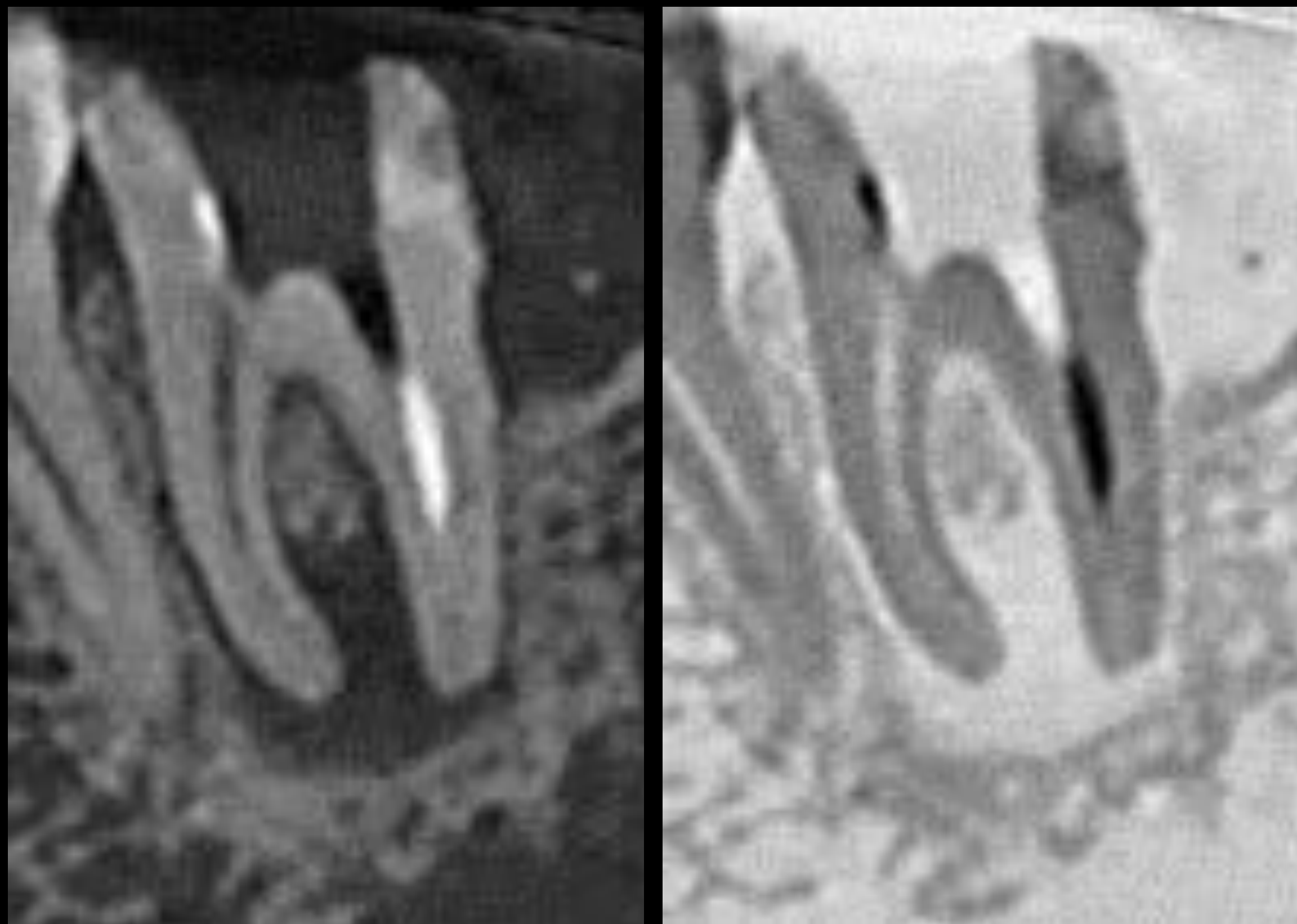

< Post-op CT >

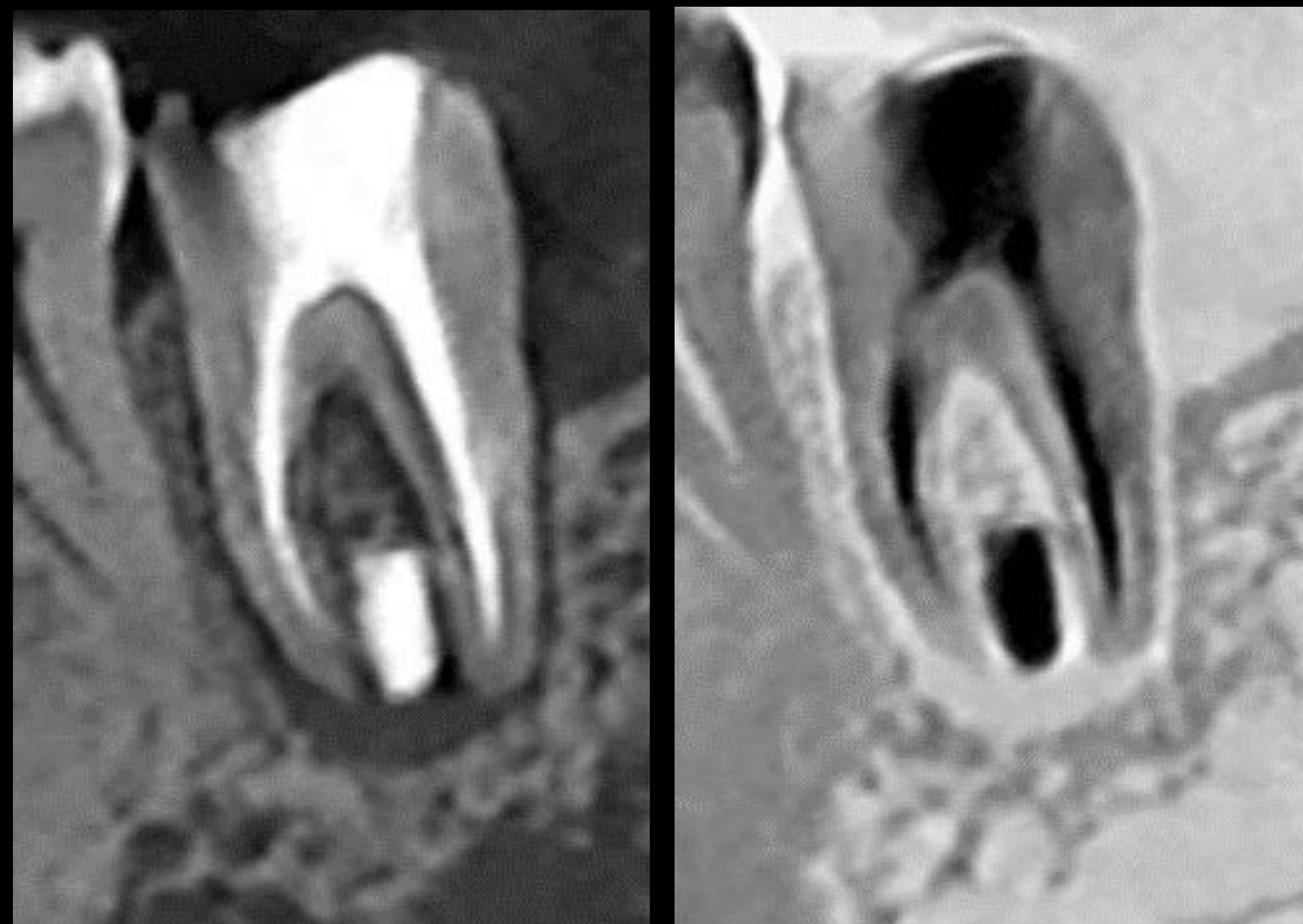



Final Periapical radiograph (2024-02-27)




Follow up : 1 month (2024-03-27)





Oral exam (2023-10-20)

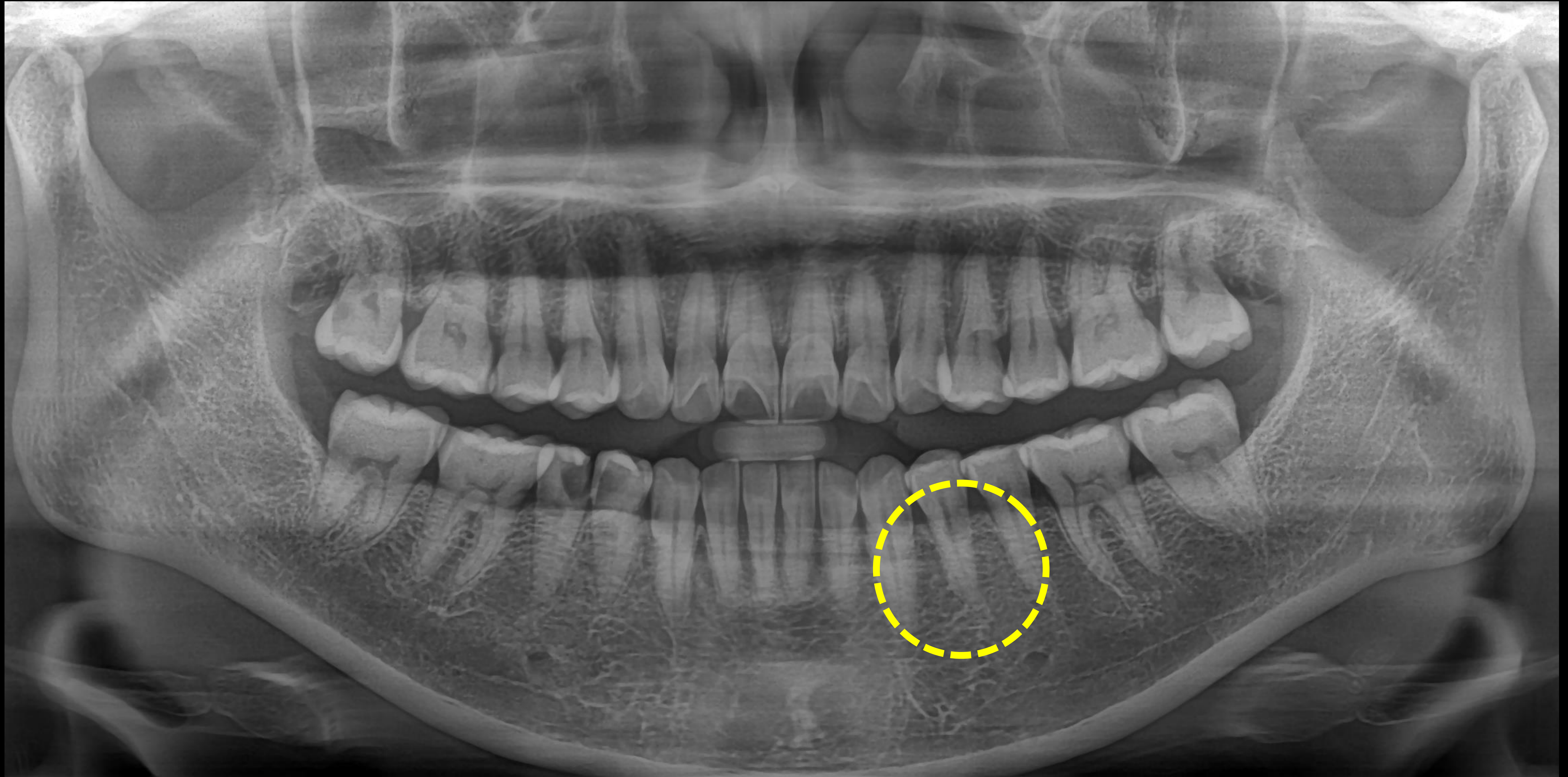



Pre-op CT (2024-03-19)



Post-op (2024-03-19)





Resin core & Abrasion building (2024-03-19)



Bright Universal Etchant







- Bright Universal Etchant
- Bight Primer
- Bright Bond Universal

- Bright High Flow
- Bright Low Flow



# Bright Flow

Light-cured Flowable Composite Resin

**High flexural strength  $\geq 140\text{MPa}$**



Excellent radiopaque

Excellent bonding strength with Bond Universal



Company A

Low flow

Company B





**Bluemoon**  
Curing Light

Dual wave LED, Glass Lens, 1 Button, All 4 modes

**3 sec**

Incomparable to  
anything

**Dentium**

판매원

서울특별시 강남구 자곡동 648, 821호 (강남메이스타워)  
Tel. 02-555-3750 / [www.dentium.co.kr](http://www.dentium.co.kr)

**GENOSS**

제조원


경기도 수원시 영통구 청량대로 256번길 56  
[www.genoss.com](http://www.genoss.com)

**Dentium**



# Compare Polymerization Depth performance with Competitors

- Bluemoon(4.1mm) ≐ Bluephase(4.1mm) > Valo(3.7mm) > Prime (3.3mm)

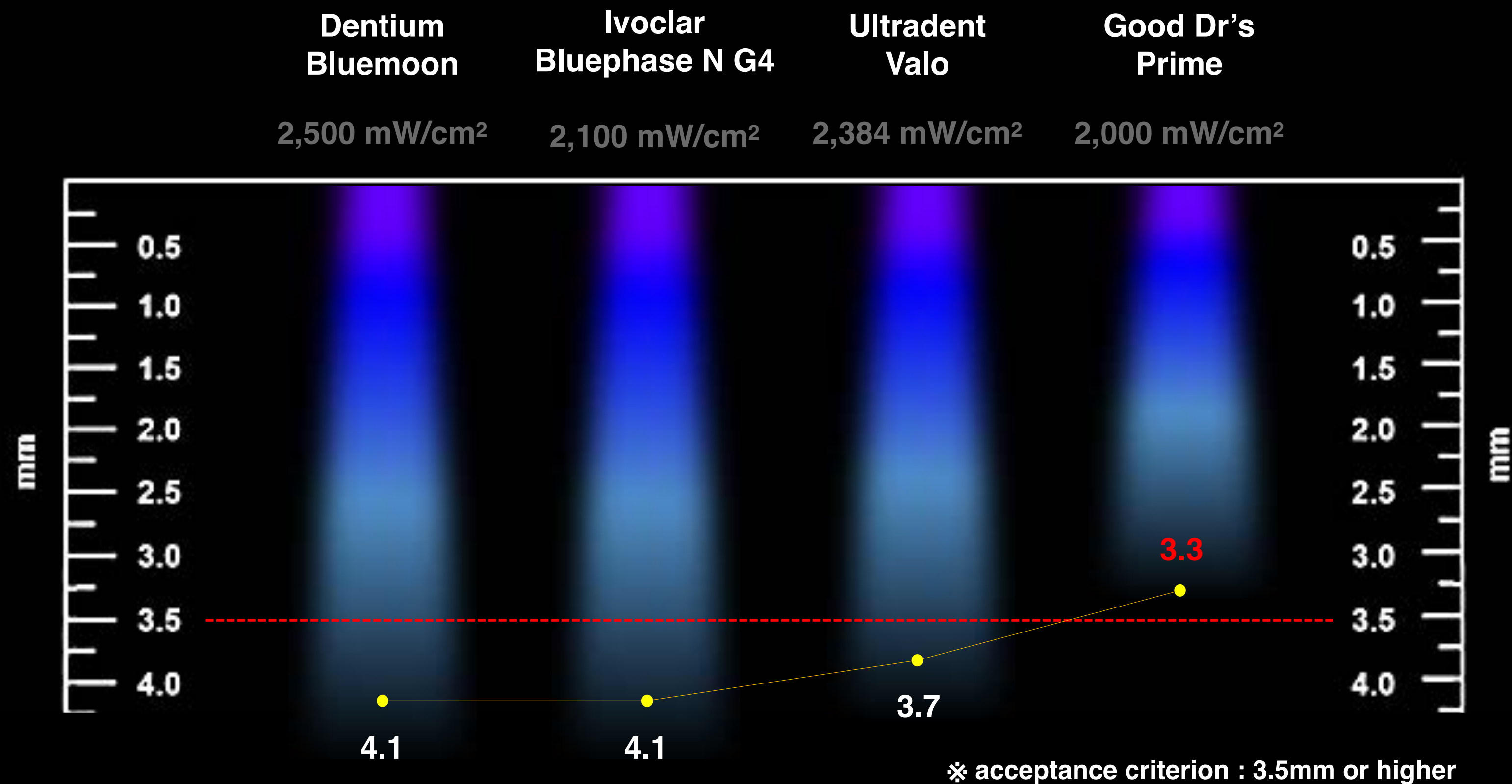



방문해 주신 여러분들께 감사합니다

자료 준비에 도움을 준  
개발팀 치과팀  
임상연구팀 마케팅팀  
노고에 감사합니다